

Operation Manual Multi-parameter Transmitter M400



Operation Manual Multi-parameter Transmitter M400

Content

| | | |
|----------|--|-----------|
| 1 | Introduction | 9 |
| 2 | Safety instructions | 10 |
| 2.1 | Definition of equipment and documentation symbols and designations | 10 |
| 2.2 | Correct disposal of the unit | 11 |
| 2.3 | Ex Classification | 12 |
| 2.4 | Ex instructions for M400 series multi-parameter Transmitters – FM Approval | 13 |
| 2.4.1 | Instructions of use to be considered under FM approval | 13 |
| 2.4.1.1 | General notes | 14 |
| 2.4.1.2 | Cautionary notes, warnings and markings | 14 |
| 2.5 | Ex instructions for M400 series multi-parameter transmitters | 16 |
| 3 | Unit Overview | 18 |
| 3.1 | M400 ½ DIN Versions | 18 |
| 3.2 | Menu Structure | 19 |
| 3.3 | Display | 20 |
| 3.4 | Operating Elements | 21 |
| 3.5 | Entry of Data | 21 |
| 3.6 | Selection Menus | 21 |
| 3.7 | "Save changes" Dialog | 22 |
| 3.8 | Security Passwords | 22 |
| 3.9 | Graphic Trend Measurement | 22 |
| 3.9.1 | Activation Trend Display Screen | 23 |
| 3.9.2 | Settings for Trend Display Screen | 24 |
| 3.9.3 | Deactivation Trend Display Screen | 24 |
| 4 | Installation instruction | 25 |
| 4.1 | Unpacking and inspection of equipment | 25 |
| 4.2 | Mounting ½ DIN Versions | 25 |
| 4.2.1 | Dimensions ½ DIN Version | 25 |
| 4.2.2 | Mounting Procedure – ½ DIN Version | 26 |
| 4.2.3 | ½ DIN – Panel Mounting | 27 |
| 4.2.4 | ½ DIN Version – Wall Mounting | 28 |
| 4.2.5 | ½ DIN Version – Pipe Mounting | 29 |
| 4.3 | Electrical Connection | 29 |
| 4.4 | Terminal Definition | 30 |
| 4.4.1 | TB1 Terminal definition – All transmitter versions | 31 |
| 4.4.2 | TB2 Terminal definition | 31 |
| 4.4.3 | TB3 Terminal definition – Analog sensors | 32 |
| 4.4.4 | TB3 Terminal definition – ISM sensors | 34 |
| 5 | Placing transmitter in, or out, of service | 35 |
| 5.1 | Placing transmitter in service | 35 |
| 5.2 | Placing transmitter out of service | 35 |
| 6 | Calibration | 36 |
| 6.1 | Sensor Calibration | 36 |
| 6.1.1 | Select the desired sensor calibration task | 36 |
| 6.1.2 | Terminate Sensor Calibration | 37 |
| 6.2 | Calibration of UniCond 2-e and UniCond 4-e Sensors (ISM Sensors only) | 37 |
| 6.2.1 | Conductivity Calibration of UniCond 2-e and UniCond 4-e Sensors | 37 |
| 6.2.1.1 | One-Point Calibration | 39 |
| 6.2.1.2 | Two-Point Calibration | 40 |
| 6.2.1.3 | Process Calibration | 41 |
| 6.2.2 | Temperature Calibration of UniCond 2-e Sensors and UniCond 4-e Sensors | 42 |
| 6.2.2.1 | One-Point Calibration | 42 |
| 6.2.2.2 | Two-Point Calibration | 43 |
| 6.3 | Calibration of Cond2e Sensors or Cond4e Sensors | 45 |
| 6.3.1 | One-Point Calibration | 45 |
| 6.3.2 | Two-Point Calibration | 46 |
| 6.3.3 | Process Calibration | 46 |
| 6.4 | pH Calibration | 47 |
| 6.4.1 | One-Point Calibration | 47 |
| 6.4.2 | Two-Point Calibration | 48 |
| 6.4.3 | Process Calibration | 48 |
| 6.5 | ORP Calibration of pH Sensors | 49 |

| | | |
|----------|---|-----------|
| 6.6 | Calibration of Amperometric Oxygen Sensors | 49 |
| 6.6.1 | One-Point Calibration | 50 |
| 6.6.2 | Process Calibration | 51 |
| 6.7 | Calibration of Optical Oxygen Sensors (ISM Sensors only) | 51 |
| 6.7.1 | One-Point Calibration | 52 |
| 6.7.2 | Two-Point Calibration | 52 |
| 6.7.3 | Process Calibration | 53 |
| 6.8 | Calibration of Dissolved Carbon Dioxide Sensors (ISM Sensors only) | 54 |
| 6.8.1 | One-Point Calibration | 54 |
| 6.8.2 | Two-Point Calibration | 55 |
| 6.8.3 | Process Calibration | 56 |
| 6.9 | Calibration of Thermal Conductivity CO ₂ (CO ₂ high) Sensors (ISM Sensors only) | 56 |
| 6.9.1 | One-Point Calibration | 57 |
| 6.9.2 | Process Calibration | 57 |
| 6.10 | Calibration of O ₃ Sensors | 58 |
| 6.10.1 | One-Point Calibration | 58 |
| 6.10.2 | Process Calibration | 59 |
| 6.11 | Calibration of a tunable diode laser (TDL) analyzer | 60 |
| 6.11.1 | One-point calibration for TDL gas sensors | 60 |
| 6.11.2 | Process calibration for TDL gas sensors | 61 |
| 6.12 | Sensor Verification | 62 |
| 6.13 | UniCond 2-e Electronics Calibration (ISM Sensor only) | 63 |
| 6.14 | Meter Calibration (Analog Sensors only) | 63 |
| 6.14.1 | Resistance (Analog Sensors only) | 64 |
| 6.14.2 | Temperature (Analog Sensors only) | 65 |
| 6.14.3 | Voltage (Analog Sensors only) | 66 |
| 6.14.4 | Current (Analog Sensors only) | 66 |
| 6.14.5 | Rg (Analog Sensors only) | 67 |
| 6.14.6 | Rr (Analog Sensors only) | 67 |
| 6.15 | Analog Output Calibration | 67 |
| 6.16 | Analog Input Calibration | 68 |
| 6.17 | Maintenance | 68 |
| 7 | Configuration | 69 |
| 7.1 | Measurement | 69 |
| 7.1.1 | Channel Setup | 69 |
| 7.1.2 | Analog sensor | 69 |
| 7.1.3 | ISM sensor | 70 |
| 7.1.4 | Parameter Related Settings | 71 |
| 7.1.4.1 | Conductivity Settings | 71 |
| 7.1.4.2 | pH Settings | 72 |
| 7.1.4.3 | Settings for Oxygen Measurement Based on Amperometric Sensors | 74 |
| 7.1.4.4 | Settings for Oxygen Measurement Based on Optical Sensors | 75 |
| 7.1.4.5 | Dissolved Carbon Dioxide Settings | 76 |
| 7.1.4.6 | Settings for Thermal Conductivity Dissolved CO ₂ Measurement (CO ₂ hi) | 77 |
| 7.1.4.7 | Settings for Tunable Diode Laser (TDL) Analyzer | 78 |
| 7.1.4.8 | Setting the correct process side purging | 79 |
| 7.1.5 | Concentration Curve Table | 80 |
| 7.2 | Temperature Source (Analog Sensors only) | 80 |
| 7.3 | Analog Outputs | 81 |
| 7.4 | Set Points | 82 |
| 7.5 | ISM Setup (ISM Sensors only) | 83 |
| 7.5.1 | Sensor Monitor | 83 |
| 7.5.2 | CIP Cycle Limit | 85 |
| 7.5.3 | SIP Cycle Limit | 85 |
| 7.5.4 | AutoClave Cycle Limit | 86 |
| 7.5.5 | DLI Stress Adjustment | 87 |
| 7.5.6 | SAN Cycle Parameters | 87 |
| 7.5.7 | Reset Counters for UniCond 2-e Sensors | 88 |
| 7.5.8 | Set Calibration Interval for UniCond 2-e Sensors | 88 |
| 7.6 | General Alarm | 88 |
| 7.7 | ISM / Sensor Alarm | 89 |
| 7.8 | Clean | 89 |
| 7.9 | Display Setup | 90 |
| 7.10 | Digital Inputs | 90 |
| 7.11 | System | 91 |
| 7.12 | PID Controller | 92 |

| | | |
|-----------|--|------------|
| 7.13 | Service | 96 |
| 7.13.1 | Set Analog Outputs | 96 |
| 7.13.2 | Read Analog Outputs | 96 |
| 7.13.3 | Set Relay | 96 |
| 7.13.4 | Read Relay | 96 |
| 7.13.5 | Read Digital Inputs | 96 |
| 7.13.6 | Memory | 97 |
| 7.13.7 | Display | 97 |
| 7.13.8 | Calibrate TouchPad | 97 |
| 7.13.9 | Channel Diagnostic | 97 |
| 7.14 | User Management | 98 |
| 7.15 | Reset | 99 |
| 7.15.1 | System Reset | 99 |
| 7.15.2 | Reset Sensor Calibration for UniCond 2-e Sensors | 99 |
| 7.16 | USB | 100 |
| 7.16.1 | Printer Output Configuration | 100 |
| 7.16.2 | USB data logging | 101 |
| 7.17 | Configuration via USB | 102 |
| 8 | ISM | 103 |
| 8.1 | iMonitor | 103 |
| 8.2 | Messages | 104 |
| 8.3 | ISM Diagnostics | 104 |
| 8.3.1 | pH/ORP, Oxygen, O ₂ , Cond4e Sensors and TDL | 105 |
| 8.3.2 | UniCond 2-e and UniCond 4-e Sensors | 105 |
| 8.4 | Calibration Data | 106 |
| 8.4.1 | Calibration Data for All ISM Sensors excluding UniCond 2-e and UniCond 4-e | 106 |
| 8.4.2 | Calibration Data for UniCond 2-e and UniCond 4-e Sensors | 107 |
| 8.5 | Sensor Info | 107 |
| 8.6 | HW / SW Version | 108 |
| 9 | Custom Key | 109 |
| 9.1 | Set Favorite | 109 |
| 10 | Maintenance | 110 |
| 10.1 | Front panel cleaning | 110 |
| 11 | Software History | 110 |
| 11.1 | M400 Type 1 | 110 |
| 11.2 | M400 Type 2 | 110 |
| 11.3 | M400 Type 3 | 110 |
| 11.4 | M400 4-wire FF | 110 |
| 12 | Troubleshooting | 111 |
| 12.1 | Conductivity (resistive) Error messages/ Warning- and Alarm list for analog sensors | 111 |
| 12.2 | Conductivity (resistive) Error messages/ Warning- and Alarm list for ISM sensors | 112 |
| 12.3 | pH Error messages/Warning- and Alarm list | 112 |
| 12.3.1 | pH, pH/pNa and Dissolved Carbon Dioxide sensors | 112 |
| 12.3.2 | ORP messages | 113 |
| 12.4 | Amperometric O ₂ Error messages/ Warning- and Alarm list | 114 |
| 12.4.1 | High level oxygen sensors | 114 |
| 12.4.2 | Low level oxygen sensors | 114 |
| 12.4.3 | Trace oxygen sensors | 115 |
| 12.5 | Warning- and Alarm Indication | 116 |
| 12.5.1 | Warning Indication | 116 |
| 12.5.2 | Alarm Indication | 117 |
| 13 | Ordering Information, Accessories and Spare Parts | 118 |
| 14 | Specifications | 119 |
| 14.1 | General specifications | 119 |
| 14.2 | Electrical specifications | 122 |
| 14.3 | FOUNDATION fieldbus specifications | 122 |
| 14.4 | Environmental specifications | 123 |
| 14.5 | Mechanical specifications | 123 |
| 15 | Warranty | 124 |

| | | |
|-----------|--|------------|
| 16 | Buffer tables | 125 |
| 16.1 | Standard pH buffers | 125 |
| 16.1.1 | Mettler-9 | 125 |
| 16.1.2 | Mettler-10 | 126 |
| 16.1.3 | NIST Technical Buffers | 126 |
| 16.1.4 | NIST standard buffers (DIN and JIS 19266: 2000-01) | 127 |
| 16.1.5 | Hach buffers | 127 |
| 16.1.6 | Ciba (94) buffers | 128 |
| 16.1.7 | Merck Titrisole, Riedel-de-Haën Fixanale | 128 |
| 16.1.8 | WTW buffers | 129 |
| 16.1.9 | JIS Z 8802 buffers | 129 |
| 16.2 | Dual membrane pH electrode buffers | |
| 16.2.1 | Mettler-pH/pNa buffers (Na+ 3.9M) | |

1 Introduction

The M400 is a 4-wire transmitter, with 4(O) to 20 mA output signal and HART or FOUNDATION Fieldbus communication capabilities, for analytical measurements. The M400 is a multi-parameter transmitter and supports the measurements listed in the parameter fit guide below.

The M400 transmitter is designed for use in the process industries.

M400 parameter fit guide

| | M400 Type 1 | | M400 Type 2/ M400 4-wire FF | | M400 Type 3 | |
|--|-------------|-----|--------------------------------|----------------------|-------------|-------|
| | Analog | ISM | Analog | ISM | Analog | ISM |
| pH/ORP | • | • | • | • | • | • |
| pH/pNa | – | • | – | • | – | • |
| UniCond 2-e/4-e | – | • | – | • | – | • |
| Conductivity 2-e | • | – | • | – | • | – |
| Conductivity 4-e | • | • | • | • | • | • |
| Amp. dissolved oxygen ppm/ppb/trace | – | – | •/• ¹⁾ /– | •/• ¹⁾ /– | •/•/• | •/•/• |
| Opt. dissolved oxygen ppm/ppb | – | – | –/– | •/• ²⁾ | –/– | •/• |
| Amp. O ₂ gas ppm/ppb/trace | – | – | –/–/– | –/–/– | •/•/• | •/•/• |
| Opt. O ₂ gas ppm | – | – | – | – | – | • |
| Dissolved ozone | – | – | • | • | • | • |
| Dissolved carbon dioxide | – | – | • | • | • | • |
| CO ₂ hi | – | – | – | – | – | • |
| GPro 500 TDL | – | – | – | – | – | • |

1) M400 4-wire FF supports Ingold Amp. DO ppb sensors.

2) Thornton high performance dissolved oxygen and pure water optical sensors only.

A black & white touch screen conveys measuring data and setup information. The menu structure allows the operator to modify all operational parameters. A menu-lockout feature, with password protection, is available to prevent the unauthorized use of the meter. The M400 Multi-parameter transmitter can be configured to use up to four analog and/or up to four relay outputs as well as HART communication protocol for process control.

The M400 Multi-parameter transmitter is equipped with a USB communication interface. This interface provides up- and download capabilities of the transmitter configuration via a Personal Computer (PC).

This description corresponds to the firmware release, version 1.0. Changes are taking place constantly, without prior notification.

2 Safety instructions

This manual includes safety information with the following designations and formats.

2.1 Definition of equipment and documentation symbols and designations



WARNING: POTENTIAL FOR PERSONAL INJURY.



CAUTION: Possible instrument damage or malfunction.



NOTE: Important operating information.



On the transmitter or in this manual text indicates: Caution and/or other possible hazard including risk of electric shock (refer to accompanying documents).

The following is a list of general safety instructions and warnings. Failure to adhere to these instructions can result in damage to the equipment and/or personal injury to the operator.

- The M400 Transmitter should be installed and operated only by personnel familiar with the transmitter and who are qualified for such work.
- The M400 Transmitter must only be operated under the specified operating conditions (see chapter “14 Specifications” on page 119).
- Repair of the M400 Transmitter must be performed by authorized, trained personnel only.
- With the exception of routine maintenance, cleaning procedures, as described in this manual, the M400 Transmitter must not be tampered with or altered in any manner.
- Mettler-Toledo accepts no responsibility for damage caused by unauthorized modifications to the transmitter.
- Follow all warnings, cautions, and instructions indicated on and supplied with this product.
- Install equipment as specified in this instruction manual. Follow appropriate local and national codes.
- Protective covers must be in place at all times during normal operation.
- If this equipment is used in a manner not specified by the manufacturer, the protection provided by it against hazards may be void.



WARNINGS:

- Installation of cable connections and servicing of this product require access to shock hazard voltage levels.
- Main power and relay contacts wired to a separate power source must be disconnected before servicing.
- Switch or circuit breaker shall be in close proximity to the equipment and within easy reach of the OPERATOR; it shall be marked as the disconnecting device for the equipment.
- Main power must employ a switch or circuit breaker as the disconnecting device for the equipment.
- Electrical installation must be in accordance with the National Electrical Code and/or any other applicable national or local codes.

**NOTE: RELAY CONTROL ACTION**

the M400 Transmitter relays will always de-energize on loss of power, equivalent to normal state, regardless of relay state setting for powered operation. Configure any control system using these relays with fail-safe logic accordingly.

**NOTE: PROCESS UPSETS**

Because process and safety conditions may depend on consistent operation of this transmitter, provide appropriate means to maintain operation during sensor cleaning, replacement, or sensor or instrument calibration.



NOTE: This is a 4-wire-product with an active 4–20 mA analog output.

Do not supply power to the analog output terminals (TB2: terminal 1 to 8).

2.2 Correct disposal of the unit

When the transmitter is finally removed from service, observe all local environmental regulations for proper disposal.

2.3 Ex Classification



NOTE: The Ex classification is valid for the transmitters M400 Type 1, M400 Type 2 and M400 Type 3.

| | | |
|------------------|-------------------------------|--------------------------|
| Standards | CSA Std C22.2 No. 213-16; | UL 60079-0-2013 |
| | CAN/CSA-C22.2 No. 60079-0-15 | UL 60079-15-2013 |
| | CAN/CSA-C22.2 No. 60079-15-16 | EN 60079-0:2012/A11:2013 |
| | ANSI/ISA-12.12.01-2016 | EN 60079-15:2010 |

Special condition of safe use

1. This equipment uses external non-metallic components, therefore may generate an ignition-capable level of electrostatic charge under certain extreme conditions. The user should ensure that the equipment is not installed in a location where it may be subjected to external conditions (such as high pressure steam) which might cause a build up of electrostatic charge on non-conducting surfaces.
2. The display has not been tested for resistance to ultraviolet light. The display shall be protected from direct light (e.g. from sunlight or luminaires).



- **WARNING**
THE EQUIPMENT IS SUITABLE FOR USE IN CLASS I, DIVISION 2, GROUPS A, B, C, D OR NON-HAZARDOUS LOCATIONS ONLY;
- **WARNING**
EXPLOSION HAZARD – DO NOT REMOVE OR REPLACE LAMPS, FUSES OR PLUG-IN MODULES (AS APPLICABLE) UNLESS POWER HAS BEEN DISCONNECTED OR THE AREA IS FREE OF IGNITIBLE CONCENTRATIONS.
- **WARNING**
EXPLOSION HAZARD. DO NOT CONNECT OR DISCONNECT WHILE THE CIRCUIT IS LIVE OR UNLESS THE AREA IS FREE OF IGNITIBLE CONCENTRATIONS;
- **WARNING**
DO NOT OPEN WHEN ENERGISED.
- **WARNING**
THIS EQUIPMENT IS DESIGNED TO BE USED IN RESTRICTED ACCESS LOCATION. ONLY SERVICE PERSONS OR TRAINED PERSONS ARE ALLOWED TO ACCESS THIS EQUIPMENT.

Mettler Toledo's M400 G2 series transmitters have been approved by FM.

Shall you need any more information, please contact:

process.service@mt.com

2.4 Ex instructions for M400 series multi-parameter Transmitters – FM Approval

2.4.1 Instructions of use to be considered under FM approval



M400 series multi-parameter transmitters are produced by Mettler-Toledo GmbH.

| US marking | |
|-----------------------------|---|
| Operating temperature range | –20 °C to +50 °C (–4 °F to +122 °F) |
| Environmental designation | Enclosure type 4X, 1P 66 |
| Nonincendive | – Class 1, Division 2, Groups A, B, C, D T4 – Class 1, Zone 2, Group 11C T4 |
| Certificate no. | FM17US0240X |
| Standards | <ul style="list-style-type: none"> – FM3810:2018 Approval Standard for Electrical Equipment for Measurement, Control and Laboratory Use. – FM3611:2018 Approval Standard for Nonincendive Electrical Equipment for Use in Class 1 & 11, Division 2, and Class 111, Division 1 & 2, Hazardous (Classified) Locations. – FM3600:2018 Approval Standard for Electrical Equipment for Use in Hazardous (Classified) Locations – General Requirements. – ANS1/IEC 60529-2004: R2011 Degrees of Protection Provided by Enclosures (1P Codes). – ANS1/UL 121201: 2017 Nonincendive Electrical Equipment for Use in Class 1 & 2, Division 2, and Class 111, Divisions 1 & 2, Hazardous (Classified) Locations. – ANS1/UL 61010-1: 2016 Safety Requirements for Electrical Equipment for Measurement, Control, and Laboratory Use. – Part 1: General Requirements – ANS1/UL 50E: 2015 Enclosures for Electrical Equipment, Environmental Considerations |

2.4.1.1 General notes

The Multi-parameter Transmitter M400 Type 1, 2, 3 are suitable for use in hazardous atmospheres of all combustible materials of explosion groups A, B, C, D for applications requiring Class 1, Division 2 and Class 1, Zone 2 instruments (National Electrical Code® [ANSI/NFPA 70 (NEC®)], Article 500.

If the Multi-parameter Transmitter M400 Type 1, 2, 3 are installed and operated in hazardous areas, the general Ex installation regulations as well as these safety instructions must be observed.

The operating instructions as well as the installation regulations and standards that apply for explosion protection of electrical systems must always be observed.

The installation of explosion-endangered systems must always be carried out by qualified personnel.

For mounting instructions on specific valves refer to the mounting instructions supplied with the mounting kit. Mounting does not affect the suitability of the transmitter for use in a potentially hazardous environment.

The equipment is not intended to be used as personal protective equipment. To prevent injury, read the manual before use.

For language translation assistance contact your local representative or email to process.service@mt.com

2.4.1.2 Cautionary notes, warnings and markings

Hazardous location notes:

1. Installations in the US shall comply with the relevant requirements of the National Electrical Code® (ANSI/NFPA 70 [NEC®]), Article 500.
2. Installations shall comply with the latest edition of the manufacturer's instruction manual.
3. Care must be taken during installation to avoid impacts or friction that could create an ignition source.
4. Extreme care should be taken with the installation of his equipment, and any problems should be resolved by consultation with the factory or the authorized representative.
5. For equipment marked with an environmental ingress protection rating, provision shall be made for maintaining that environmental ingress protection rating whether the cable assembly is connected to the equipment or not, within the enclosure and within the connector body.
6. Any plugs and sockets incorporated as part of the equipment apparatus shall be capable of being connected to a wiring method as permitted by the National Electrical Code® (ANSI/NFPA 70 [NEC®]) for the involved hazardous (classified) location in accordance with the associated location restrictions.
7. The internal grounding terminal shall be used as the primary equipment grounding means and the external grounding terminal is only for a supplemental (secondary) bonding connection where local authorities permit or require such a connection.
8. Tighten cover screws to 2.5 N·m (22 lb·in.) maximum. Over-torquing may cause enclosure breakage.
9. The minimum tightening torque for M4 (No. 6) binding screw protective conductor terminals is 1.2 N·m (10.6 lb·in.) or greater, as specified.
10. Use copper, copper-clad aluminum or aluminum conductors only.

11. For ambient temperatures above +40 °C (+104 °F) use installation wiring connection suitable for maximum ambient temperatures, as prescribed by the manufacturer.
12. The multi-parameter transmitter must be connected to limited output NEC Class 2 circuits, as outlined in the National Electrical Code® (ANSI/NFPA 70 [NEC®]) only. If the devices are connected to a redundant power supply (two separate power supplies), both must meet this requirement
13. The Class I, Zone 2 certifications are based on Division evaluations and the marking acceptance of Article 505 of the National Electrical Code® (ANSI/NFPA 70 [NEC®]).
14. Tampering and replacement with non-factory components may adversely affect the safe use of the system.
15. The multi-parameter transmitter is intended for servicing or maintenance operation. Malfunctioning units operating out of manufacturer's specification should be returned to the authorized service center for repair. On-site repairs are not permitted.
16. If the equipment is installed via panel-mount configuration within an ultimate enclosure, the inner service temperature of the enclosure corresponds to the ambient temperature of the module.
17. If the panel-mount configuration module is operated in an ambient temperature between +40 °C and +50 °C, the temperature of the module housing may be higher than +50 °C. The device must therefore be installed to that it is only accessible to service personnel or users that are aware of the reason for restricted access and the required safety measures at an ambient temperature of +40 °C to +50 °C.
18. Insertion or withdrawal of removable electrical connectors or modules is to be accomplished only when the area is known to be free of flammable vapors.
- 19. WARNING** – POTENTIAL ELECTROSTATIC CHARGING HAZARD – SEE INSTRUCTIONS.
- 20. WARNING** – SUBSTITUTION OF COMPONENTS MAY IMPAIR SUITABILITY FOR DIVISION 2.
- 21. WARNING** – DO NOT REMOVE OR REPLACE WHILE CIRCUIT IS LIVE WHEN A FLAMMABLE OR COMBUSTIBLE ATMOSPHERE IS PRESENT.
- 22. WARNING** – EXPLOSION HAZARD, DO NOT DISCONNECT EQUIPMENT WHEN A FLAMMABLE OR COMBUSTIBLE ATMOSPHERE IS PRESENT.
- 23. WARNING** – FOR CONNECTION ONLY TO NON-FLAMMABLE PROCESSES.
- 24. WARNING** – SUBSTITUTION OF COMPONENTS MAY IMPAIR SUITABILITY OF THE EQUIPMENT.
- 25. CAUTION** – TO PREVENT INJURY, READ THE MANUAL BEFORE USE.
- 26. WARNING** – TO MAINTAIN THE ENCLOSURE TYPE AND INGRESS PROTECTION RATING, THE COVER MUST BE CLOSED AND SECURED.

2.5 Ex instructions for M400 series multi-parameter transmitters

M400 series multi-parameter transmitters are produced by Mettler-Toledo GmbH. It has passed the inspection of IECEx and conforms to following standards:

- **IEC 60079-0: 2017 Edition:**
7.0 Explosive atmospheres – Part 0: General requirements
- **IEC 60079-11: 2011 Edition:**
6.0 Explosive atmospheres – Part 11: Equipment protection by intrinsic safety "i"
- **IEC 60079-15: 2017 Edition:**
5.0 Explosive atmospheres – Part 15: Equipment protection by type of protection "n"
- **IEC 60079-7: 2015 Edition:**
5.0 Explosive atmospheres – Part 7: Equipment protection by increased safety "e"

Ex Marking: Ex ec ic nC IIC T4 Gc
Certificate No.: IECEx NEP 19.0008X

Rated ambient temperature range: –20 ~ +50 °C
Um = 253 Vac

Special Conditions of use (X-marking in the Certificate Number):.

1. Avoid electrostatic discharge on enclosure surface, use wet cloth only for cleaning.
2. The display shall be protected from direct light (e.g. from sunlight or luminaires).
3. Take protective measure to avoid risk of mechanical danger "high" on the display.
4. When installation in explosive atmosphere, cable gland separately certified according to IEC 60079-0:2017 and IEC 60079-7:2015 with marking Ex ec IIC IP66 shall be adopted.
5. This equipment shall only be used in an area of at least pollution degree 2, as defined in IEC 60664-1.
6. Observe the warnings:
DO NOT CONNECT OR DISCONNECT WHILE CIRCUIT IS LIVE UNLESS AREA IS KNOWN TO BE NON-HAZARDOUS! DO NOT OPEN WHILE ENERGIZED!
Potential electrostatic charging hazard — see instructions!
7. While installation, use and maintenance, IEC 60079-14 shall be observed.
8. The equipment is provided with external earthing facilities (M4) at the bottom, which is suitable for connection lug.



Earth connection cable with lug

The Multi-parameter Transmitter M400 Type 1, 2, 3 non-incentive version, bears the following label marking :

METTLER TOLEDO

Do not connect or disconnect USB while live unless an explosive atmosphere is shown to be absent.

Analog Sensor Input,
See Manual

| | | | | | | | | | | | | | | | | | |
|---|---|---|---|---|---|---|---|---|----|----|----|----|----|----|----|----|----|
| 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 | 13 | 14 | 15 | 16 | 17 | 18 |
|---|---|---|---|---|---|---|---|---|----|----|----|----|----|----|----|----|----|

TB3 (Sensor Input 1)

1-Wire
GND
RS485B
RS485A
5V
GND_24V
24V

Always Connect Protective Earth!

L(+) N(-) PE

Power

| | | | | | | | | | |
|------|-----|------|-----|-----|------|-----|-----|------|-----|
| COM4 | NO4 | COM3 | NO3 | NC2 | COM2 | NO2 | NC1 | COM1 | NO1 |
| 10 | 9 | 8 | 7 | 6 | 5 | 4 | 3 | 2 | 1 |

TB1 (Relays)

MODEL: M400 Type *
 UL File # E480774
 PWR: 100-240VAC 50/60 Hz, 10W
 20-30VDC, 10W
 T_[amb]: -20 to +50 °C

| | | | | | | | | | | | | | | | |
|----|----|----|----|----|----|----|---|---|---|---|---|---|---|---|---|
| 16 | 15 | 14 | 13 | 12 | 11 | 10 | 9 | 8 | 7 | 6 | 5 | 4 | 3 | 2 | 1 |
|----|----|----|----|----|----|----|---|---|---|---|---|---|---|---|---|

TB2 (I/O Terminal)

Ain+ Ain- D1+ D1-/D2-
D1+ A04- A04+ A03- A03+ A02- A02+ A01- /HART-
A01+ /HART+

P/N: S/N:

UL US LISTED
Measuring Equipment
26MR

MODEL: M400 Type *
 UL File # E480774
 PWR: 100-240VAC 50/60 Hz, 10W
 20-30VDC, 10W
 T_[amb]: -20 to +50 °C ENCL: IP66 / NEMA 4X

P/N: S/N:

METTLER TOLEDO
 www.mt.com/pro
 Mettler-Toledo GmbH
 Im Hackacker 15
 8902 Urdorf, Switzerland

UL US LISTED
Measuring Equipment
26MR

Made by METTLER TOLEDO in China

IECEx

SP
c us

APPROVED
-20 °C ≤ T_a ≤ +50 °C

Ex ec to nC IIC T4 Gc
 IECEx NEP 19.0008X
 Um=253VAC
 Sira 16 ATEX 4341X
 II 3G Ex nA nC IIC T4 Gc
 CSA File 18.70100658X ENCL: IP66 Type 4X
 Class I, Division 2, Groups A, B, C and D T4;
 Ex nA nC IIC T4 Gc;
 Class I, Zone 2, AEx nA nC IIC T4 Gc
 FM17US0240X ENCL IP66 Type 4X
 Class I, Division 2, Groups A, B, C and D T4;
 Class I, Zone 2, Group IIC T4

WARNING
 EXPLOSION HAZARD - DO NOT CONNECT OR DISCONNECT WHILE CIRCUIT IS LIVE UNLESS AREA IS KNOWN TO BE NON-HAZARDOUS
 POTENTIAL ELECTROSTATIC CHARGING HAZARD - SEE INSTRUCTIONS
 DO NOT OPEN WHILE ENERGIZED
AVERTISSEMENT
 RISQUE D'EXPLOSION. NE PAS DEBRANCHER TANT QUE LE CIRCUIT EST SOUS TENSION, A MOINS QU'IL NE S'AGISSE D'UN EMPLACEMENT NON - DANGEREUX
 NE PAS OUVRIR SOUS TENSION
 DANGER POTENTIEL DE CHARGES ELECTROSTATIQUES - VOIR INSTRUCTIONS

Entry thread: Metric, 4xM20, 1xM25
 NOTE:
 1. Conduit Hub/Fittings Entry Thread;
 2. Must use minimum Class I, Division 2, Groups A, B, C, D, Type 4X and IP66 suitable Hub/Fittings & Cable Glands to fulfill the complete FM certification.
 Operation Manual No. 30413330

3 Unit Overview

The M400 transmitter is available as ½ DIN version.

For dimensions refer to chapter "13 Ordering Information, Accessories and Spare Parts" on page 118.

3.1 M400 ½ DIN Versions

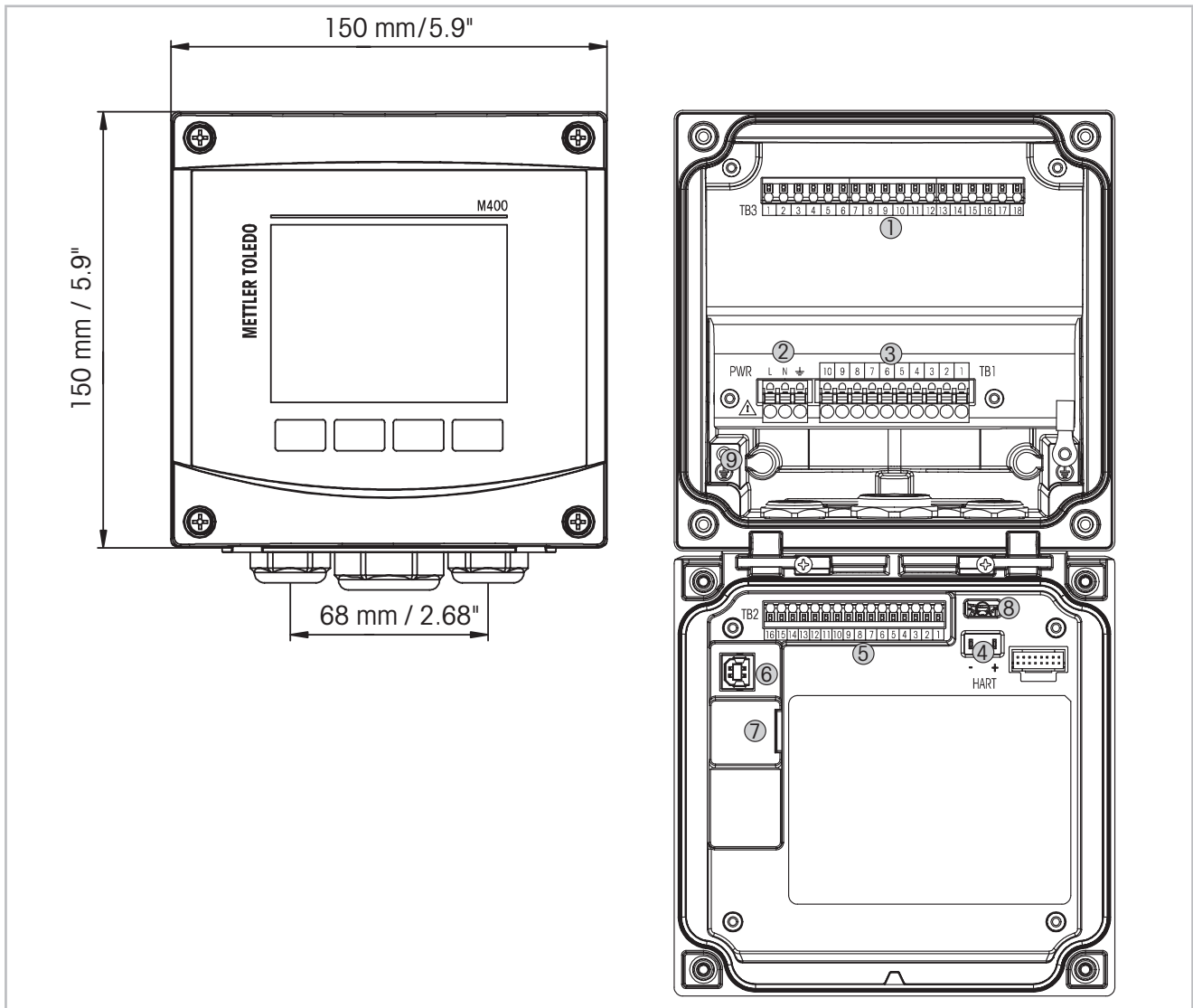


Fig. 1: M400 ½ DIN versions

- 1 TB3 – Terminal block for sensor connection
- 2 Terminals for supply voltage
- 3 TB1 – Terminal block for relay outputs
- 4 HART, for M400 Type 1, Type 2, Type 3 only
- 5 TB2 – Terminal block for analog output and digital input signals
- 6 USB Device – Software update interface
- 7 USB Host – Printer connection, data logging, loading and saving configuration
- 8 **Warning! Do not disconnect the internal ground wire between the front and back modules.**
- 9 **Warning! Tightly secure an earth wire to the internal PE screw terminal:**
 ↓ (Protective Conductor Terminal).
The cross-section of the PE wire must be above 18 AWG (0.8 mm).



3.2 Menu Structure

Below is the structure of the M400 menu tree:

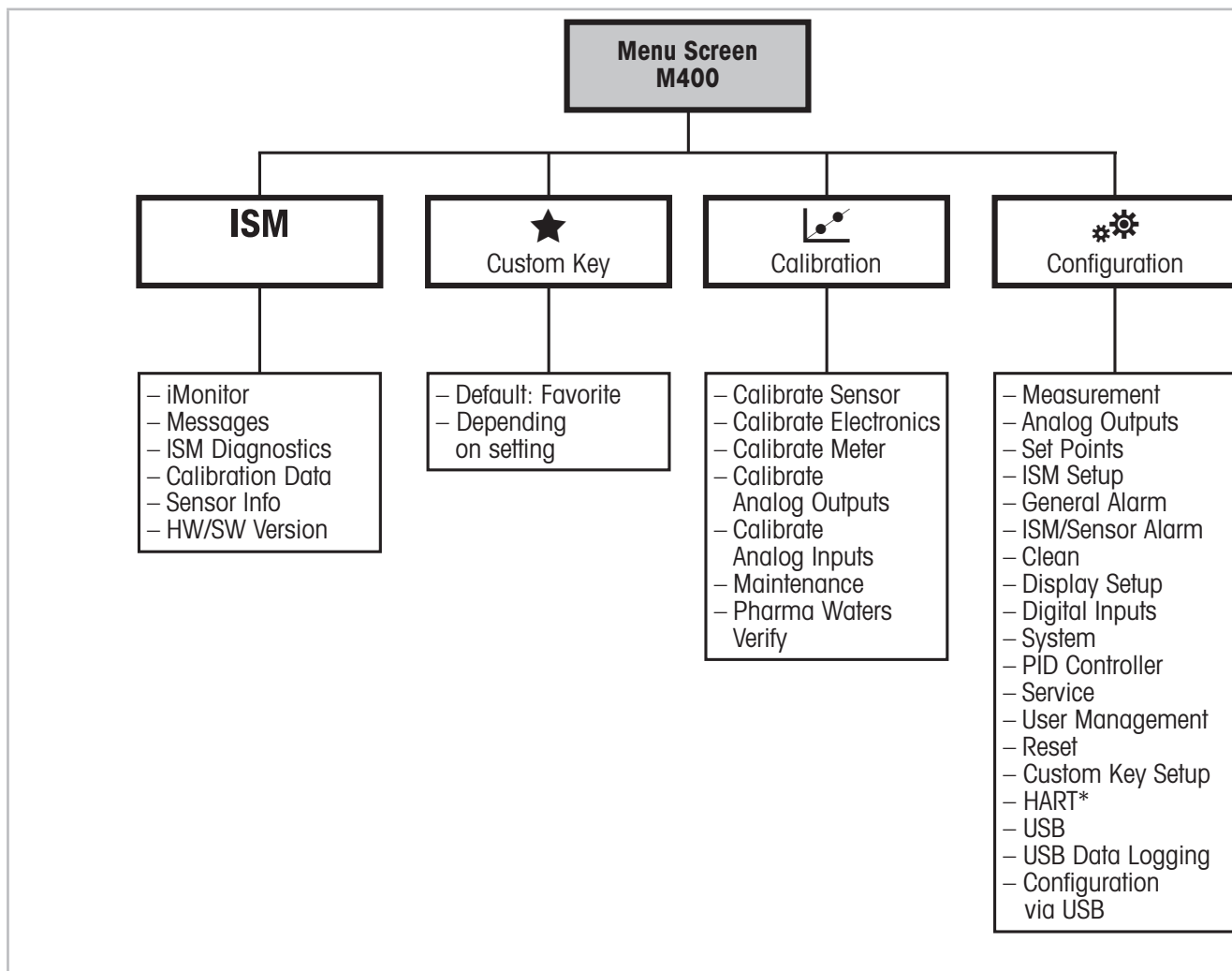


Fig. 2: Menu overview

* For Type 1, Type 2 and Type 3 only

3.3 Display

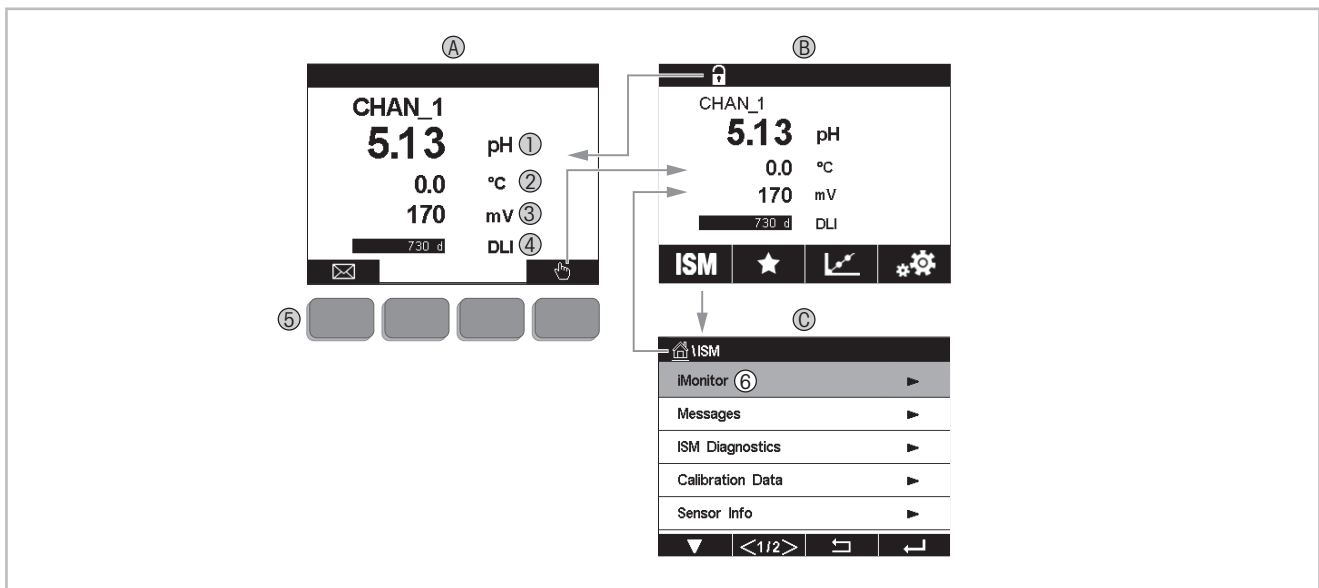


Fig. 3: M400 Display, navigation

A Start screen (example)

- 1 1st line, standard configuration
- 2 2nd line, standard configuration
- 3 3rd line, depends on configuration
- 4 4th line, depends on configuration
- 5 Soft key with indicated functions on the screen
- 6 Cursor, indicates the current item for soft key operation

B Menu screen (example)

C ISM Menu screen














NOTE: In the event of an alarm or other error conditions the M400 Transmitter will display a symbol in the head line of the display. This head line is blinking until the condition that caused it has been cleared (see chapter "12.5 Warning- and Alarm Indication" on page 116).




NOTE: During calibrations, clean, Digital In with Analog Output/Relay/USB in HOLD state, a flashing "H" (HOLD) will appear in the upper right corner of the display for the corresponding channel. This symbol will remain for 20 sec., after end of calibration. This symbol will remain for 20 seconds until after the calibration or clean is completed. This symbol will also disappear when Digital In is deactivated.

3.4 Operating Elements

| Operating element | Description |
|---|--|
|  | Enter Messages menu |
|  | Enter Menu screen |
|  | Enter Start screen |
| ISM | Enter ISM menu |
|  | Enter Favorite menu |
|  | Enter Calibration menu |
|  | Enter Configuration menu |
|  | Return to Menu screen |
|  | Enter next-lower menu level, here e.g. iMonitor, Messages or ISM Diagnostics |
|  | Return to next-higher menu level |
|  | Navigate menu for soft key operation |
|  | Enter selected menu or item for soft key operation |

3.5 Entry of Data

The M400 displays a keypad for modifying values. Press the  button and the transmitter will store the value. Press the ESC button to exit the keypad without changing data.



NOTE: For some values, the units can be modified. In this case the keypad shows a button with a U. To select another unit for the entered value on the keypad press the U button. To return again press the 0–9 button.



NOTE: For some entries letters and/or numbers can be used. In this case the keypad shows a button 'A,a,0'. Press this button to change between capital letters, small letters and numbers on the keypad.

3.6 Selection Menus

Some menus require a selection of a parameter / data. In this case the transmitter displays a pop up window. Press the according field to select the value. The pop-up window will be closed and the selection will be stored.

3.7 "Save changes" Dialog

If the M400 brings up the "Save changes" dialog there are the following options. No will discard the entered values, Yes will save changes made and Cancel will bring you back to continue configuring.

3.8 Security Passwords

The M400 Transmitter allows a security lock-out of various menus. If the security lock-out feature of the transmitter has been enabled, a security password must be entered to allow access to the menu. See chapter "7.14 User Management" on page 98.


3.9 Graphic Trend Measurement

Any single measurement may be displayed as a trend measurement over time. Measurement values will be indicated by a value on the Y-axis and time elapsed on the X-axis of the graph displayed. An actual measurement for the selected value will also be displayed numerically above the graphic trend display. The measurement value is refreshed once per second.

Graphic trending will only display the data within maximum/minimum range. Out of range values or invalid values will not be displayed. The Y-axis will display the maximum value unit with its range; X-axis unit uses "mins" for minutes for measurements less than one hour and "hrs" for one day. 4 scales for X/Y-axis. The maximum value on Y-axis is one decimal place.


3.9.1 Activation Trend Display Screen

While the M400 is displaying the Menu Screen, touch any measurement value line of the display screen once to activate the trend display for that measurement. Or you can use custom key setup to access this function when operating with tactile keys. (See chapter "9 Custom Key" on page 109.)

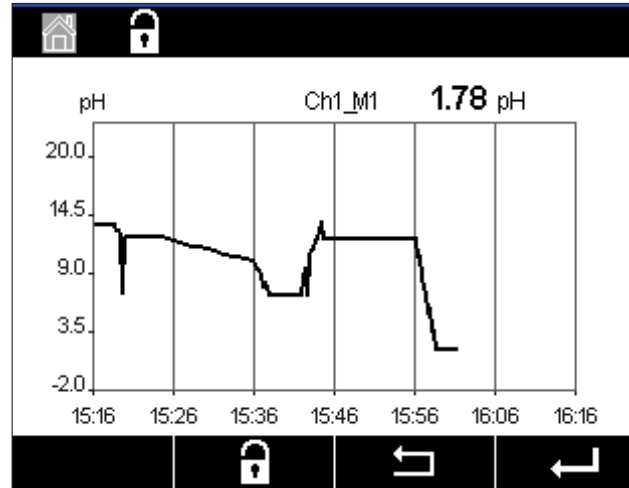

When a sensor is disconnected/connected a pop-up window will appear; after closing the window, the display will return to the Menu Screen.

The top line will display any message that occurs during trending. "H", "P" will display when this channel is in hold or process.

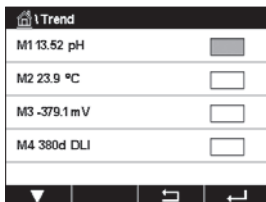


When using custom key setup to access trend display, press second left soft key after setting trend as the custom key. (See chapter "9 Custom Key" on page 109)

Use ▼ and ↵ to select the measurement.

3.9.2 Settings for Trend Display Screen

For setting configurations, touch any area of the graphic trend display to go to the pop-up window of this measurement parameter. Settings are at the default values. However, these settings may be changed when options are available, as needed.

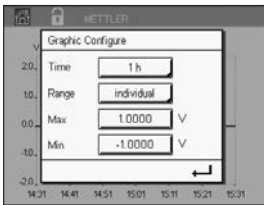

Time: Option button. For graphic display time (X-axis)
 1-h (default value)
 1-day

NOTE: 1 h means: 1 meas storage/15 seconds, totally 240 measurements for 1h. 1 day means: 1 meas storage/6 minutes, totally 240 measurements for 1 day;

Range: Option button
 Default(default value)
 Individual

When "Default" modes are set for the maximum or minimum value, this indicates the full measurement range for this unit. A Max or Min button is not displayed. If setting is selectable, the user can set maximum and minimum settings manually.

Max: Edit button.
 Maximum value of this unit on Y-axis. xxxxxx, floating decimal point.

Min: Edit button.
 Minimum value of this unit on Y-axis. xxxxxx, floating decimal point.
 Max Value > Min Value



NOTE: Settings for Y-and Y-axis and the corresponding measurement values are stored the transmitters memory. A power down returns to default settings.

3.9.3 Deactivation Trend Display Screen

Press  in activated graphic trend screen to return to Menu Screen.



NOTE: If a sensor is disconnected/connected a pop-up window come up; after closing the window, it will go back to the Menu Screen.

4 Installation instruction

4.1 Unpacking and inspection of equipment

Inspect the shipping container. If it is damaged, contact the shipper immediately for instructions. Do not discard the box.

If there is no apparent damage, unpack the container. Be sure all items shown on the packing list are present.

If items are missing, notify Mettler-Toledo immediately.

4.2 Mounting ½ DIN Versions

4.2.1 Dimensions ½ DIN Version

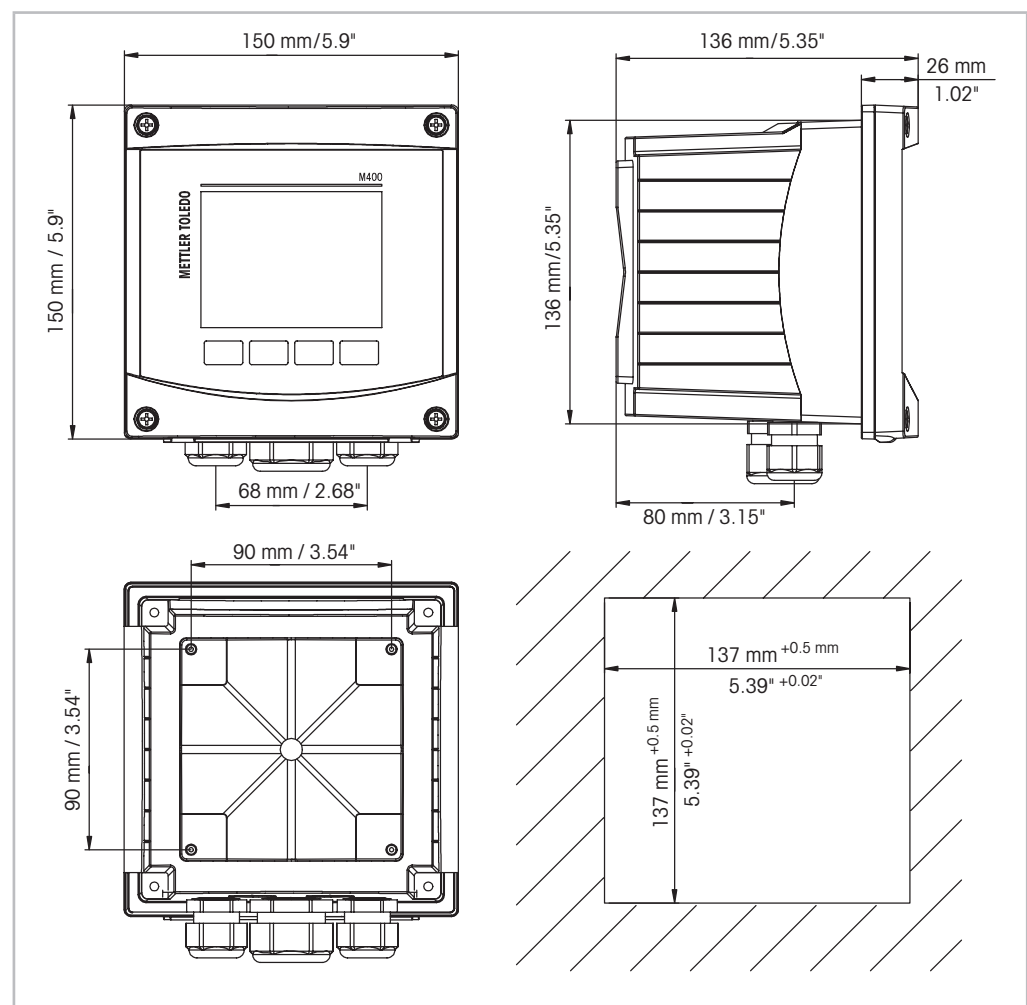


Fig. 4: Dimensions ½ DIN version

1 Dimensions for panel cutout

4.2.2 Mounting Procedure – ½ DIN Version

½ DIN versions transmitters are designed for the following mounting versions: panel mount, wall mount or pipe mount. For wall mount the integral rear cover is used.

Optional hardware accessories are available that allow for panel- or pipe-mount. Refer to section "13 Ordering Information, Accessories and Spare Parts" on page 118.

Assembly:

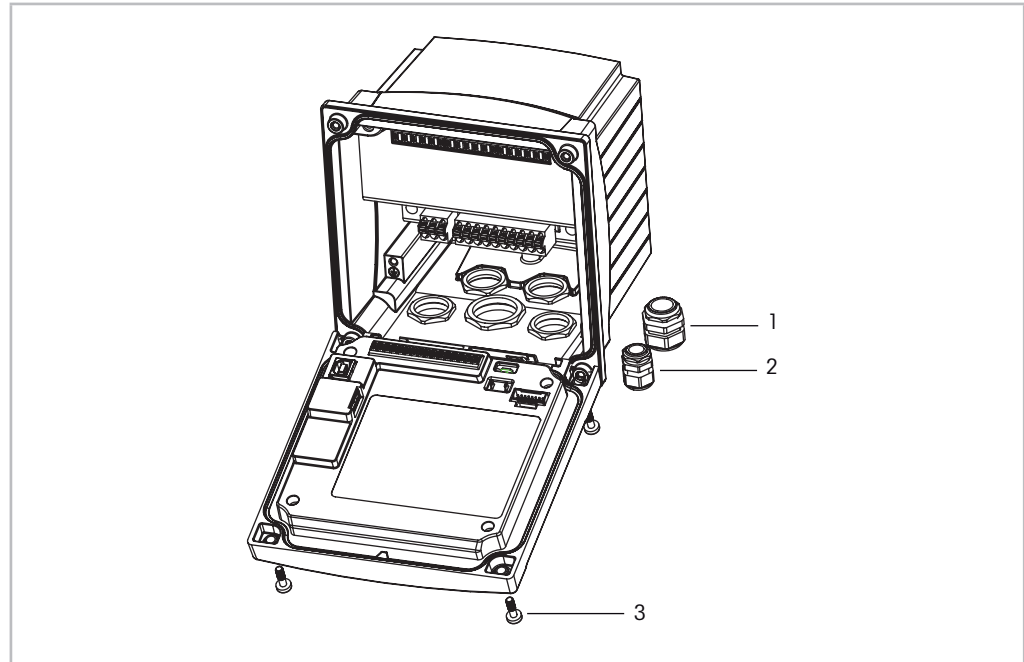


Fig. 5: Assembly

- 1 1 piece M25 x 1.5 cable gland
- 2 4 pieces M20 x 1.5 cable glands
- 3 4 pieces screws

General:

- Orient the transmitter so that the cable grips face downward.
- Wiring routed through the cable grips shall be suitable for use in wet locations.
- In order to provide IP66 enclosure ratings, all cable glands must be in place. Each cable gland must be filled using a cable.
- Tighten the screws of the front panel with a tightening torque of 2 Nm to 2.5 Nm.

4.2.3 ½ DIN – Panel Mounting

To insure a good seal, the panel or door must be flat and have a smooth finish. Textured or rough surfaces are not recommended and may limit the effectiveness of the gasket seal provided.

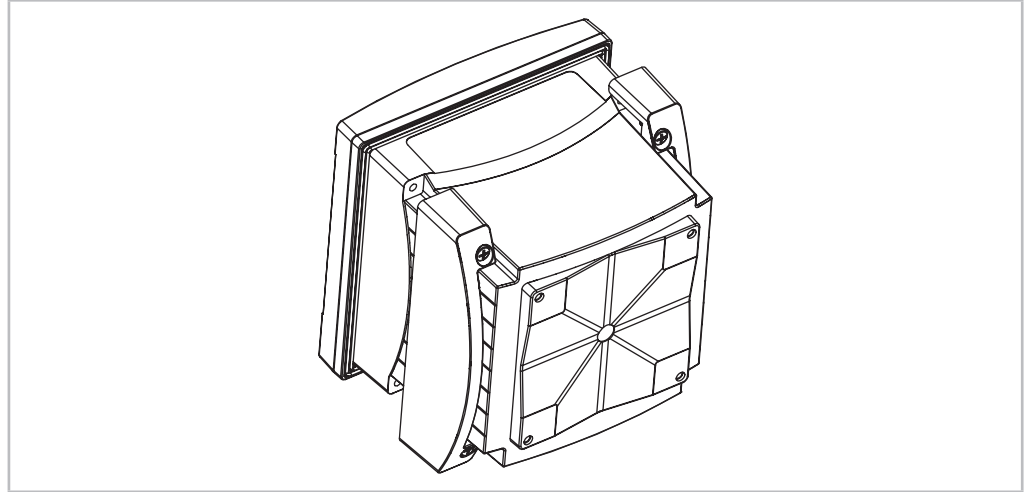


Fig. 6: Panel mounting

1. Make cutout in panel. For dimensions refer to chapter "4.2.1 Dimensions ½ DIN Version" on page 25.
 - Be sure surface surrounding cutout is clean, smooth and free of burrs.
2. Slide face gasket around transmitter from the back of the unit.
3. Place transmitter into cutout hole. Be sure there are no gaps between the transmitter and panel surface.
4. Place the two mounting brackets on either side of the transmitter as shown.
5. While holding transmitter firmly into the cutout hole, push the mounting brackets toward the backside of panel.
6. Once secure, use a screwdriver to tighten the brackets against the panel. In order to provide IP66 environmental enclosure rating, the two clamps provided shall be securely tightened to create an adequate seal between the panel enclosure and transmitter.
 - Face gasket will compress between transmitter and panel.

4.2.4 ½ DIN Version – Wall Mounting



DANGER! Mortal danger by electric shock or risk of electrical shock: The maximum screw-in depth of the mounting holes in the housing is 12 mm (0.47 inch). Do not exceed maximum screw-in depth.

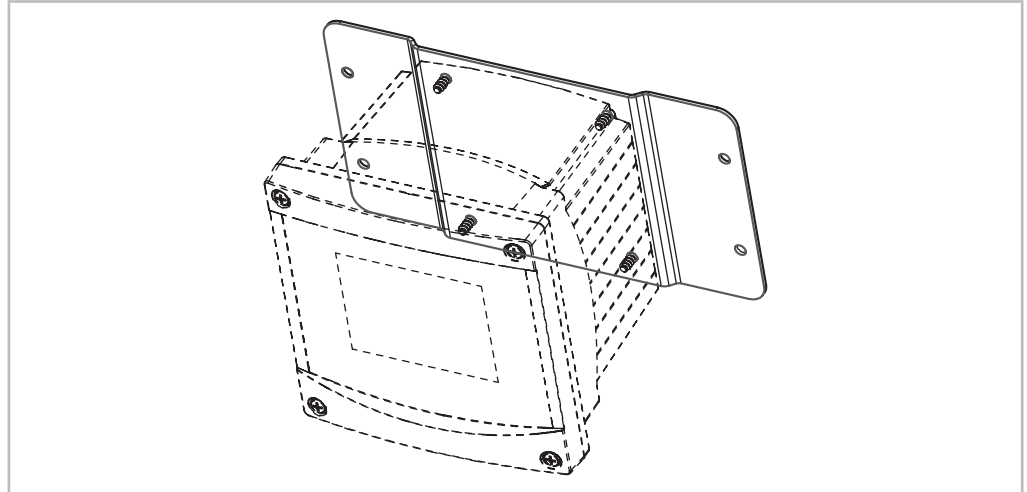


Fig. 7: Wall mounting with wall mounting kit

1. Mount wall mounting kit to the housing. Do not exceed maximum screw-in depth.
2. Mount wall mounting kit with the housing to the wall.
Attach to wall using appropriate mounting hardware for wall surface. Be sure it is level and securely fastened and the installation adheres to any and all clearance dimensions required for transmitter service and maintenance. Orient the transmitter so that the cable grips are facing downward.

4.2.5 ½ DIN Version – Pipe Mounting

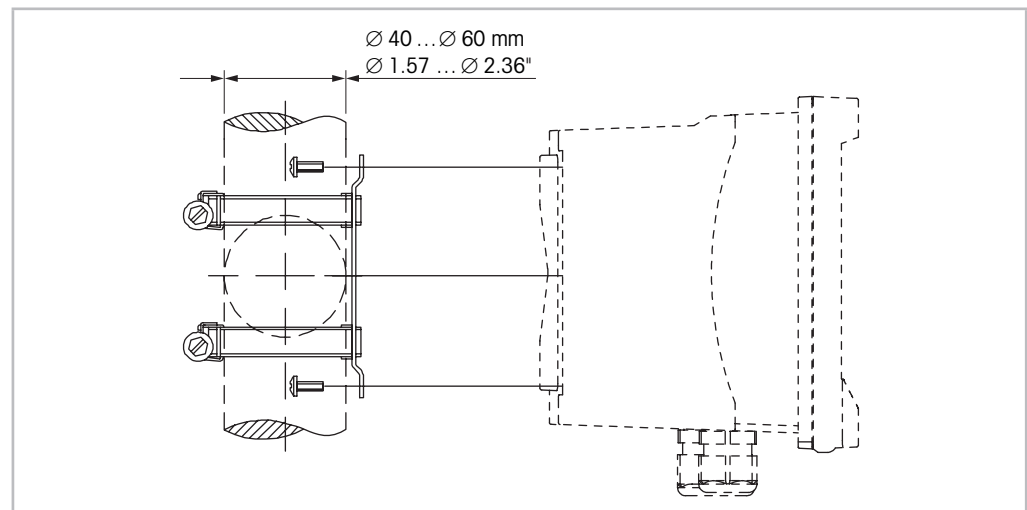


Fig. 8: Pipe mounting ½ DIN version

- Use only manufacturer-supplied components for pipe-mounting the M400 transmitter. See chapter "13 Ordering Information, Accessories and Spare Parts" on page 118 for ordering information.
- Tighten the fixing screws with a tightening torque of 2 to 3 Nm.

4.3 Electrical Connection



DANGER! Mortal danger by electric shock: Power off instrument during electrical connection.

1. Switch off supply voltage.
2. Connect mains supply to the terminals L, N, and \downarrow (Ground).
3. Connect sensor to terminal block TB3.
4. Connect analog output, analog input and digital input signals to terminal block TB2.
5. Connect relay output signals to terminal block TB1.
6. Connect the HART modem to AO1+/HART+ and AO1-/HART- for HART communication (communication load 230–500 ohm). Notice polarity.
7. Connect FOUNDATION fieldbus to FF+ and FF- for FF communication. Notice polarity.

For terminal definitions refer to the Operation Manual.

This is a 4-wire-product with an active 4–20 mA analog output. Do not supply power to the analog output terminals. M400 4-wire FF version has no analog outputs.



Warning! Do not disconnect the internal ground wire between the front and back modules.

Warning! Tightly secure an earth wire to the internal PE screw terminal:

\downarrow (Protective Conductor Terminal).

The cross-section of the PE wire must be above 18 AWG (0.8 mm).

4.4 Terminal Definition

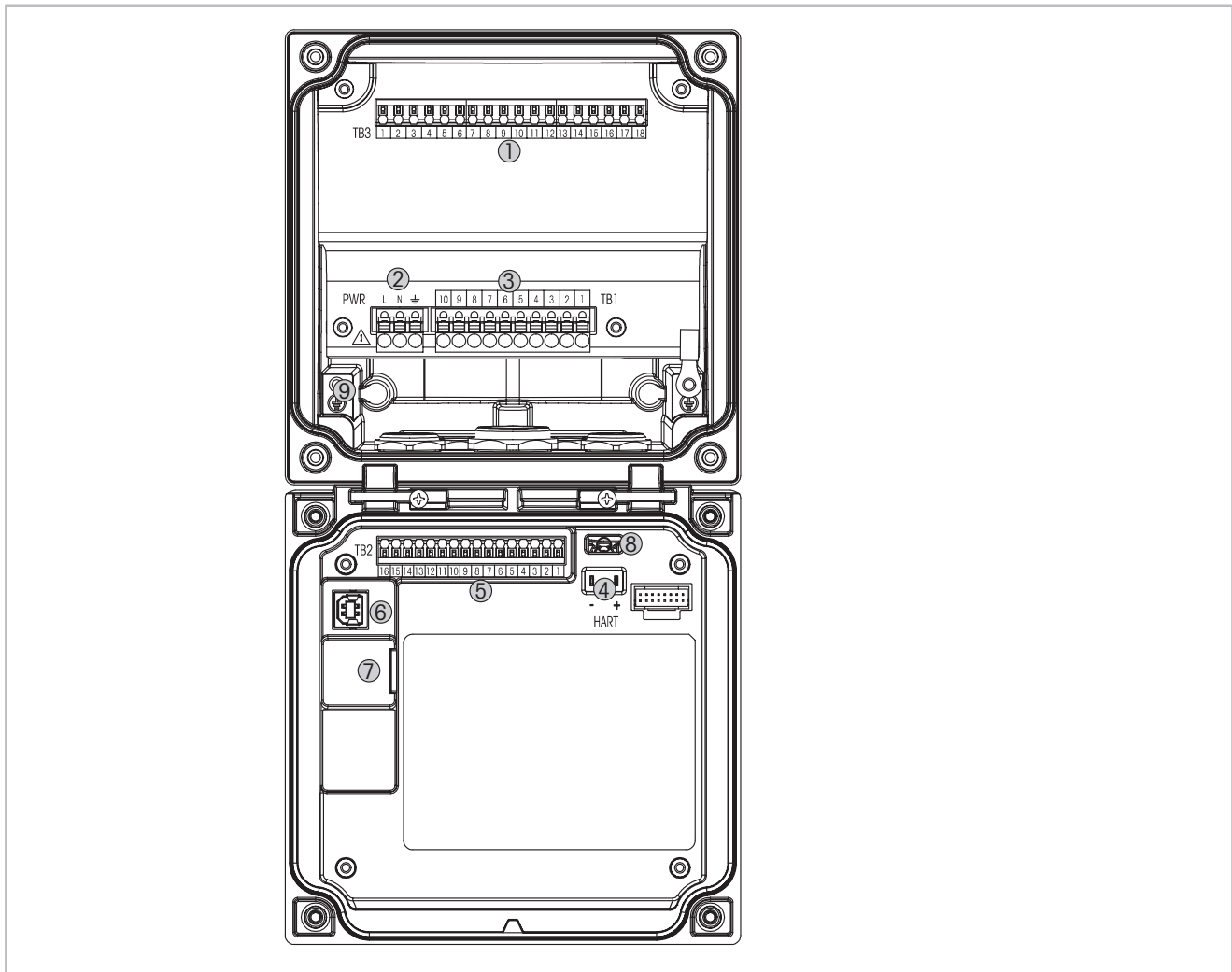


Fig. 9: M400 1/2 DIN versions

- 1 TB3 – Terminal block for sensor connection
- 2 Terminals for supply voltage
- 3 TB1 – Terminal block for relay outputs
- 4 HART, for M400 Type 1, Type 2, Type 3 only
- 5 TB2 – Terminal block for analog output and digital input signals
- 6 USB Device – Software update interface
- 7 USB Host – Printer connection, data logging, loading and saving configuration
- 8 **Warning! Do not disconnect the internal ground wire between the front and back modules.**
- 9 **Warning! Tightly secure an earth wire to the internal PE screw terminal:**
 ↓ (Protective Conductor Terminal).
 The cross-section of the PE wire must be above 18 AWG (0.8 mm).



4.4.1 TB1 Terminal definition – All transmitter versions

| Terminal | Description | Contact rating |
|----------|-------------|------------------------------|
| 1 | NO 1 | 250VAC or 30VDC, 3 A |
| 2 | COM | |
| 3 | NC 1 | |
| 4 | NO 2 | 250VAC or 30VDC, 3 A |
| 5 | COM | |
| 6 | NC 2 | |
| 7 | NO 3 | 250VAC or 30VDC, 0.5 A, 10 W |
| 8 | COM | |
| 9 | NO 4 | 250VAC or 30VDC, 0.5 A, 10 W |
| 10 | COM | |

4.4.2 TB2 Terminal definition

| Type 1, 2, 3 | | FF version | |
|--------------|---------------|------------|-------------|
| Terminal | Description | Terminal | Description |
| 1 | AO 1 +/HART + | 1 | FF+ |
| 2 | AO 1 –/HART – | 2 | FF- |
| 3 | AO 2 + | 3 | FF+ |
| 4 | AO 2 – | 4 | FF– |
| 5 | AO 3 + | 5 | Not used |
| 6 | AO 3 – | 6 | Not used |
| 7 | AO 4 + | 7 | Not used |
| 8 | AO 4 – | 8 | Not used |
| 9 | DI 1 + | 9 | DI1+ |
| 10 | DI 1 –/DI 2 – | 10 | DI1–/DI2– |
| 11 | DI 2 + | 11 | DI2+ |
| 12 | AI + | 12 | AI+ |
| 13 | AI – | 13 | AI– |
| 14 to 16 | Not used | 14 to 16 | Not used |

4.4.3 TB3 Terminal definition – Analog sensors

Conductivity 2-e/4-e

| Terminal | Function | Color |
|----------|---------------------------------|-------------|
| 1 | Cnd inner 1 ¹⁾ | White |
| 2 | Cnd outer 1 ¹⁾ | White/blue |
| 3 | Cnd outer 1 | – |
| 4 | Not used | – |
| 5 | Cnd outer 2 | – |
| 6 | Cnd inner 2 ²⁾ | Blue |
| 7 | Cnd outer 2 (GND) ²⁾ | Black |
| 8 | Not used | – |
| 9 | RTD ref/GND | Bare shield |
| 10 | RTD sense | Red |
| 11 | RTD | Green |
| 12 to 18 | Not used | – |

1) For third party Conductivity 2-e sensors a jumper between 1 and 2 may be required.

2) For third party Conductivity 2-e sensors a jumper between 6 and 7 may be required.

pH/ORP, Dissolved carbon dioxide

| Terminal | pH/dissolved carbon dioxide (InPro 5000) | | Redox (ORP) | |
|----------|--|---------------------|----------------------------|--------------|
| | Function | Color ¹⁾ | Function | Color |
| 1 | Glass | Transparent | Platinum | Transparent |
| 2 | Not used | – | – | – |
| 3 | Not used | – | – | – |
| 4 | Not used | – | – | – |
| 5 | Reference | Red | Reference | Red |
| 6 | Reference ²⁾ | – | Reference ²⁾ | – |
| 7 | Solution GND ²⁾ | Blue ³⁾ | Solution GND ²⁾ | – |
| 8 | Not used | – | – | – |
| 9 | RTD ref/GND | White | – | – |
| 10 | RTD sense | – | – | – |
| 11 | RTD | Green | – | – |
| 12 | Not used | – | – | – |
| 13 | Shield (GND) | Green/yellow | Shield (GND) | Green/yellow |
| 14 to 18 | Not used | – | – | – |

1) Grey wire not used.

2) Install jumper between 6 and 7 for ORP sensors and pH electrodes without SG.

3) Blue wire for electrode with SG.

Amperometric oxygen and Ozone – Analog sensors (continued)

| Terminal | Function | Oxygen | | | Ozone | |
|----------|---------------|-----------------|-----------------|-------------|-----------------------|-----------------|
| | | InPro 6800 | InPro 6900 | InPro 6950 | Hi Performance Oxygen | InPro 6510 |
| | | Color | Color | Color | Color | Color |
| 1 | Not used | – | – | – | – | – |
| 2 | Anode | Red | Red | Red | Red | Red |
| 3 | Anode | – ¹⁾ | – ¹⁾ | – | – ¹⁾ | – ¹⁾ |
| 4 | Reference | – ¹⁾ | – ¹⁾ | Blue | – ¹⁾ | – ¹⁾ |
| 5 | Not used | – | – | – | – | – |
| 6 | Not used | – | – | – | – | – |
| 7 | Guard | – | Grey | – | – | – |
| 8 | Cathode | Transparent | Transparent | Transparent | Grey | Grey |
| 9 | NTC ref (GND) | Green | Green | Green | Green | Green |
| 10 | Not used | – | – | – | – | – |
| 11 | NTC | White | White | White | White | White |
| 12 | Not used | – | – | – | – | – |
| 13 | Shield (GND) | Green/yellow | Green/yellow | | | Green/yellow |
| 14 to 18 | Not used | – | – | – | – | – |

1) Install jumper between 3 and 4

ISFET¹⁾

| Terminal | Function | Color |
|----------|-------------------------|-------------------|
| 1 | FET | Coax inner / pink |
| 2 | – | – |
| 3 | – | – |
| 4 | – | – |
| 5 | Reference | Yellow |
| 6 | Reference ²⁾ | – |
| 7 | Reference ²⁾ | – |
| 8 | – | – |
| 9 | RTD ref/GND | White |
| 10 | – | – |
| 11 | RTD | Grey |
| 12 | – | – |
| 13 | GND/Shield | Green |
| 14 | – | – |
| 15 | – | – |
| 16 | +5V | Brown |
| 17 to 18 | – | – |

1) When using InPro3300 sensor with special 5V cable [52300404]

2) Install jumper between 6 and 7 for ISFET sensors

4.4.4 TB3 Terminal definition – ISM sensors

pH/ORP, amperometric oxygen, dissolved ozone, conductivity 4-e, dissolved CO₂ low

| Terminal | Function | Color |
|----------|----------|--------------------------|
| 1 to 11 | Not used | – |
| 12 | 1-wire | Transparent (cable core) |
| 13 | GND | Red (shield) |
| 14 | RS485-B | – |
| 15 | RS485-A | – |
| 16 | 5V | – |
| 17 | GND 24V | – |
| 18 | 24V | – |

UniCond 2-e, UniCond 4-e

| Terminal | Function | Color |
|----------|----------|-------|
| 1 to 12 | Not used | – |
| 13 | GND | White |
| 14 | RS485-B | Black |
| 15 | RS485-A | Red |
| 16 | 5V | Blue |
| 17 to 18 | Not used | – |

Optical oxygen, dissolved CO₂ hi (InPro 5500 i), GPro 500 TDL

| Terminal | Function | Optical Oxygen, CO ₂ hi (InPro 5500 i) | | GPro 500 TDL |
|----------|-----------|---|-------------------------|--------------|
| | | VP8 cables wire color | 5-pin cables wire color | Color |
| 1 to 12 | Not used | – | – | – |
| 13 | GND | Green/Yellow | Green/Yellow | Brown |
| 14 | RS485-B | Brown | Blue | Yellow |
| 15 | RS485-A | Pink | White | Green |
| 16 | 5V | – | – | – |
| 17 | GND (24V) | Blue | Black | Blue |
| 18 | 24V | Grey | Brown | Red |

5 Placing transmitter in, or out, of service

5.1 Placing transmitter in service



After connecting the transmitter to power supply circuit, it will be active as soon as the circuit is powered.

5.2 Placing transmitter out of service

First disconnect the unit from the main power source, then disconnect all remaining electrical connections. Remove the unit from the panel. Use the installations instruction in this manual as reference for dis-assembling mounting hardware.

All transmitter settings stored in memory are non volatile.

6 Calibration

For the menu structure refer to chapter "3.9 Graphic Trend Measurement" on page 22.

PATH:  \ Cal



NOTE: During calibration, the outputs for the corresponding channel will default to be held at their current values until 20 seconds after the calibration menu is exited. A flashing H appears in the upper right corner of the display while outputs are held. Refer to chapter "7.3 Analog Outputs" on page 81 and chapter "7.4 Set Points" on page 82 to change the HOLD output status.

6.1 Sensor Calibration

PATH:  \ Cal \ Calibrate Sensor

6.1.1 Select the desired sensor calibration task

For analog sensors depending on sensor type, the following choices are available:

| Analog sensor | Calibration task |
|---------------|--|
| pH | pH, mV, Temperature, Edit, Verify |
| Conductivity | Conductivity, Resistivity, Temperature, Edit, Verify |
| Amp. Oxygen | Oxygen, Temperature, Edit, Verify |
| Ozone | Ozone, Temperature, Edit, Verify |

For ISM (digital) sensors depending on sensor type, the following choices are available:

| ISM sensor | Calibration task |
|----------------|-----------------------------------|
| pH | pH, ORP, Verify |
| Conductivity | Conductivity, Resistivity, Verify |
| Amp. Oxygen | Oxygen, Verify |
| Ozone | Ozone, Verify |
| Opt. Oxygen | Oxygen, Verify |
| Carbon Dioxide | Carbon dioxide, Verify |
| GPro 500 TDL | Cal, Verify |

6.1.2 Terminate Sensor Calibration

After every successful calibration different options are available. If "Adjust", "SaveCal" or "Calibrate" is chosen, the message "Calibration saved successfully! Reinstall sensor" is displayed. Press "Done" to return to the measuring mode.

| Option | Analog sensors | ISM (digital) sensors |
|--|---|--|
| Analog sensors: SaveCal ISM sensors: Adjust | Calibration values are stored in the transmitter and used for the measurement. Additionally, the calibration values are stored in the calibration data. | Calibration values are stored in the sensor and used for the measurement. Additionally, the calibration values are stored in the calibration history. |
| Calibrate | The function "Calibrate" is not applicable for analog sensors. | Calibration values are stored in the calibration history for documentation, but not be used for the measurement. The calibration values from the last valid adjustment are further used for the measurement. |
| Cancel | Calibration values are discarded. | Calibration values are discarded. |

6.2 Calibration of UniCond 2-e and UniCond 4-e Sensors (ISM Sensors only)

6.2.1 Conductivity Calibration of UniCond 2-e and UniCond 4-e Sensors

The M400 provides the ability to perform a one-point, two-point or process conductivity or resistivity calibration for 2-e-sensors and 4-e-sensors.




NOTE: When performing calibration on a conductivity sensor, results will vary depending on the method, calibration apparatus and/or quality of reference standards used to perform the calibration.



NOTE: For measuring tasks the temperature compensation for the application as defined through the parameter settings for conductivity will be considered and not the temperature compensation selected through the calibration procedure (see also chapter "7.1.4.1 Conductivity Settings" on page 71; PATH: \ CONFIG \ Meas \ Parameter Setting).

Enter the menu Calibrate Sensor (see chapter "6.1 Sensor Calibration" on page 36; PATH: \ Cal \ Calibrate Sensor) and choose the desired channel for calibration.

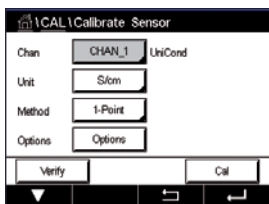

The following menus can be called up:

Unit: Choose between the units for conductivity (S/cm) and resistivity (Ω -cm).

Method: Select the desired calibration procedure. Available are 1-point, 2-point or process calibration.

Options: The desired compensation mode for the calibration process can be selected.

Choices are "None", "Standard", "Light 84", "Std 75 °C", "Linear 25°C", "Linear 20°C", "Glycol.5", "Glycol1", "Cation", "Alcohol" and "Ammonia".

None does not make any compensation of the measured conductivity value. The uncompensated value will be displayed and proceeded.

Standard compensation includes compensation for non-linear high purity effects as well as conventional neutral salt impurities and conforms to ASTM standards D1125 and D5391.

Light 84 compensation matches the high purity water research results of Dr. T.S. Light published in 1984. Use only if your institution has standardized on that work.

Std 75 °C compensation is the Standard compensation algorithm referenced to 75 °C. This compensation may be preferred when measuring Ultrapure Water at an elevated temperature (Resistivity of ultrapure water compensated to 75 °C is 2.4818 Mohm-cm.)

Linear 25 °C compensation adjusts the reading by a coefficient or factor expressed as %/°C (deviation from 25 °C). Use only if the solution has a well-characterized linear temperature coefficient. The factory default setting is 2.0%/°C. (2.4818 Mohm-cm.)

Linear 20 °C compensation adjusts the reading by a coefficient or factor expressed as %/°C (deviation from 20 °C). Use only if the solution has a well-characterized linear temperature coefficient. The factory default setting is 2.0%/°C.

Glycol.5 compensation matches the temperature characteristics of 50% ethylene glycol in water. Compensated measurements using this solution may go above 18 Mohm-cm.

Glycol1 compensation matches the temperature characteristics of 100% ethylene glycol. Compensated measurements may go well above 18 Mohm-cm.

Alcohol compensation provides for the temperature characteristics of a 75% solution of isopropyl alcohol in pure water. Compensated measurements using this solution may go above 18 Mohm-cm.




NOTE: If compensation mode "Linear 25 °C" or "Linear 20 °C" has been chosen, the coefficient for the adjustment of the reading can be modified. In this case an additional input field will be displayed.


The changes are valid until the calibration mode has been escaped. After the values defined in the configuration menu are valid again.

6.2.1.1 One-Point Calibration


Select calibration procedure 1-Point (see chapter “6.2.1 Conductivity Calibration of UniCond 2-e and UniCond 4-e Sensors” on page 37). With 2-e-sensors or 4-e-sensors a one-point calibration is always performed as a slope calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.

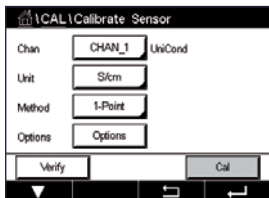

Press the button Cal for starting calibration.

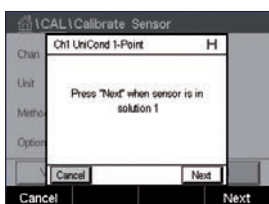

Place the electrode in the reference solution and press Next button.

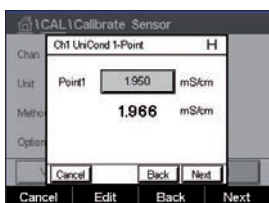


The second value displayed on the screen is the value being measured by the transmitter and sensor in units selected by the user.

Press the input field for **Point1** to enter the value for the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button and the transmitter will take over the value.




NOTE: To select another unit for the entered value on the keypad press the U button. To return again press the 0-9 button.

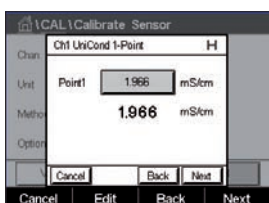

The screen shows the entered value for the reference solution (1st line) and the measured value of the M400 (2nd line).

Press the Next button to start the calculation of the calibration results.

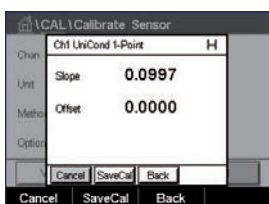

The display shows the value for the slope and the offset as the result of the calibration.


The calibration values are stored in the calibration history and taken over (press button SaveCal) or discarded (press button Cancel).

Use the Back button to go one step back in the calibration procedure.



6.2.1.2 Two-Point Calibration

Select calibration procedure 2-Point. With 4-e-sensors a two-point calibration is always performed as an offset and slope calibration. The following procedure shows the calibration with a 4-e-sensor.

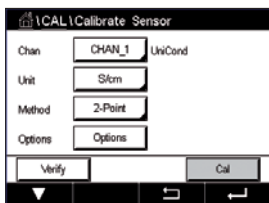


Press the button Cal for starting calibration.



Place the electrode in the first reference solution and press Next button.

CAUTION: Rinse sensors with a high-purity water solution between calibration points to prevent contamination of the reference solutions.

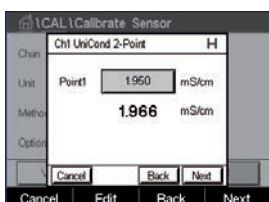


The second value displayed on the screen is the value being measured by the transmitter and sensor in the units selected by the user.

Press the input field for **Point1** to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button to accept the value.



NOTE: To select another unit for the entered value on the keypad press the U button. To return again press the 0-9 button.

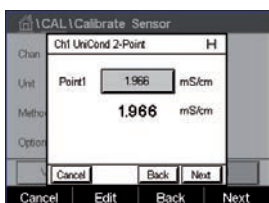


The screen shows the entered value for the first reference solution (1st line) and the measured value of the M400 (2nd line).

Press the Next button to go on with the calibration.



Place the electrode in the second reference solution and press Next button.

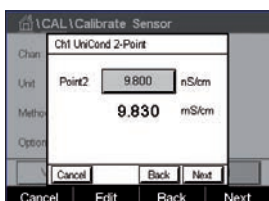


The second value displayed on the screen is the value being measured by the transmitter and sensor in the units selected by the user.

Press the input field for **Point2** to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button to accept the value.




NOTE: To select another unit for the entered value on the keypad press the U button. To return again press the 0-9 button.

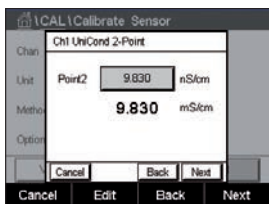

The screen shows the entered value for the second reference solution (1st line) and the measured value of the M400 (2nd line).


Press the Next button to start the calculation of the calibration results.

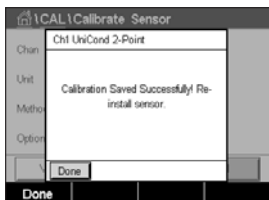

The display shows the value for the slope and the offset as the result of the calibration.


The calibration values are stored in the calibration history. To save (press button SaveCal) or to discard (press button Cancel).

Use the Back button to go one step back in the calibration procedure.

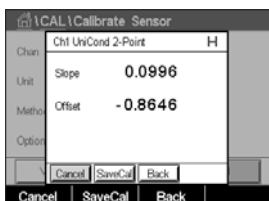


6.2.1.3 Process Calibration

Select calibration procedure Process (see chapter “6.2.1 Conductivity Calibration of UniCond 2-e and UniCond 4-e Sensors” on page 37). With 2-e-sensors or 4-e-sensors a process calibration is always performed as a slope calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.

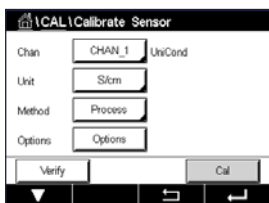


Press the button Cal for starting calibration.

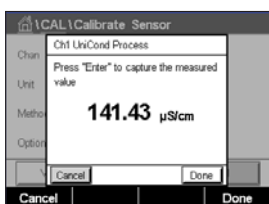


Take a sample and press the DONE button to store the current measuring value. To show the ongoing calibration process, P is blinking in the Start and Menu screen if the related channel is selected in the display.



After determining the conductivity value of the sample, press the calibration icon in the Menu Screen again.

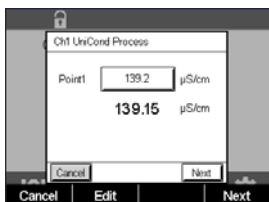

Press the input field for **Point1** and enter the conductivity value of the sample. Press the Next button to start the calculation of the calibration results.

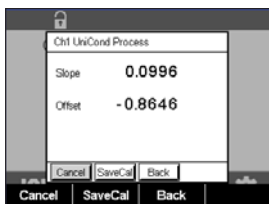

The display shows the value for the slope and the offset as the result of the calibration.

The calibration values are stored in the calibration history. To save (press button SaveCal) or to discard (press button Cancel).


Use the Back button to go one step back in the calibration procedure.



6.2.2 Temperature Calibration of UniCond 2-e Sensors and UniCond 4-e Sensors

The M400 provides the ability to perform a one-point or two-point calibration for the temperature sensor of the UniCond 2-e and UniCond 4-e.

Enter the menu Calibrate Sensor (see chapter "6.1 Sensor Calibration" on page 36; PATH: \ Cal \ Calibrate Sensor).

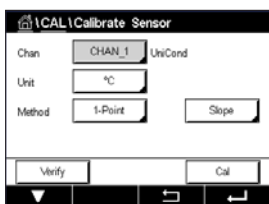

The following menus can be called up:


- Unit:** Choose between the units °C and °F.
- Method:** Select the desired calibration procedure. Available are 1-point and 2-point calibration.

6.2.2.1 One-Point Calibration


Select calibration procedure 1-Point. With 2-e-sensors or 4-e-sensors a one-point temperature calibration can be performed as a slope or offset calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.




Press the right input field for the parameter **Method**. Choose Slope or Offset calibration through pressing the corresponding field.

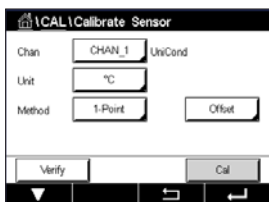

Press the button Cal for starting calibration.

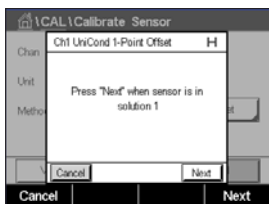


Place the electrode in the reference solution and press Next button.

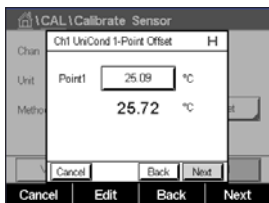


The second value displayed on the screen is the value being measured by the transmitter and sensor.

Press the input field or EDIT button for **Point1** to enter the value for the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button to accept the value.

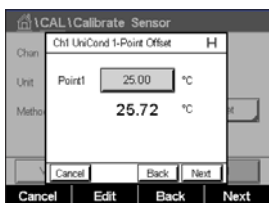

The screen shows the entered value for the reference solution (1st line) and the measured value of the M400 (2nd line).


Press the Next button to start the calculation of the calibration results.

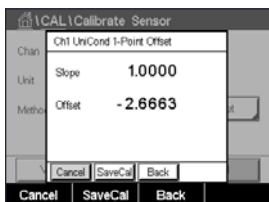

The display shows the value for the slope and the offset as the result of the calibration.


The calibration values are stored in the calibration history. To save (press button SaveCal) or to discard (press button Cancel).

Use the Back button to go one step back in the calibration procedure.

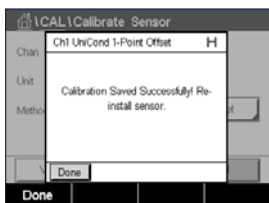

6.2.2.2 Two-Point Calibration


Select calibration procedure 2-Point (see chapter "6.2.2 Temperature Calibration of UniCond 2-e Sensors and UniCond 4-e Sensors" on page 42). With 2-e-sensors or 4-e-sensor a two-point calibration is always performed as an offset and slope calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.

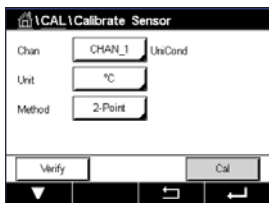

Press the button Cal for starting calibration.




Place the electrode in the first reference solution and press Next button.

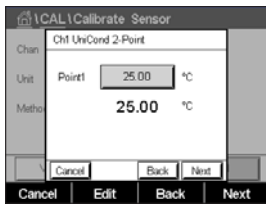


The second value displayed on the screen is the value being measured by the transmitter and sensor in the units selected by the user.

Press the input field for **Point1** to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button and to accept the value.

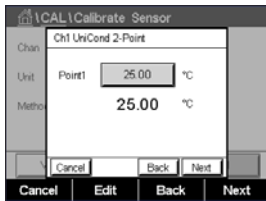


The screen shows the entered value for the first reference solution (1st line) and the measured value of the M400 (2nd line).

Press the Next button to go on with the calibration.

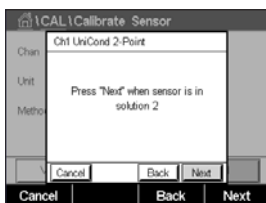


Place the electrode in the second reference solution and press Next button.

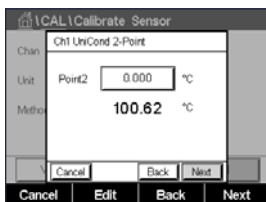


The second value displayed on the screen is the value being measured by the transmitter and sensor in the units selected by the user.

Press the input field for **Point2** to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button and to accept the value.

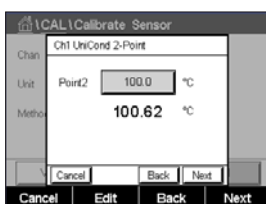

The screen shows the entered value for the second reference solution (1st line) and the measured value of the M400 (2nd line).


Press the Next button to start the calculation of the calibration results.

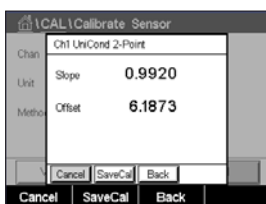

The display shows the value for the slope and the offset as the result of the calibration.

The calibration values are stored in the calibration history. To save (press button SaveCal) or to discard (press button Cancel).

Use the Back button to go one step back in the calibration procedure.

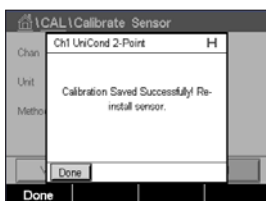

6.3 Calibration of Cond2e Sensors or Cond4e Sensors

PATH:  \ Cal \ Calibrate Sensor


The M400 provides the ability to perform a one-point, two-point or process conductivity or resistivity calibration for 2-e-sensors and 4-e-sensors.



NOTE: When performing calibration on a conductivity sensor, results will vary depending on the method, calibration apparatus and/or quality of reference standards used to perform the calibration.



NOTE: For measuring tasks the temperature compensation for the application as defined through the parameter settings for conductivity will be considered and not the temperature compensation selected through the calibration procedure (see also chapter "7.1.4.1 Conductivity Settings" on page 71).

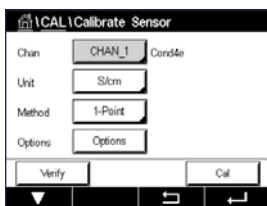

The following menus can be called up:

Unit: Between the units for conductivity and resistivity can be chosen.

Method: Select the desired calibration procedure, 1-point, 2-point or process calibration.

Options: Select the desired temperature compensation mode for the calibration process.




NOTE: If compensation mode "Linear 25 °C" or "Linear 20 °C" has been chosen, the coefficient for the adjustment of the reading can be modified.

The changes are valid until the calibration mode has been exited. After the values defined in the configuration menu are valid again.

6.3.1 One-Point Calibration

With 2-e-sensors or 4-e-sensors a one-point calibration is always performed as a slope calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.

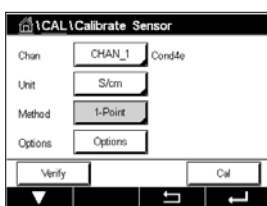

Press the button Cal for starting calibration.

Place the electrode in the reference solution and press Next button.

Enter the value for the calibration point (**Point1**).


Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See chapter "6.1.2 Terminate Sensor Calibration" on page 37.

6.3.2 Two-Point Calibration

With 2-e-sensors or 4-e-sensors a two-point calibration is always performed as an offset and slope calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.

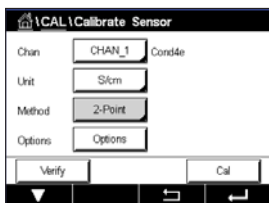

Press the button Cal for starting calibration.

Place the electrode in the first reference solution and press Next button.

CAUTION: Rinse sensors with a high-purity water solution between calibration points to prevent contamination of the reference solutions.

Enter the value for the first calibration point (**Point1**).

Press the Next button to go on with the calibration.

Place the electrode in the second reference solution and press Next button.

Enter the value for the second calibration point (**Point2**).


Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

6.3.3 Process Calibration

With 2-e-sensors or 4-e-sensors a process calibration is always performed as a slope calibration. The following procedure shows the calibration with a 2-e-sensor. The calibration with a 4-e-sensor works respectively.

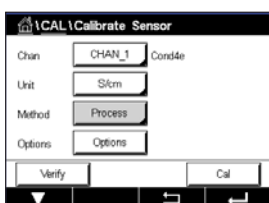

Press the button Cal for starting calibration.

Take a sample and press the \leftarrow button to store the current measuring value. To show the ongoing calibration process, P is blinking in the Start and Menu screen if the related channel is selected in the display.

After determining the conductivity value of the sample, press the calibration icon in the Menu Screen again.

Enter the conductivity value of the sample. Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.


6.4 pH Calibration

PATH:  \ Cal \ Calibrate Sensor

For pH sensor include ISFET, the M400 Transmitter features one-point, two-point or process calibration with 9 preset buffer sets or manual buffer entry. Buffer values refer to 25 °C. To calibrate the instrument with automatic buffer recognition, you need a standard pH buffer solution that matches one of these values. Please select the correct buffer table before using automatic calibration (see chapter "16 Buffer tables" on page 125). The stability of the sensor signal during calibration can be checked by the user or automatically by the transmitter (see chapter "7.1.4.2 pH Settings" on page 72).



NOTE: For dual membrane pH electrodes (pH/pNa) only buffer Na+ 3.9M (see chapter "16.2.1 Mettler-pH/pNa buffers (Na+ 3.9M)" on page 130) is available.

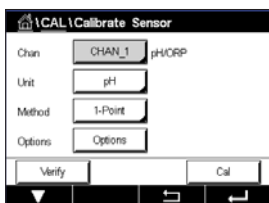

The following menus can be called up:

Unit: Select pH.

Method: Select the desired calibration procedure, 1-point, 2-point or process calibration.

Options: The buffer used for the calibration and the required stability of the sensor signal during the calibration can be selected (see also chapter "7.1.4.2 pH Settings" on page 72). The changes are valid until the calibration mode has been escaped. After the values defined in the configuration menu are valid again.



NOTE: When do ISFET one-point calibration, select unit as "mV".

6.4.1 One-Point Calibration

With pH sensors a one-point calibration is always performed as an offset calibration.

Press the button Cal for starting calibration.

Place the electrode in the buffer solution and press the Next button.

The display shows the buffer the transmitter has recognized **Point 1** and the measured value.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.



NOTE: When measuring with an ISFET sensor, the nominal zero point must be adjusted each time a new sensor is connected or after CIP. Immerse sensor in a zero point buffer (pH7). Make a mV calibration and enter for point 1 the value 00.00 mV.

The transmitter shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

6.4.2 Two-Point Calibration

With pH sensors a two-point calibration is always performed as calibration of slope and offset.

Press the Cal button to start calibration.

Place the electrode in buffer solution 1 and press Next button.

The display shows the buffer the transmitter has recognized **Point 1** and the measured value.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.

NOTE: If option Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.

The transmitter prompts you to place the electrode in the second buffer solution.

Press the Next button to proceed with the calibration.

The display shows the buffer the transmitter has recognized **Point 2** and the measured value.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.

NOTE: If option Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.

The transmitter shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See chapter "6.1.2 Terminate Sensor Calibration" on page 37.

6.4.3 Process Calibration

With pH sensors a process calibration is always performed as an offset calibration.

Press the Cal button to start calibration.



Take a sample and press the \leftarrow button to store the current measuring value. To show the ongoing calibration process, P is blinking in the Start and Menu Screen if the related channel is selected in the display.

After determining the pH value of the sample, press the calibration icon in the Menu Screen again.

Enter the pH value of the sample. Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See chapter "6.1.2 Terminate Sensor Calibration" on page 37.

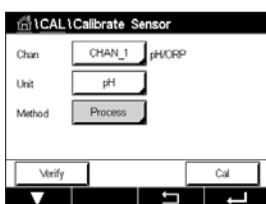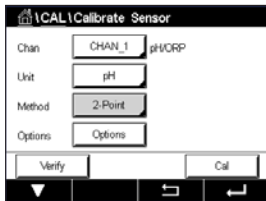


6.5 ORP Calibration of pH Sensors

PATH:  \ Cal \ Calibrate Sensor

For pH sensors with solution ground based on ISM technology the M400 Transmitter gives the option to make, in addition to the pH calibration, an ORP calibration.



NOTE: In case of choosing ORP calibration the parameters defined for pH (see chapter "7.1.4.2 pH Settings" on page 72) will not be considered. For pH sensors, the M400 Transmitter features one-point calibration for ORP.

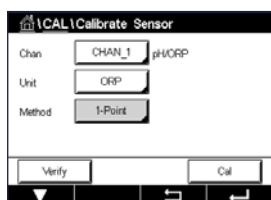

The following menus can be called up:

Unit: Select ORP through pressing the corresponding field.

Method: 1-Point calibration is displayed.

Press the button Cal for starting calibration.

Enter the value for calibration point 1 (**Point1**).

Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.


6.6 Calibration of Amperometric Oxygen Sensors

PATH:  \ Cal \ Calibrate Sensor

The M400 provides the ability to perform a one-point or process calibration for amperometric oxygen sensors.



NOTE: Before air calibration, for highest accuracy, enter the barometric pressure and relative humidity, as described in chapter "7.1.4.3 Settings for Oxygen Measurement Based on Amperometric Sensors" on page 74.

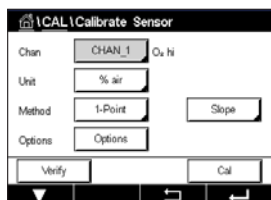

The following menus can be called up:


Unit: Between several units for Dissolved Oxygen can be chosen.

Method: Select the desired calibration procedure, 1-point or process calibration.

Options: In case the method 1-point has been chosen the calibration pressure, relative humidity and - for slope calibration - the stability mode for the sensor signal during the calibration can be selected. For the method Process the values for the process pressure, calibration pressure and the parameter ProcCalPress can be modified. See also chapter "7.1.4.3 Settings for Oxygen Measurement Based on Amperometric Sensors" on page 74. The changes are valid until the calibration mode has been escaped. After the values defined in the configuration menu are valid again.

6.6.1 One-Point Calibration

A one-point calibration of oxygen sensors is always either a one-point slope (i.e. with air) or a zero (offset) calibration. A one-point slope calibration is done in air and a one-point offset calibration is done at 0 ppb oxygen. A one-point zero dissolved oxygen calibration is available but not normally recommended since zero oxygen is very hard to achieve. A zero-point calibration is only recommended if high accuracy at low oxygen level (below 5% air) is needed.

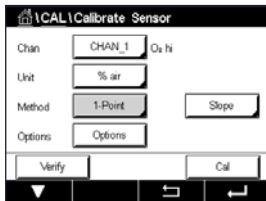

Choose Slope or Offset calibration through pressing the corresponding field.

Press the button Cal for starting calibration.



NOTE: If the polarization voltage for the measuring mode and calibration mode is different, the transmitter will wait 120 seconds before starting the calibration. In this case the transmitter will also go after the calibration for 120 seconds to the HOLD Mode, before returning to the measuring mode again.

Place the sensor in air or the calibration gas and press Next button

Enter the value for the calibration point (**Point1**).

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.




NOTE: For an offset calibration the Auto mode is not available. If Auto mode has been chosen and afterwards slope calibration has been changed to offset calibration, the transmitter will perform the calibration in Manual mode.

The transmitter shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

6.6.2 Process Calibration

A process calibration of oxygen sensors is always either a slope or an offset calibration.

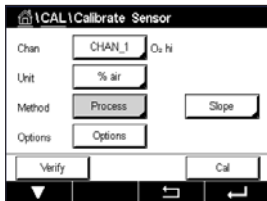

Choose Slope or Offset calibration through pressing the corresponding field.

Press the Cal button to start calibration.

Take a sample and press the \leftarrow button to store the current measuring value. To show the ongoing calibration process, P is blinking in the Start and Menu screen if the related channel is selected in the display.

After determining the oxygen value of the sample, press the calibration icon in the Menu Screen again.

Enter the oxygen value of the sample. Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.


6.7 Calibration of Optical Oxygen Sensors (ISM Sensors only)

PATH: \ Cal \ Calibrate Sensor

Oxygen calibration for optical sensors can be performed as a two-point, process or, depending on the sensor model connected to the transmitter, also as a one-point calibration.



NOTE: Before air calibration, for highest accuracy, enter the barometric pressure and relative humidity, as described in chapter "7.1.4.4 Settings for Oxygen Measurement Based on Optical Sensors" on page 75.

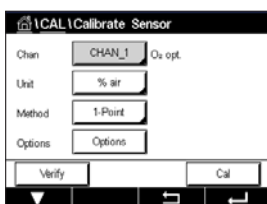


The following menus can be called up:

- Unit:** Between several units can be chosen. The units are displayed during the calibration.
- Method:** Select the desired calibration procedure, 1-point, 2-point or process calibration.
- Options:** In case the method 1-point has been chosen the calibration pressure, relative humidity and the stability mode for the sensor signal during the calibration can be selected. For the method Process the values for the process pressure, calibration pressure, the parameter ProcCalPress and the mode of the process calibration can be modified. See also chapter "7.1.4.4 Settings for Oxygen Measurement Based on Optical Sensors" on page 75. The changes are valid until the calibration mode has been escaped. After the values defined in the configuration menu are valid again.

6.7.1 One-Point Calibration

Typically a one-point calibration is done in air. Nevertheless other calibration gases and solutions are possible.

The calibration of an optical sensor is always a calibration of the phase of the fluorescence signal towards the internal reference. During a one-point calibration the phase in this point is measured and extrapolated over the measuring range.

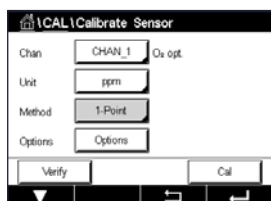

Press the button Cal for starting calibration.

Place the sensor in air or the calibration gas and press Next button

Enter the value for the calibration point (**Point1**).

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.


The transmitter shows the value for the phase of the sensor at 100% air (P100) and at 0% air (P0) as the result of the calibration.

Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed. In either case you will see the message "Please re-install sensor".

6.7.2 Two-Point Calibration

The calibration of an optical sensor is always a calibration of the phase of the fluorescence signal towards the internal reference. A two-point calibration is a combination of first a calibration in air (100%) where a new phase P100 is measured and then a calibration in nitrogen (0%) where a new phase P0 is measured. This calibration routine gives the most accurate calibration curve over the whole measuring range.

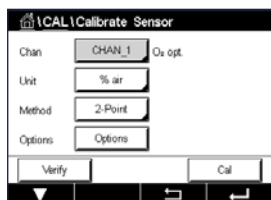

Press the Cal button to start calibration.

Place the sensor in air or the calibration gas and press Next button

Enter the value for the first calibration point (**Point1**).

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.

The transmitter prompts you to change the gas.

Press the Next button to proceed with the calibration.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.




NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.

The transmitter shows the value for the phase of the sensor at 100% air (P100) and at 0% air (P0) as the result of the calibration.

Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed. In either case you will see the message "Please re-install sensor".

6.7.3 Process Calibration

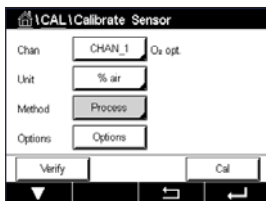

Press the Cal button to start calibration.

Take a sample and press the \leftarrow button to store the current measuring value. To show the ongoing calibration process, P is blinking in the start and Menu Screen if the related channel is selected in the display.

After determining the oxygen value of the sample, press the calibration icon in the Menu Screen.

Enter oxygen value of the sample. Press the Next button to start the calculation of the calibration results.

The display shows now the values for the phase of the sensor at 100% air (P100) and at 0% (P0) air.

Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.




NOTE: If for process calibration Scaling has been chosen (see chapter "7.1.4.4 Settings for Oxygen Measurement Based on Optical Sensors" on page 75) the calibration values are not stored in the calibration history.

If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed.

6.8 Calibration of Dissolved Carbon Dioxide Sensors (ISM Sensors only)

For dissolved carbon dioxide (CO₂) sensors, the M400 Transmitter features one-point, two-point or process calibration. For the one-point or two-point calibration the solution with pH = 7.00 and/or pH = 9.21 of the Mettler – 9 standard buffer can be used (see also chapter "7.1.4.5 Dissolved Carbon Dioxide Settings" on page 76) or the buffer value can be entered manually.

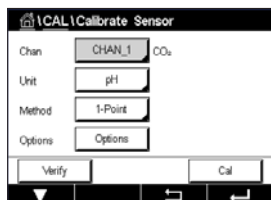

The following menus can be called up:


Unit: Between several units for partial pressure, and dissolved carbon dioxide can be selected.

Method: Select the desired calibration procedure, 1-point or process calibration.

Options: The buffer used for the calibration and the required stability of the sensor signal during the calibration can be selected (see also chapter "7.1.4.5 Dissolved Carbon Dioxide Settings" on page 76). The changes are valid until the calibration mode has been escaped. After the values defined in the configuration menu are valid again.

6.8.1 One-Point Calibration

With CO₂ sensors a one-point calibration is always performed as an offset calibration.

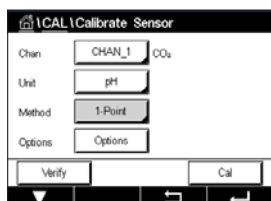

Press the button Cal for starting calibration.

Place the electrode in the buffer solution and press the Next button.

The display shows the buffer the transmitter has recognized **Point 1** and the measured value.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.


The transmitter shows the value for the slope and the offset as the result of the calibration.

Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed. In either case you will get the message "Please re-install sensor".

6.8.2 Two-Point Calibration

With CO₂ sensors a two-point calibration is always performed as calibration of slope and offset.

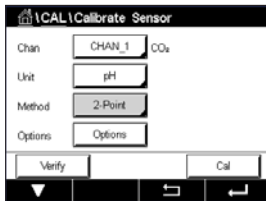

Press the Cal button to start calibration.

Place the electrode in buffer solution 1 and press Next button.

The display shows the buffer the transmitter has recognized **Point 1** and the measured value.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.

The transmitter prompts you to place the electrode in the second buffer solution.

Press the Next button to proceed with the calibration.

The display shows the buffer the transmitter has recognized **Point 2** and the measured value.

The M400 checks the stability of the measuring signal and proceeds as soon as the signal is sufficiently stable.



NOTE: If **option** Stability is set to Manual press 'Next' after the measuring signal is stable enough to go on with the calibration.


The transmitter shows the value for the slope and the offset as the result of the calibration.

Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed. In either case you will see the message "Please re-install sensor".

6.8.3 Process Calibration

With CO₂ sensors a process calibration is always performed as an offset calibration.

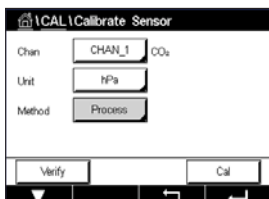

Press the Cal button to start calibration.

Take a sample and press the \leftarrow button to store the current measuring value. To show the ongoing calibration process, P is blinking in the Start and Menu Screen if the related channel is selected in the display.

After determining the corresponding value of the sample, press the calibration icon in the Menu Screen again.

Enter the value of the sample. Press the Next button to start the calculation of the calibration results.

The display shows the value for the slope and the offset as the result of the calibration.

Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration. If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed.

6.9 Calibration of Thermal Conductivity CO₂ (CO₂ high) Sensors (ISM Sensors only)


PATH: \ Cal \ Calibrate Sensor

The M400 provides the ability to perform a one-point calibration using a reference gas (CO₂) with a known carbon dioxide partial pressure value. It offers also to perform a process calibration based on an analyzed process sample.



NOTE: The sensor is designed to measure CO₂ partial pressure or concentration values accurately in liquid phase only! In gas phase the sensor will only show correct CO₂ gas partial pressure values in the 1-point calibration menu.

The following menus can be called up:

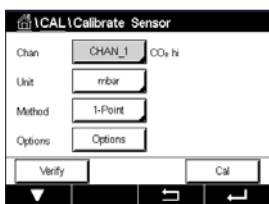

Unit: Between the units of CO₂ pressure or concentration can be chosen.

Method/Options: Select the desired calibration procedure, 1-point or process calibration and stability option (manual/auto).

In case the method 1-point has been chosen only the calibration pressure and the option stability mode for the sensor signal during the calibration can be selected (Sensor expects to be in a calibration gas).

For the method Process only concentration values can be chosen as pressure or concentration values (Sensor expects to be in liquids).




NOTE: With reference Gas (CO₂) use 1-point calibration. With liquids use process calibration. When changing MembraCap always first perform a 1-point gas calibration. The changes are valid until the calibration mode has been exited. After the values defined in the Configuration menu are valid again.

6.9.1 One-Point Calibration




With the thermal conductivity sensor a one-point calibration is always performed as a slope calibration. Press the button Cal for starting calibration.

Expose the TC-Sensor to a reference gas of a known CO₂ concentration and press Next button. Enter the value for the calibration point (Point1) in mbar or hPa.

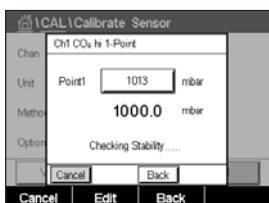


Press the Next button to start the calculation of the calibration results.

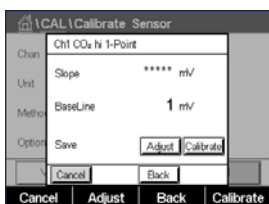

The display shows the value for the slope and the baseline as the result of the calibration. Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.


If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed.

6.9.2 Process Calibration

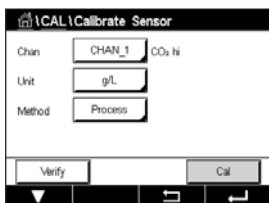

With the thermal conductivity sensor a process calibration is always performed as a slope calibration.


Choose process calibration and desired unit in the calibration menu. Press the button Cal for starting calibration.

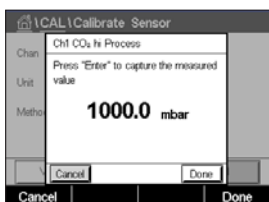

Take a sample and press the \leftarrow button to store the current measuring value. To show the ongoing calibration process, P is blinking in the Start and Menu screen if the related channel is selected on the display.

After determining the CO₂ value of the sample, press the calibration icon in the Menu Screen again. Enter the CO₂ value of the sample.

Press the Next button to start the calculation of the calibration results.

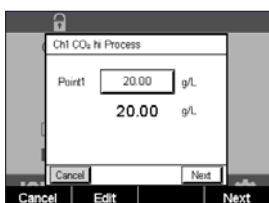

The display shows the value for the slope and the baseline as the result of the calibration.


Press the Adjust button to perform the calibration and store the calculated values in the sensor. Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Calibration Saved Successfully!" is displayed.

6.10 Calibration of O₃ Sensors

The M400 provides the ability to perform a 1-Point or process calibration for O₃ sensors. Dissolved Ozone must be performed quickly because O₃ decays rapidly into oxygen, especially at warm temperatures.

Enter the menu Calibrate Sensor (see chapter "6.1 Sensor Calibration" on page 36; PATH: \ Cal \ Calibrate Sensor) and choose the desired channel for calibration.

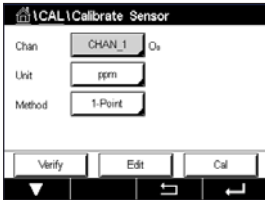


The following menus can be called up:

Unit: Several units for dissolved O₃ can be chosen.

Method: Select the desired calibration procedure, 1-Point or process calibration.

6.10.1 One-Point Calibration


Select the 1-Point calibration method. A one-point calibration of O₃ sensors is always a zero (offset) calibration

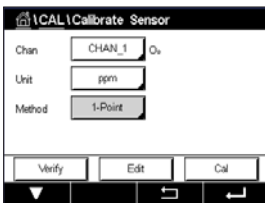

Press the button Cal for starting calibration.




Place the sensor in the calibration gas, such as air, and press the Next button.

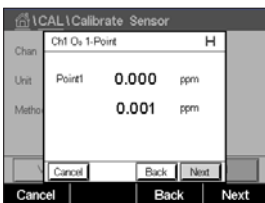

The second value displayed on the screen is the value being measured by the transmitter and sensor in the units selected by the user.

Press the input field for **Point1** to enter the value for the calibration point. The M400 displays a keypad for modifying the value. Press the button to accept the value.


When the measuring signal is stable, press Next to continue with the calibration

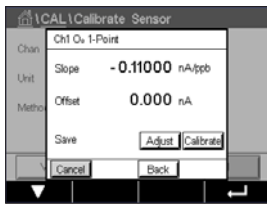

The display shows the value for the slope and the offset as result of the calibration.


For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

Use the Back button to go one step back in the calibration procedure

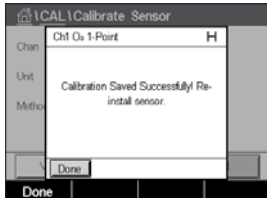


6.10.2 Process Calibration


Select the Process calibration method. A Process calibration of O₃ sensors can be performed as a slope or offset calibration.

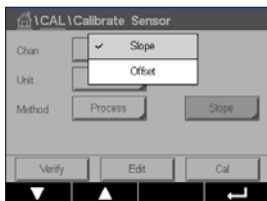


Select the desired calibration **Method**.

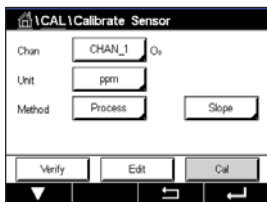

Press Cal to start the calibration.

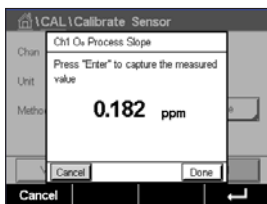

Take a sample and press the \leftarrow button to store the current measuring value. "P" will blink in the measurement screen indicating a Process calibration is active.

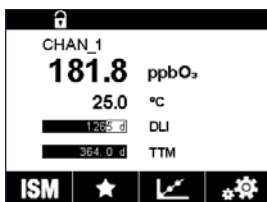

After determining the O₃ value of the sample, press the calibration icon to complete the Process calibration.



Press the input field for **Point1** and enter the O₃ value of the sample. Press the \leftarrow button to accept the value.

Press the Next button to start the calculation of the calibration results.



The display shows the value for the slope and the offset as the result of the calibration.

For ISM (digital) sensors select "Adjust", "Calibrate" or "Cancel" to finish calibration. For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

Use the Back button to go one step back in the calibration procedure.




6.11 Calibration of a tunable diode laser (TDL) analyzer

PATH: \ Cal \ Calibrate Sensor

Calibration for a TDL sensor is performed as a one-point or process calibration.

The following menus can be called up:

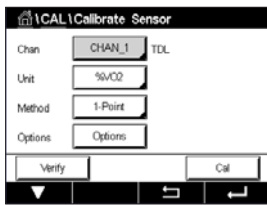

Unit: One of several units can be chosen. The units are displayed during the calibration.


Method: Select the desired calibration procedure, 1-point or process calibration.

Options: If the 1-point method has been chosen the calibration pressure, temperature and the path length for the sensor signal during the calibration can be edited.


See also chapter 7.1.5.4 "Settings for Oxygen Measurement Based on Optical Sensors". The changes are valid until the calibration mode has been exited. After, the values defined in the configuration menu are valid again.


6.11.1 One-point calibration for TDL gas sensors

A one-point calibration of gas sensors is always a slope (i.e. with air) calibration. A one point slope calibration is done in air or any other calibration gas with defined gas concentration.

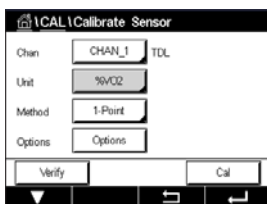


In the case of a dual gas (for example CO and CO₂) TDL selects the gas to be calibrated.

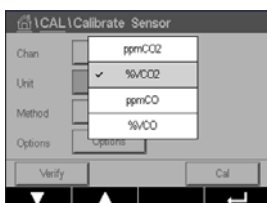

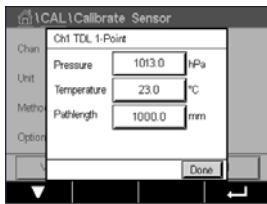

Adjust calibration pressure and temperature, which are applied during calibration.

Adjust the optical path length for your individual system.

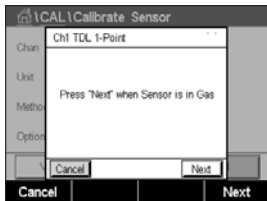

Press the button Cal for starting calibration

Place the sensor in the calibration gas (e.g. air). Press NEXT.

Enter the value for the calibration point then press Next to start the calculation.

The M400 checks the deviation of the measuring signal and proceeds as soon as the signal is sufficiently stable.


The display shows the value of the sensor as the result of calibration.

Press the adjust button to perform the calibration and store the calculated values in the sensor.


Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Adjustment Saved Successfully!" or "Calibration Saved Successfully!" is displayed. In either case you will see the message "Please re-install sensor".

6.11.2 Process calibration for TDL gas sensors

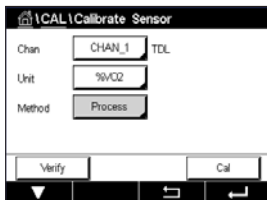


A process calibration of gas sensors is always a slope calibration.

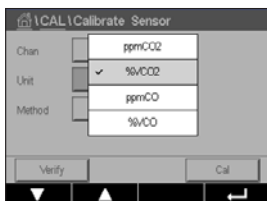


In the case of a dual gas (for example CO and CO₂) TDL selects the gas to be calibrated.

Press the Cal button to start calibration.

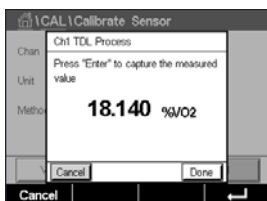


Take a sample and press the [ENTER] key again to store the current measuring value. To show the ongoing calibration process, P is blinking in the start and menu screen.

After determining the concentration value of the sample press the calibration icon in the menu screen to proceed with the calibration.

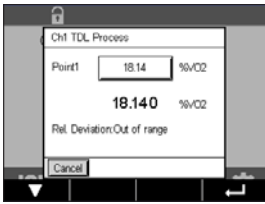

Enter the value for the calibration point then press Next to start the calculation.

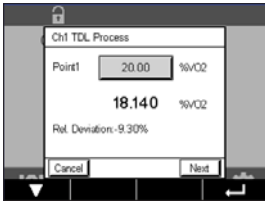

The M400 checks the deviation of the measuring signal and proceeds as soon as the signal is sufficiently stable.

The display shows the value of the sensor as the result of calibration.


Press the adjust button to perform the calibration and store the calculated values in the sensor.

Press the Calibrate button to store the calculated values in the sensor. Calibration is not performed. Press the Cancel button to terminate the calibration.

If "Adjust" or "Calibrate" are chosen, the message "Adjustment Saved Successfully!" or "Calibration Saved Successfully!" is displayed. In either case you will see the message "Please re-install sensor."

6.12 Sensor Verification

Enter the menu Calibrate Sensor (see chapter "6.1 Sensor Calibration" on page 36; PATH: \ Cal \ Calibrate Sensor) and choose the desired channel for verification

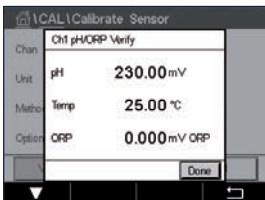


Press the Verify button to start verification.

The measured signal of the primary and the secondary measurement in basic (mostly electrical) units are shown. The meter calibration factors are used when calculating these values.

Press the button and the transmitter returns to the calibration menu.

6.13 UniCond 2-e Electronics Calibration (ISM Sensor only)

The M400 provides the ability to calibrate or verify the electronic circuits of UniCond 2-e conductivity sensors. UniCond 2-e sensors have 3 resistance range circuits that require individual calibration. These measuring circuits are calibrated using the THORNTON ISM Conductivity Sensor Calibration Module part number 58 082 305 and supplied Y-connector. Before calibration, remove the sensor from the process, rinse with deionized water and allow to completely dry. Power the transmitter and sensor at least 10 minutes prior to calibration to assure stable operating temperature of the circuitry.

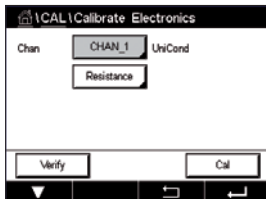

Press the Cal button.

Enter menu Calibrate Electronics.

Press the Chan_x button and select the desired channel for calibration.


Choose **Verify** or **Cal**.

Reference THORNTON ISM Conductivity Sensor Calibration Module (part number 58 082 305) for detailed calibration and verification instructions.

6.14 Meter Calibration (Analog Sensors only)

Although it is not normally necessary to perform meter re-calibration unless extreme conditions cause an out of spec operation shown by Calibration Verification, periodic verification/re-calibration may be necessary to meet Q.A. requirements. The frequency calibration requires a two-point calibration. It is recommended that point one be at the low end of the frequency range and point two at the high end.

Press the Cal button.

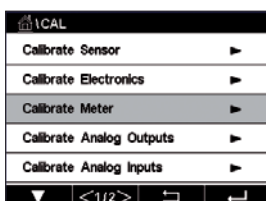


Enter menu Calibrate Meter.

6.14.1 Resistance (Analog Sensors only)

The meter is equipped with five (5) internal ranges of measurement. Each resistance range and temperature is calibrated separately, with each resistance range consisting of a two-point calibration.


Below is a table showing the resistance values for all calibration ranges.

| Range | Point 1 | Point 2 | Point 4 |
|---------------|-------------|-------------|----------|
| Resistivity 1 | 1.0 Mohms | 10.0 Mohms | – |
| Resistivity 2 | 100.0 Kohms | 1.0 Mohms | – |
| Resistivity 3 | 10.0 Kohms | 100.0 Kohms | – |
| Resistivity 4 | 1.0 Kohms | 10.0 Kohms | – |
| Resistivity 5 | 100 Ohms | 1.0 Kohms | – |
| Temperature | 1000 Ohms | 3.0 Kohms | 66 Kohms |

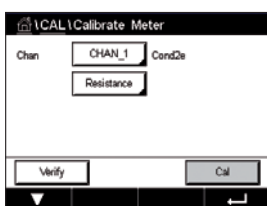

Press the input field in the second line to select Resistance.

Press the Cal button.

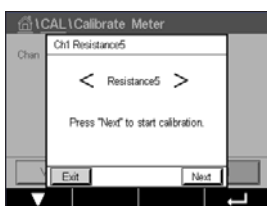


Press the Next button to start the calibration process.




Connect source 1 to input terminals. Each resistance range consists of a two-point calibration.

Press the Next button to continue.

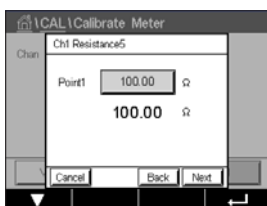


Press input field for Point1 to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button and the transmitter will take over the value.

The second line shows the current value.

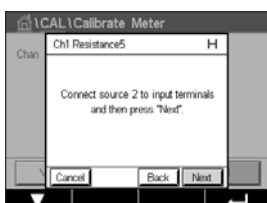


Connect source 2 to input terminals.

Press the Next button to continue.

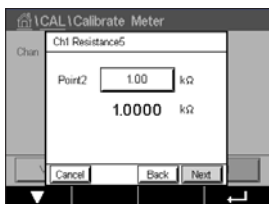

Press input field for Point2 to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button to accept the value.

The second line shows the current value.

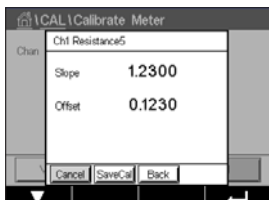


The display shows the value for the slope and the offset as the result of the calibration.

Select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

Use the Back button to go one step back in the calibration procedure.


6.14.2 Temperature (Analog Sensors only)

Temperature is performed as a three point calibration. The table in chapter "6.14.1 Resistance (Analog Sensors only)" on page 64 shows the resistance values of these three points.

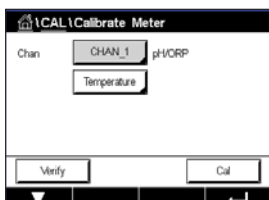


Press the input field in the second line to select Temperature.

Press the Cal button.

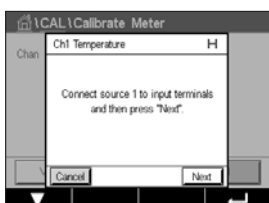


Connect source 1 to input terminals. Press the Next button to start the calibration process.

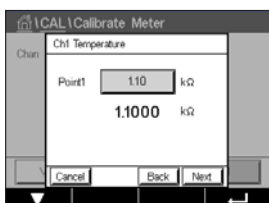

Press input field for Point1 to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button and the transmitter will take over the value.


The second line shows the current value.

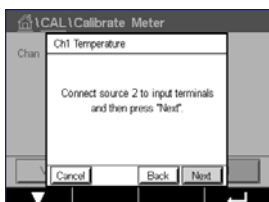

Connect source 2 to input terminals.

Press the Next button to continue.

Repeat the calibration procedure for Point2 and Point3 as for Point1.

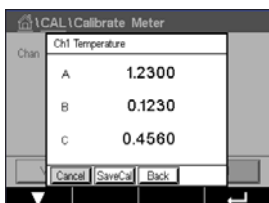


The display shows the result of the calibration.

Select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

Use the Back button to go one step back in the calibration procedure.


6.14.3 Voltage (Analog Sensors only)

Voltage calibration is performed as a two-point calibration.

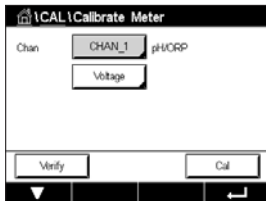


Press the input field in the second line to select Temperature.

Press the Cal button.

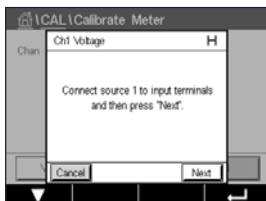


Connect source 1 to input terminals. Press the Next button to start the calibration process.

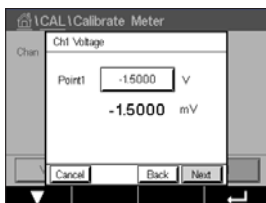

Press input field for Point1 to enter the calibration point. The M400 displays a keypad for modifying the value. Press the \leftarrow button to accept the value.


The second line shows the current value.

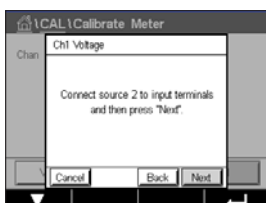

Connect source 2 to input terminals.

Press the Next button to continue.

Repeat the calibration procedure for Point2 and Point3 as for Point1.

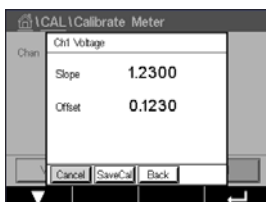

The display shows the result of the calibration.

For Analog sensors select "SaveCal" or "Cancel" to finish calibration. See chapter "6.1.2 Terminate Sensor Calibration" on page 37.

Use the Back button to go one step back in the calibration procedure.

6.14.4 Current (Analog Sensors only)

Current calibration is performed as a two-point calibration.

Perform current calibration according to chapter "6.14.3 Voltage (Analog Sensors only)" on page 66.

6.14.5 Rg (Analog Sensors only)

Rg Diagnostic calibration is performed as a two-point calibration.

Perform current calibration according to chapter "6.14.3 Voltage (Analog Sensors only)" on page 66.


6.14.6 Rr (Analog Sensors only)

Rr Diagnostic calibration is performed as a two-point calibration.

Perform current calibration according to chapter "6.14.3 Voltage (Analog Sensors only)" on page 66.

6.15 Analog Output Calibration

PATH:  \ CAL \ Calibrate Analog Outputs

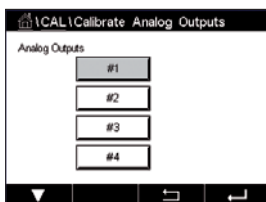

Each analog output can be calibrated at 4 and 20 mA. Select the desired output signal for calibration by pressing the #1 button for output signal 1, #2 for output signal 2, etc.

Connect an accurate milliamp meter to the analog output terminals and then adjust the 5-digit number in the display until the milliamp meter reads 4.00 mA and repeat for 20.00 mA.

As the 5-digit number is increased the output current increases and as the number is decreased the output current decreases. Thus coarse changes in the output current can be made by changing the thousands or hundreds digits and fine changes can be made by changing the tens or ones digits.


After adjusting both values press the Next button to start the calculation of the calibration results.

The display shows the calibration slope and zero point as the result of the output signal calibration.

Select "SaveCal" or "Cancel" to finish calibration. See "6.1.2 Terminate Sensor Calibration" on page 37.

6.16 Analog Input Calibration

PATH:  \ CAL \ Calibrate Analog Inputs

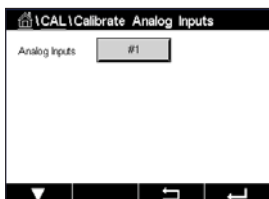

Analog input can be calibrated at 4 and 20 mA by pressing the #1 button.

Connect an 4 mA signal to the analog input terminals. Press the Next button.

Enter the right value for the input signal (**Point1**).

Press the Next button to go on with the calibration.

Connect an 20 mA signal to the analog input terminals. Press the Next button.

Enter the right value for the input signal (**Point2**).

Press the Next button to go on with the calibration.

The display shows the calibration slope and zero point as the result of the input signal calibration.


Selecting Cancel will discard the entered values. Pressing SaveCal will making the entered values the current ones.

If "SaveCal" is chosen, "Calibration Saved Successfully" is displayed.

6.17 Maintenance

PATH:  \ CAL \ Maintenance

The different channels of the M400 Transmitter can be switched manually into HOLD state. Furthermore a cleaning cycle can be started/stopped manually.

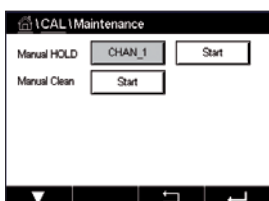

Press Start button for **Manual HOLD** to activate the HOLD state for the selected channel. To deactivate the HOLD state again, press the Stop button, which is now displayed instead of the Start button.

Press the Start button for **Manual Clean** to switch the cleaning relay to the state for starting a cleaning cycle. To switch back the relay press the Stop button, which is now displayed instead of the Start button.

7 Configuration


For the menu structure refer to chapter “3.2 Menu Structure” on page 19.

7.1 Measurement

PATH:  \ CONFIG \ Meas

7.1.1 Channel Setup

PATH:  \ CONFIG \ Meas \ Channel Setup

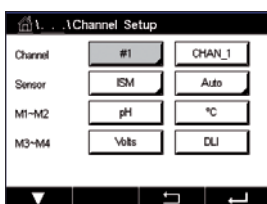

Press the right input field in the line of the setting for **the transmitter**. A parameter for the corresponding channel is chosen through pressing the according field.

If Auto is selected, M400 Transmitter automatically recognizes the ISM sensor type. The channel can also be fixed to a certain measurement parameter, depending on the type of transmitter.

7.1.2 Analog sensor

Select sensor type Analog.

Available measurement types are (depends on transmitter type):

| | M400 Type 1 | M400 Type 2/ M400 4-wire FF | M400 Type 3 |
|--|-------------|--------------------------------|-------------|
| | Analog | Analog | Analog |
| pH/ORP | • | • | • |
| pH/pNa | – | – | – |
| UniCond 2-e/4-e | – | – | – |
| Conductivity 2-e | • | • | • |
| Conductivity 4-e | • | • | • |
| Amp. dissolved oxygen ppm/ppb/trace | – | •/• ¹⁾ /– | •/•/• |
| Opt. dissolved oxygen ppm/ppb | – | –/– | –/– |
| Amp. O ₂ gas ppm/ppb/trace | – | –/–/– | •/•/• |
| Opt. O ₂ gas ppm | – | – | – |
| Dissolved ozone | – | • | • |
| Dissolved carbon dioxide | – | • | • |
| CO ₂ hi | – | – | – |
| GPro 500 TDL | – | – | – |

1) M400 4-wire FF supports Ingold Amp. DO ppb sensors.

7.1.3 ISM sensor

Select sensor type ISM.

If an ISM sensor is connected, the transmitter automatically (Parameter = Auto) Recognizes the type of sensor. You can also fix the transmitter to a certain measurement parameter e.g. "pH", depending on the type of transmitter you have.

| | M400 Type 1 | M400 Type 2/ M400 4-wire FF | M400 Type 3 |
|--|-------------|--------------------------------|-------------|
| | ISM | ISM | ISM |
| pH/ORP | • | • | • |
| pH/pNa | • | • | • |
| UniCond 2-e/4-e | • | • | • |
| Conductivity 2-e | – | – | – |
| Conductivity 4-e | • | • | • |
| Amp. dissolved oxygen ppm/ppb/trace | – | •/• ¹⁾ /– | •/•/• |
| Opt. dissolved oxygen ppm/ppb | – | •/• ²⁾ | •/• |
| Amp. O ₂ gas ppm/ppb/trace | – | –/–/– | •/•/• |
| Opt. O ₂ gas ppm | – | – | • |
| Dissolved ozone | – | • | • |
| Dissolved carbon dioxide | – | • | • |
| CO ₂ hi | – | – | • |
| GPro 500 TDL | – | – | • |

1) M400 4-wire FF supports Ingold Amp. DO ppb sensors.

2) Thornton high performance dissolved oxygen and pure water optical sensors only.

Enter the name with a maximum length of 6 characters for the channel through pressing the input field in the line **Descriptor**. The name of the channel will always be displayed. The name will also be displayed on the Start Screen and Menu Screen.

Choose one of the measurements **M1 to M4** (e.g. for measuring value M1 the left button, for measuring M2 the right button in the corresponding line).

Select in the input field for **Measurement** the desired parameter to show.



NOTE: Beside the parameters pH, O₂, T, etc. also the ISM values DLI, TTM and ACT can be linked to the measurements.

Choose **Range factor** of the measuring value. Not all parameters allow a modification of the range.

The menu **Resolution** allows the setting of the resolution for the measurement. The accuracy of the measurement is not effected by this setting. Possible setting are 1, 0.1, 0.01, 0.001.


Selected the menu **Filter**. The averaging method (noise filter) for the measurement can be selected. The options are None (default), Low, Medium, High, Special and Custom.

| Option | Description |
|---------|--|
| None | No averaging or filtering |
| Low | Equivalent to a 3 point moving average |
| Medium | Equivalent to a 6 point moving average |
| High | equivalent to a 10 point moving average |
| Special | Averaging depending on signal change (normally High averaging, but Low averaging for large changes in input signal) |
| Custom | 1 point to 15 points moving average selection |

7.1.4 Parameter Related Settings

PATH:  \ CONFIG \ Meas \ Parameter Setting


Measuring and calibration parameters can be set for the parameters pH, conductivity and oxygen.

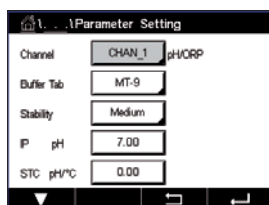

Depending on the selected channel and assigned sensor the measuring and calibration parameters are displayed.

See the following explanation to get more details about the different parameter settings.

7.1.4.1 Conductivity Settings

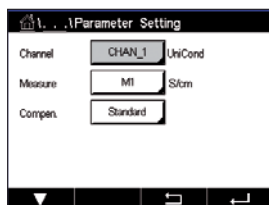

Select measurement (M1-M4). For more information regarding measurements see chapter "7.1.1 Channel Setup" on page 69.

If the selected measurement can be temperature compensated, the compensation method may be selected.



NOTE: During calibration, the compensation method must also be selected. (see chapter "6.2 Calibration of UniCond 2-e and UniCond 4-e Sensors (ISM Sensors only)" on page 37 and chapter "6.3 Calibration of Cond2e Sensors or Cond4e Sensors" on page 45).

Press **Compen.** to select the desired temperature compensation method. Choices are "None", "Standard", "Light 84", "Std 75 °C", "Linear 25°C", "Linear 20°C", "Glycol.5", "Glycol1", "Cation", "Alcohol" and "Ammonia".

None does not make any compensation of the measured conductivity value. The uncompensated value will be displayed and proceeded.

Standard compensation includes compensation for non-linear high purity effects as well as conventional neutral salt impurities and conforms to ASTM standards D1125 and D5391.

Light 84 compensation matches the high purity water research results of Dr. T.S. Light published in 1984. Use only if your institution has standardized on that work.

Std 75 °C compensation is the Standard compensation algorithm referenced to 75 °C. This compensation may be preferred when measuring Ultrapure Water at an elevated temperature (Resistivity of ultrapure water compensated to 75 °C is 2.4818 Mohm-cm.)

Linear 25 °C compensation adjusts the reading by a coefficient or factor expressed as %/°C (deviation from 25 °C). Use only if the solution has a well-characterized linear temperature coefficient. The factory default setting is 2.0%/°C.

Linear 20 °C compensation adjusts the reading by a coefficient or factor expressed as %/°C (deviation from 20 °C). Use only if the solution has a well-characterized linear temperature coefficient. The factory default setting is 2.0%/°C.

Glycol.5 compensation matches the temperature characteristics of 50% ethylene glycol in water. Compensated measurements using this solution may go above 18 Mohm-cm.

Glycol1 compensation matches the temperature characteristics of 100% ethylene glycol. Compensated measurements may go well above 18 Mohm-cm.

Cation compensation is used in power industry applications measuring the sample after a cation exchanger. It takes into account the effects of temperature on the dissociation of pure water in the presence of acids.

Alcohol compensation provides for the temperature characteristics of a 75% solution of isopropyl alcohol in pure water. Compensated measurements using this solution may go above 18 Mohm-cm.


Ammonia compensation is used in power industry applications for specific conductivity measured on samples using ammonia and/or ETA (ethanolamine) water treatment. It takes into account the effects of temperature on the dissociation of pure water in the presence of these bases.



NOTE: If compensation mode “Linear 25 °C” or “Linear 20 °C” has been chosen, the coefficient for the adjustment of the reading can be modified. In this case an additional input field will be displayed.

Press the input field for **Coef.** and adjust the coefficient or factor for the compensation.

7.1.4.2 pH Settings

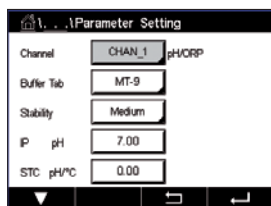

If a pH sensor is connected while during the channel setup (see chapter “7.1.1 Channel Setup” on page 69) Auto has been chosen the parameters Buffer Tab, Stability, IP, STC and calibration temperature as well as the displayed units for slope and/or zero point can be set or adjusted. The same parameters will be displayed if during the channel setup not Auto but pH/ORP has been set.

Select the buffer through the parameter **Buffer Tab**.

For automatic buffer recognition during calibration, select the buffer solution set that will be used: Mettler-9, Mettler-10, NIST Tech, NIST Std = JIS Std, HACH, CIBA, MERCK, WTW, JIS Z 8802 or None. See chapter “16 Buffer tables” on page 125 for buffer values. If the auto buffer feature will not be used or if the available buffers are different from those above, select None.



NOTE: For dual membrane pH electrodes (pH/pNa) buffer Na+ 3.9M (see chapter “16.2.1 Mettler-pH/pNa buffers (Na+ 3.9M)” on page 130.

Select the required **Stability** of the measuring signal during the calibration procedure. Choose manual if the user will decide when a signal is stable enough to complete the calibration. Select Low, Medium or Strict if an automatic stability control of the sensor signal during calibration through the transmitter should be done.

If the parameter stability is set to medium (default) the signal deviation has to be less than 0.8 mV over a 20 second interval to be recognized by the transmitter as stable. The calibration is done using the last reading. If the criteria is not met within 300 seconds then the calibration times out and the message “Calibration Not Done” is displayed.

Adjust the parameter **IP pH**.

IP is the isothermal point value (Default = 7.000 for most applications). For specific compensation requirements or non standard inner buffer value, this value can be changed.


Adjust the value of the parameter **STC pH/°C**.

STC is the solution temperature coefficient in units of pH/°C referenced to the defined temperature. (Default = 0.000 pH/°C for most applications). For pure waters, a setting of -0.016 pH/°C should be used. For low conductivity power plant samples near 9 pH, a setting of -0.033 pH/°C should be used.

If the value for STC is $\neq 0.000$ pH/°C an additional input field for the reference temperature will be displayed.

The value for **pH Ref Temperature** indicates to which temperature the solution temperature compensation is referenced. The displayed value and the output signal is referenced to this temperature. Most common reference temperature is 25°C.

7.1.4.3 Settings for Oxygen Measurement Based on Amperometric Sensors

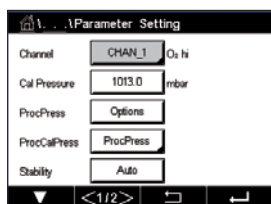

If an amperometric oxygen sensor is connected while during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) Auto has been chosen the parameters CalPressure, ProcPressure, ProcCalPress, Stability, Salinity, RelHumidity, UpolMeas and UpolCal can be set or adjusted. The same parameters will be displayed if during the channel setup not Auto but O₂ hi, or O₂ lo has been set.

Enter the value for the calibration pressure through the parameter **CalPressure**.

NOTE: For a modification of the unit for the calibration pressure press U on the displayed keypad.

Press the Option button for the parameter **ProcPressure** and select the how to get applying process pressure through choosing the **Type**.


The applied process pressure can be entered by choosing Edit or measured over the analog input of the M400 by choosing Ain_1.

If Edit has been chosen an input field for entering the value manually is displayed on the screen. In case that Ain_1 has been selected two input fields are displayed to enter the start value (4 mA) and the end value (20 mA) of the range for the 4 to 20 mA input signal.

For the algorithm of the process calibration the applied pressure has to be defined. Select the pressure through the parameter **ProcCalPress**. For the process calibration the value of the process pressure (ProcPress) or the calibration pressure (CalPress) can be used.

Select the required **Stability** of the measuring signal during the calibration procedure. Choose Manual if the user will decide when a signal is stable enough to complete the calibration. Select Auto and an automatic stability control of the sensor signal during calibration through the transmitter will be done.

Additional settings can be done by navigating to the next page of the menu.

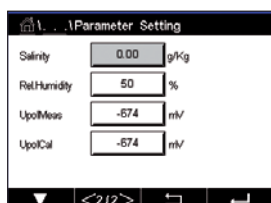

The **Salinity** of the measured solution can be modified.

In addition the relative humidity (button **Rel.Humidity**) of the calibration gas can also be entered. The allowed values for relative humidity are in the range 0% to 100%. When no humidity measurement is available, use 50% (default value).


The polarization voltage of amperometric oxygen sensors in the measuring mode can be modified through the parameter **UpolMeas**. For entered values 0 mV to -550 mV the connected sensor will be set to a polarization voltage of -500mV. If the entered value is less then -550 mV, the connected sensor will set to a polarization voltage of -674 mV.

The polarization voltage of amperometric oxygen sensors for calibration can be modified through the parameter **UpolCal**. For entered values 0 mV to -550 mV the connected sensor will be set to a polarization voltage of -500mV. If the entered value is less then -550mV, the connected sensor will set to a polarization voltage of -674mV.

NOTE: During a process calibration, the polarization voltage UpolMeas, defined for the measuring mode, will be used.

NOTE: If a one-point calibration is executed, the transmitter sends the polarization voltage, valid for the calibration, to the sensor. If the polarization voltage for the measuring mode and calibration mode is different, the transmitter will wait 120 seconds before starting the calibration. In this case the transmitter will also go after the calibration for 120 seconds to the HOLD Mode, before returning to the measuring mode again.

7.1.4.4 Settings for Oxygen Measurement Based on Optical Sensors

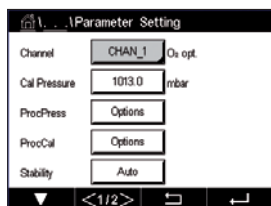

If an optical oxygen sensor is connected while during the channel setup (see chapter “7.1.1 Channel Setup” on page 69) Auto has been chosen the parameters CalPressure, ProcPressure, ProcCalPress, Stability, Salinity, RelHumidity, Sample Rate, LED Mode and Toff can be set or adjusted. The same parameters will be displayed if during the channel setup not Auto but Optical O₂ has been set.

Enter the value for the calibration pressure through the parameter **CalPressure**.

Press the button Option for the parameter **ProcPress** and select the how to get applying process pressure through pressing the according button in the line **Type**.


The applied process pressure can be entered by choosing Edit or measured over the analog input of the M400 by choosing AIN_1.

If Edit has been chosen an input field for entering the value manually is displayed on the screen. In case that AIN_1 has been selected two input fields are displayed to enter the start value (4mA) and the end value (20 mA) of the range for the 4 to 20 mA input signal.

For the algorithm of the process calibration the applied pressure has to be defined. Select the pressure through the parameter **ProcCal**. For the process calibration the value of the process pressure (ProcPress) and the value of the calibration pressure (CalPress) can be used. Select between Scaling and Calibration for the process calibration. If Scaling has been chosen, the calibration curve of the sensor will be untouched, but the output signal of the sensor will be scaled. In case of calibration value <1%, the offset of the sensor output signal will be modified during scaling, for value >1% the slope of the sensor output will be adjusted. For further information about scaling refer to the sensor manual.

Selecting the required **Stability** of the measuring signal during the calibration procedure. Choose Manual if the user will decide when a signal is stable enough to complete the calibration. Select Auto and an automatic stability control of the sensor signal during calibration through the transmitter will be done.

Additional settings can be done by navigating to the next page of the menu.

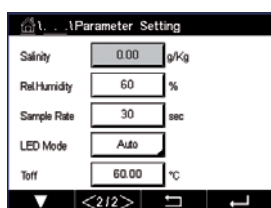

The **Salinity** of the measured solution can be modified.

In addition the relative humidity (button **Rel.Humidity**) of the calibration gas can also be entered. The allowed values for relative humidity are in the range 0% to 100%. When no humidity measurement is available, use 50% (default value).

Adjust the required **Sample Rate** of the optical sensor during measurement. The time interval from one measuring cycle of the sensor to the next can be adjusted i.e. adapted to the application. A higher value will increase the life time of the OptoCap of the sensor.

Select the **LED Mode** of the sensor. There are the following options.

Off: LED is permanently switched off.

On: LED is permanently switched on.

Auto: The LED is switched on as long as the measured media temperature is smaller then Toff (see next value) or switched off through a digital input signal (see chapter "7.10 Digital Inputs" on page 90).



NOTE: If the LED is switched off, no oxygen measurement is performed.


Enter the limit for the measuring temperature to switch off the LED of the sensor automatically for the M400 through the parameter **Toff**.

If the media temperature is higher then Toff, the LED will be switched off. The LED will be switched on as soon as the media temperature falls below Toff -3 K. This function give the option to increase the lifetime of the OptoCap by switching off the LED during SIP or CIP cycles.



NOTE: This function is only active if the LED Mode is set to "Auto".

7.1.4.5 Dissolved Carbon Dioxide Settings

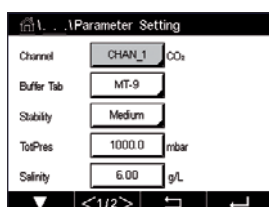

If an dissolved carbon dioxide sensor is connected while during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) Auto or CO₂ has been chosen, the buffer used for calibration and the parameters stability, salinity, HCO₃, TotPres can be set resp. adjusted.


Select the buffer through the parameter **Buffer Tab**. For automatic buffer recognition during calibration, select buffer solution Mettler-9 if it will be used. If the auto buffer feature will not be used or if the available buffer are different from Mettler-9 select None.

Select the required **Stability** of the measuring signal during the calibration procedure. Choose manual if the user will decide when a signal is stable enough to complete the calibration. Select Low, Medium or Strict if an automatic stability control of the sensor signal during calibration through the transmitter should be done.

If the unit for the measured dissolved carbon dioxide is %sat, the pressure during the calibration resp. measurement has to be considered. This will be done by setting the parameter **TotPres**. If another unit then %sat has been selected, the result will not be influenced by this parameter.

The **Salinity** describes the total amount of solved salts in the CO₂ electrolyte of the sensor connected to the transmitter. It is a sensor specific parameter. The default value (28.00 g/L) is valid for the InPro5000i. Do not change this parameter if the InPro 5000 i will be used.


Additional settings can be done by navigating to the next page of the menu.

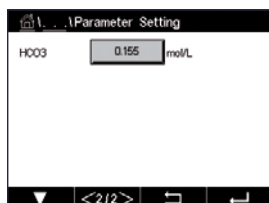

The parameter **HCO₃** describes the concentration of hydrogen carbonate in the CO₂ electrolyte of the sensor connected to the transmitter. It is also a sensor specific parameter. The default value 0.050 Mol/L is valid for the InPro 5000 i. Do not change this parameter if the InPro 5000 i will be used.

7.1.4.6 Settings for Thermal Conductivity Dissolved CO₂ Measurement (CO₂ hi)

If during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) the parameter CO₂ Hi has been chosen, the parameters stability (manual/auto) and CO₂ solubility (CO₂-solubility and Temperature Factor), be set resp. adjusted.

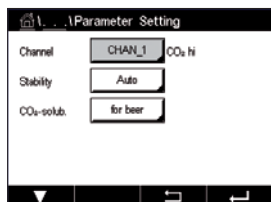

Select the required **Stability** of the measuring signal during the calibration procedure. Choose manual if the user will decide when a signal is stable enough to complete the calibration. Select Auto if an automatic stability control of the sensor signal during calibration through the transmitter should be done.

The sensor offers a choice of CO₂ **Solubility's** for measurement in beer, water and cola. The cola setting is to be used with carbonated soft drinks. For other beverages the user has the possibility to enter individual values for CO₂ solubility and temperature factors.

Default values for measurement in beer (valid for temperatures – 5 ... 50 °C):

CO₂ solubility (A): 1.420 g/L
Temp. factor (B): 2485

Values for pure water:


CO₂ solubility (A): 1.471 g/L
Temp. factor (B): 2491

Values for cola:

CO₂ solubility (A): 1.345 g/L
Temp. factor (B): 2370



NOTE: The sensor is delivered factory calibrated and is set up to measure in beer as the default.

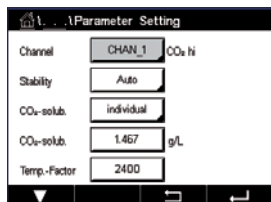

For beverages where the user knows the exact CO₂ solubility and the temperature factor the values can be changed **individually**.

If the user desires to evaluate the solubility (**CO₂-solub.**) and temperature factors (**Temp.-Factor**) they can be evaluated with the following formulas.

$$\text{HCO}_2 = A * \exp (B * (1 / T - 1 / 298.15))$$

$$c\text{CO}_2 = \text{HCO}_2 * p\text{CO}_2$$

HCO₂ : Calculated CO₂ Solubility (Henry constant) at measured process temp.


A: Solubility of CO₂ (g / L at 25 °C)

B: Temperature factor (valid for – 5 ... 50 °C)


cCO₂: Calculated CO₂ concentration in g/l or V/V

7.1.4.7 Settings for Tunable Diode Laser (TDL) Analyzer

(PATH:  \ Config \ Measurement \ TDL quick setup)

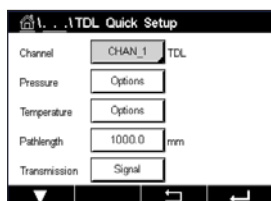

If a TDL analyzer is connected, while during the channel setup (see chapter “7.1.1 Channel Set-up” on page 69) Auto has been chosen, the parameters Pressure, Temperature and Path length can be set or adjusted. The same parameters will be displayed if during the channel set-up not Auto but TDL has been set.

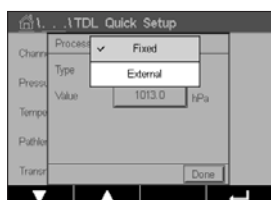

Press the button for Pressure.

- External: current external pressure value coming from a pressure transducer of 4.. 20 mA analog output
 - Fixed: pressure compensation uses a fixed value to be set manually.
- Note:** if this pressure compensation mode is selected, a considerable gas concentration measurement error resulting from a non- realistic pressure value can take place.

If External compensation is selected, then the minimum (4 mA) and maximum (20 mA) analog output signals from the pressure transducer must be mapped to the corresponding Analog input of the TDL. Key in the minimum and maximum values of the pressure in the following units:

- hPa – mmHg – mbar
- psi – kPa


In general, METTLER TOLEDO recommends the use of absolute pressure transducers for more accurate signal compensation over a broad pressure range.

If, however, small pressure variations around atmospheric pressure are to be expected, relative pressure sensors will produce better results; but the variations of the underlying barometric pressure will be ignored.

For relative pressure sensors, the minimum and maximum values must be mapped so that the TDL can interpret the analog pressure signal as “absolute”, i.e. a fixed barometric pressure of 1013 mbar (for example) has to be added to the mapped values.

If Fixed compensation is selected, the fixed pressure value with which the measurement signal will be calculated has to be keyed in manually. For the fixed pressure, the following units can be used:


- hPa – mmHg – mbar
- psi – kPa

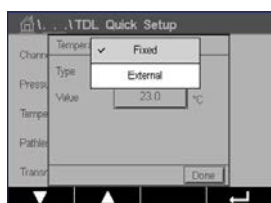

Press the button for Temperature.

If External compensation is selected, then the minimum (4 mA) and maximum (20 mA) analog output signals from the temperature transducer must be mapped to the corresponding Analog input of the TDL. Key in the minimum and maximum values of the temperature in °C.

If Fixed compensation is selected, the fixed temperature value with which the measurement signal will be calculated has to be keyed in manually. For the fixed temperature, only °C can be used.

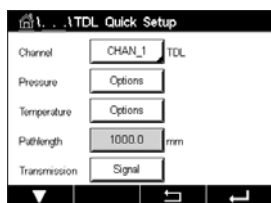

Last, select the initial optical path length corresponding to the probe length installed:


- 290 mm probe: 200 mm
- 390 mm probe: 400 mm
- 590 mm probe: 800 mm

This initial value is valid when instrument purging on the instrument and on the process side is running. Depending on the process conditions and after the optimum of the process purging flow has been found (see next chapter), this value may have to be slightly adapted.

7.1.4.8 Setting the correct process side purging

The flow rate of the purging will affect the effective path length and consequently the measurement value.

Therefore the following procedure should be used. Start with a very high flow rate and gradually decrease it. The measurement value will then start at a low value and increase with decreasing purge flow. At some point it will level out and stay constant for a while and then again start increasing. Choose a purge flow in the middle of the constant region.

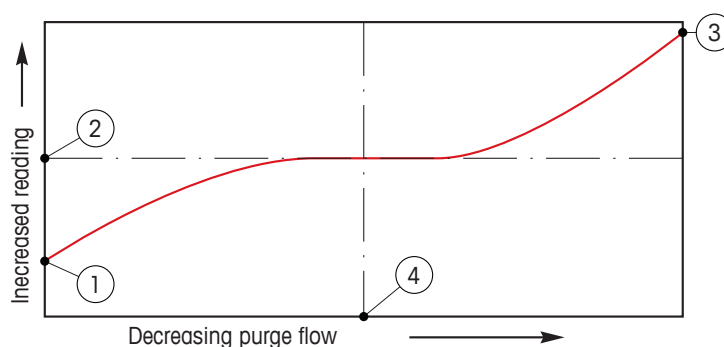

Optimizing the purge flow

On the x-axis there is purge flow and on the y-axis there is the instrument concentration reading.

- 1 Concentration reading with high purge flow. The path length is now shorter than the effective path length since the purge tubes is completely filled with purging gas and some of the purging gas is flowing into the measurement path.
- 2 Concentration reading with optimized purge flow. The path length is now equal to the effective path length since the purge tubes are completely filled with purge gas. See the illustration below.
- 3 Concentration reading with no purge flow. The path length is now equal to the nominal path length since the probe is completely filled with process gas.
- 4 The optimized purge flow.




WARNING: Always start purging at maximum flow before starting the process.



WARNING: Purging must always be switched on in order to avoid dust deposition onto the optical surfaces.

7.1.5 Concentration Curve Table


To specify a concentration curve for customer-specific solutions, up to 5 concentration values can be edited in a matrix together with up to 5 temperatures. To do so the desired values are edited under the concentration curve table menu. Beside the temperature values, the conductivity and concentration values for the corresponding temperature are edited. The concentration curve can be selected resp. used in combination with conductivity sensors.

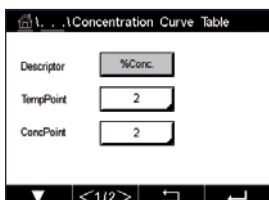

Enter the name with a maximum length of 6 characters for the concentration curve through pressing the input field in the line **Descriptor**.

Enter the amount of desired temperature points (**TempPoint**) and concentration points (**ConcPoint**).

The different values can be entered by navigating to the next page of the menu.

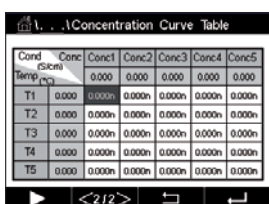

Enter the values for temperature (**T1...T5**), concentration (**Conc1...Conc5**) and the corresponding conductivity through pressing the according input field. The unit for the value of the conductivity can be adjusted as well in the according input field.



NOTE: The values for the temperature have to increase from T1 to T2 to T3, etc. The values for the concentration have to increase from Conc1 to Conc2 to Conc3, etc.



NOTE: The conductivity values at the different temperatures have to increase or decrease from Conc1 to Conc2 to Conc3, etc. Maxima and/or minima are not permitted. If the conductivity values at T1 are increasing with the different concentrations, they have to increase also at the other temperatures. If the conductivity values at T1 are decreasing with the different concentrations, they have to decrease also at the other temperatures.

7.2 Temperature Source (Analog Sensors only)

PATH: \ CONFIG \ Meas \ Temperature Source


Source: Auto(default), Pf100, Pf1000, NTC22k, Fixed

The third line shows the related temperature setting. Range: -40 to 200 °C, Default: 25 °C

7.3 Analog Outputs

PATH:  \ CONFIG \ Analog Outputs

See the following explanation to get more details about the different settings for the analog outputs.

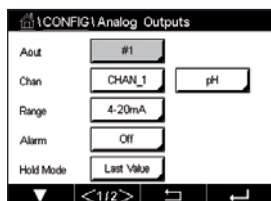

Press the input field in the line of the setting for **Aout** and select the desired output signal for configuration by pressing button #1 for output signal 1, #2 for output signal 2 etc. Press the related button for the assignment of the channel (**Chan**). Select the channel, which has to be linked to the output signal.

Press the button for the assignment of the measuring parameter – based on the selected channel – that has to be linked to the output signal.



NOTE: Besides the measurement values pH, O₂, T, etc. also the ISM values DLI, TTM and ACT can be linked to the output signal.

Select the **Range** for the output signal.

To adjust the value for the analog output signal if an alarm occurs, press the input field in the line for the setting of **Alarm**. Off means, that an alarm has now influence on the output signal.




NOTE: Not only the alarms occurred on the assigned channel will be considered, but every alarm coming up on the transmitter.

The value for the output signal if the transmitter goes into HOLD mode can be defined. It can be chosen between the last value (i.e. the value before the transmitter switched to the HOLD mode) or an fixed value.

Press the input field in the line for the setting of the **HOLD Mode** and select the value. If a fixed value is chosen, the transmitter shows an additional input field.

Additional settings can be done by navigating to the next page of the menu.

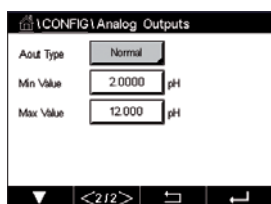

The **Aout Type** can be Normal, Bi-Linear, Auto-Range or Logarithmic. The range can be 4–20 mA or 0–20 mA. Normal provides linear scaling between the minimum and maximum scaling limits and is the default setting. Bi-Linear will also prompt for a scaling value for the mid-point of the signal and allows two different linear segments between the minimum and maximum scaling limits.

Press the button for the **Min Value**, that corresponds with start point of the analog output range.

Press the button for the **Max Value**, that corresponds with end point of the analog output signal.

Depending on the chosen Aout type additional values can be entered.

Bi-Linear will also prompt for a scaling value for the Mid Value of the signal and allows two different linear segments between the defined Min and Max Values.

Auto-Range scaling provides two ranges of output. It is designed to work with a PLC to provide a wide measurement range at the high end of the scale, and a narrower range with high resolution at the low end. Two separate settings are used, one for the maximum limit of the high range and one for the maximum limit of the low range, for the single 0/4-20 mA signal.


Max1 is the maximum limit of the low range on auto-range. The maximum value for the high range on auto-range is set with the Max Value. Both ranges have the same minimum value that is set through Min Value. If the input value is higher then value of Max1, the transmitter switches automatically to the second range. To indicate the currently valid range a relay can be assigned. The relay will be switched if the transmitter changes from on range to the other.

If **Logarithmic** Range was selected, it will prompt for the Max Value and also for the number of decades.

7.4 Set Points

PATH: \ CONFIG \ Set Points

See the following explanation to get more details about the different settings for the set points.

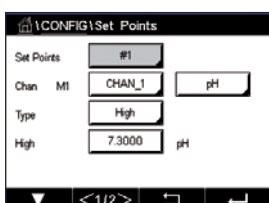

Press the input field in the line of the setting for **Set Point** and select the desired set point for configuration through pressing the button #1 for set point 1, #2 for set point 2 etc..

Press the related button for the assignment of the channel (**Chan**). Select the channel, which has to be linked to the set point.

Press the button for the assignment of the measuring parameter – based on the selected channel – that has to be linked to the set point.

Mx in the display indicates the measurement assigned to the set point. (see chapter 7.1.1 "Channel Setup").



NOTE: Beside the parameters pH, O₂, T, mS/cm, %EP WFI etc. also the ISM values DLI, TTM and ACT can be linked to the set point.


The **Type** of the setpoint can be High, Low, Between, Outside or Off. An "Outside" setpoint will cause an alarm condition whenever the measurement goes above its high limit or below its low limit. A "Between" setpoint will cause an alarm condition to occur whenever the measurement is between its high and low limits.



NOTE: If the type of set point is not Off additional settings can be done. See the following description.

According to the selected type of setpoint, value(s) regarding the limit(s) can be entered.

Additional settings can be done by navigating to the next page of the menu.

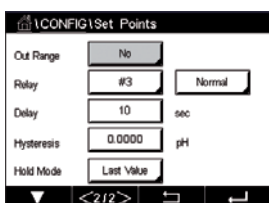

Once configured a relay could be activated if a sensor **Out of Range** condition is detected on the assigned input channel.

To select the desired relay that will be activated if the defined conditions are reached press the input field in the line for the setting of **SP Relay**. If the chosen relay is used for another task, the transmitter shows the message on the screen that there is a Relay Conflict.

The operation mode of the relay can be defined.

Relay contacts are in normal mode until the associated setpoint is exceeded, then the relay is activated and the contact states change. Select Inverted to reverse the normal operating state of the relay (i.e. normally open contacts are in a closed state, and normally closed contacts are in an open state, until the setpoint is exceeded).

Enter the **Delay** time in seconds. A time delay requires the setpoint to be exceeded continuously for the specified length of time before activating the relay. If the condition disappears before the delay period is over, the relay will not be activated.


Enter the value for the **Hysteresis**. A hysteresis value requires the measurement to return within the setpoint value by a specified percentage before the relay is deactivated.

For a high setpoint, the measurement must decrease more than the indicated percentage below the setpoint value before the relay is deactivated. With a low setpoint, the measurement must rise at least this percentage above the setpoint value before the relay is deactivated. For example, with a high setpoint of 100, when this value is exceeded, the measurement must fall below 90 before the relay is deactivated.

Enter the relay **HOLD Mode** of "Off", "Last Value" or "On". This is the state of the relay during HOLD status.

7.5 ISM Setup (ISM Sensors only)

PATH:  \ CONFIG \ ISM Setup

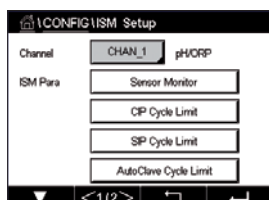


See the following explanation to get more details about the different parameter settings for the ISM Setup.

7.5.1 Sensor Monitor

If a ISM sensor is connected while during the channel setup (see "7.1.1 Channel Setup" on page 69) Auto has been chosen the parameter Sensor Monitor can be set or adjusted. The menu Sensor Monitor will also be displayed if during the channel setup not Auto but one of the mentioned sensors has been set.

Press the button Sensor Monitor.

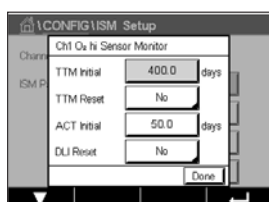

Enter the value for the initial Time To Maintenance interval (**TTM Initial**) in days. The initial value for TTM can be modified according to the application experience.

For pH/ORP sensor the timer estimates when the next cleaning cycle should be performed to keep the best possible measurement performance. The timer is influenced by significant changes on the DLI parameters.

For amperometric oxygen and ozone sensors, the time to maintenance indicates a maintenance cycle for the membrane and electrolyte.

Press the input field for **TTM Reset**. Select Yes if Time To Maintenance (TTM) for the sensor should be reset to the initial value.

Time To Maintenance needs to be reset after the following operations.

pH sensors: manual maintenance cycle on the sensor.

Oxygen or ozone sensor: manual maintenance cycle on the sensor or exchanging of the membrane of the sensor



NOTE: By connecting a sensor, the actual value for TTM of the sensor is read out from the sensor.

Enter the **ACT Initial** value in days. The new value will be loaded down to the sensor after saving the changes.

The Adaptive Calibration Timer (ACT) estimates when the next calibration should be performed to keep the best possible measurement performance. The timer is influenced by significant changes on the DLI parameters. The ACT will be reset to its initial value after a successful calibration. The initial value for the ACT can be modified according to the application experience and loaded down to the sensor.



NOTE: By connecting a sensor, the actual value for the ACT of the sensor is read out from the sensor.

Press the input field for **DLI Reset**. Select Yes if Dynamic Lifetime Indicator (DLI) for the sensor should be reset to the initial value. The reset will be done after saving the changes.

The DLI allows an estimation, when the pH electrode, the inner body of an amperometric oxygen or ozone sensor is at the end of his lifetime, based on the actual stress he is exposed to. The sensor permanently takes the averaged stress of the past days into consideration and is able to increase/decrease the lifetime accordingly.

The following parameters affect the lifetime indicator:

Dynamic parameter

- Temperature
- pH or oxygen value
- Glass impedance (only pH)
- Reference impedance (only pH)

Static parameters

- Calibration history
- Zero and Slope
- CIP/SIP/Autoclaving cycles

The sensor keeps the information stored in the built in electronics and can be retrieved via a transmitter or the iSense asset management suite.

For amperometric oxygen sensors, the DLI is related to the inner-body of the sensor. After exchanging the inner-body perform DLI Reset.



NOTE: By connecting a sensor, the actual values for the DLI of the sensor are read out from the sensor.




NOTE: The menu DLI Reset for pH sensors not available. If the actual value for the DLI of a pH sensor is 0 the sensor has to be replaced.

7.5.2 CIP Cycle Limit

If a pH/ORP, oxygen or conductivity sensor is connected during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) Auto has been chosen the parameter CIP Cycle Limit can be set or adjusted. The menu CIP Cycle Limit will also be displayed if during the channel setup not Auto but one of the mentioned sensors has been set.

Press the button CIP Cycle Limit.

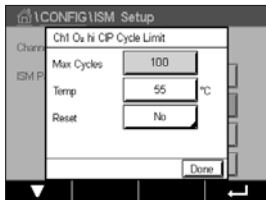

Press the button in the input field for the parameter **Max Cycles** and enter the value for the maximum CIP cycles. The new value will be written to the sensor after saving the changes.

The CIP cycles are counted by the transmitter. If the limit (value for Max Cycles) is reached, an alarm can be indicated and set to a certain output relays.

If the Max Cycles setting is on 0, the counter functionality is turned off.

Press the button in the input field for the parameter **Temp** and enter the temperature, which has to be exceeded, that the a CIP cycle will be counted.

CIP Cycles will be automatically recognized by the transmitter. Since CIP cycles will vary in intensity (duration and temperature) for each application the algorithm of the counter recognizes an increase of the measurement temperature above the level defined through the value for Temp. If the temperature does not decrease below the defined temperature level -10 °C within the next 5 minutes after the first temperature was reached, the counter in question will be incremented by one and also locked for the next two hours. In the case the CIP would last longer than two hours the counter would be incremented by one once more.

Press the input field for **Reset**. Select Yes if CIP counter for the sensor should be reset to 0. The reset will be done after saving the changes.

If an oxygen sensor is connected, the reset should be performed after the following operations.
amperometric sensor: exchanging of the inner-body of the sensor.




NOTE: For pH/ORP sensor the menu Reset is not available. A pH/ORP sensor should be replaced if the number for Max Cycles has been exceeded.

7.5.3 SIP Cycle Limit

If a pH/ORP, oxygen or conductivity sensor is connected during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) Auto has been chosen the parameter SIP Cycle Limit can be set or adjusted. The menu SIP Cycle Limit will also be displayed if during the channel setup not Auto but one of the mentioned sensors has been set.

Press the button SIP Cycle Limit.

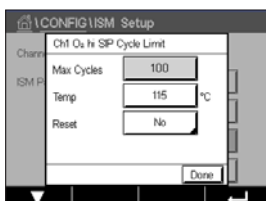

Press the button in the input field for the parameter **Max Cycles** and enter the value for the maximum SIP cycles. The new value will be written to the sensor after saving the changes.

The SIP cycles are counted by the transmitter. If the limit (value for Max Cycles) is reached, an alarm can be indicated and set to a certain output relays.

If the Max Cycles setting is on 0, the counter functionality is turned off.

Press the button in the input field for the parameter **Temp** and enter the temperature, which has to be exceeded, that the a SIP cycle will be counted.

SIP Cycles will be automatically recognized by the transmitter. Since SIP cycles will vary in intensity (duration and temperature) for each application the algorithm of the counter recognizes an increase of the measurement temperature above the level defined through the value for Temp. If the temperature does not decrease below the defined temperature level - 10°C within the next 5 minutes after the first temperature was reached, the counter in question will be incremented by one and also locked for the next two hours. In the case the SIP would last longer than two hours the counter would be incremented by one once more.

Press the input field for **Reset**. Select Yes if SIP counter for the sensor should be reset to 0. The reset will be done after saving the changes.

If an oxygen sensor is connected, the reset should be performed after the following operations.
Amperometric sensor: exchanging of the inner-body of the sensor.




NOTE: For pH/ORP sensor the menu Reset is not available. A pH/ORP sensor should be replaced if the number for Max Cycles has been exceeded.

7.5.4 AutoClave Cycle Limit

If a pH/ORP, amperometric oxygen is connected during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) Auto has been chosen the parameter AutoClave Cycle Limit can be set or adjusted. The menu AutoClave Cycle Limit will also be displayed if during the channel setup not Auto but one of the mentioned sensors has been set.

Press the button AutoClave Cycle Limit.

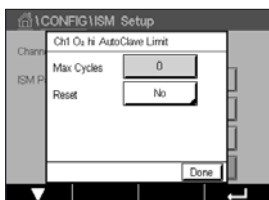

Press the button in the input field for the parameter **Max Cycles** and enter the value for the maximum AutoClave cycles. The new value will be written to the sensor after saving the changes.

If the Max Cycles setting is on 0, the counter functionality is turned off.

Since during the autoclaving cycle the sensor is not connected to the transmitter, you will be asked after every sensor connection, whether the sensor was autoclaved or not. According to your selection, the counter will be incremented or not. If the limit (value for Max Cycles) is reached, an alarm can be indicated and set to a certain output relay. Press the input field for **Reset**. Select Yes if the AutoClave counter for the sensor should be reset to 0. The reset will be done after saving the changes.


If an oxygen sensor is connected, the reset should be performed after the following operations.
Amperometric sensor: exchanging of the inner-body of the sensor.



NOTE: For pH/ORP sensor the menu Reset is not available. A pH/ORP sensor should be replaced if the number for Max Cycles has been exceeded.

7.5.5 DLI Stress Adjustment

If a pH/ORP is connected during the channel setup (see chapter "7.1.1 Channel Setup" on page 69) Auto has been chosen the parameter DLI Stress Adjustment can be adjusted. With this setting the user can adjust the sensor sensitivity to the stress of his specific application for the DLI calculation.

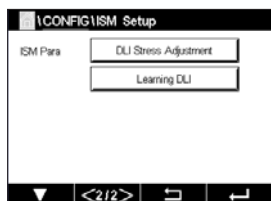

Browse to page 2 of "ISM Setup".

Press the button **DLI Stress Adjustment**.

Select between low / medium / high for the **Type** of DLI Stress Adjustment.

LOW: DLI extended (-30% sensitivity)

MEDIUM: standard DLI (default)


HIGH: DLI reduced (+30% sensitivity)

Press \leftarrow to accept the setting.

7.5.6 SAN Cycle Parameters

If an ozone sensor is connected, values for the following SAN Cycle Parameters can be set, Max Cycles (the maximum number of sanitization cycles), Conc. Max (the maximum allowed O_3 concentration), Conc. Min (the minimum allowed O_3 concentration), Cycle Time (length of cycle), and Reset.

Press the button SAN Cycle Parameters.

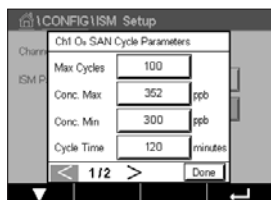

Press the input field next to Max Cycles and enter the value for the maximum SAN cycles. Press \leftarrow to accept the value. The new value will be written to the sensor after saving the changes.

The SAN cycles are counted by the transmitter. If the limit (value for Max Cycles) is reached, an alarm can be configured. If the Max Cycles setting = 0, the counter functionality is turned off.

Press the input field next to Conc. Max and enter the ozone concentration above which a sanitization cycle is to be detected. Press \leftarrow to accept the value.

Press the input field next to Conc. Min. Enter the value for the ozone concentration below which a sanitization cycle is no longer detected. Press \leftarrow to accept the value

Press the input field next to Cycle Time. Enter the value for the time, the ozone concentration has to be higher than the Conc. Min value after the Conc. Max value has been exceeded to count a sanitization cycle. Press \leftarrow to accept the value.


Press the input field next to Reset. Select Yes to reset the sanitization counter to zero. This is typically performed after sensor replacement. The reset will be done after saving the changes

Press \leftarrow to exit the menu SAN Cycle Parameters.

7.5.7 Reset Counters for UniCond 2-e Sensors

For UniCond 2-e sensors, the following counters can be reset: High Temp and High Conductivity.

Press the button Reset Counters.

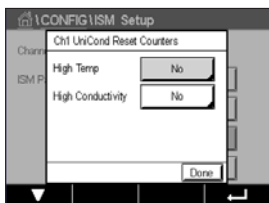


Select Yes for the desired counter to be reset and press enter. The reset will be done after saving the changes.

Press \leftarrow to exit the menu Reset Counters.

7.5.8 Set Calibration Interval for UniCond 2-e Sensors

For UniCond 2-e sensor the Cal Interval (calibration interval) can be set.

Press the button Cal Interval.

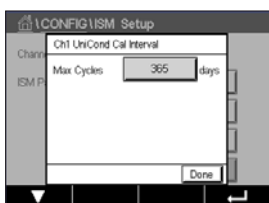

Press the input field next to **Cal Interval** and enter the value for the calibration interval. Based on this value the Time To Calibration (TTCal) will be calculated by the transmitter. Press \leftarrow to accept the value. The new value will be written to the sensor after saving the changes.

Press \leftarrow to exit the menu Cal Interval.

7.6 General Alarm

PATH: \leftarrow \ CONFIG \ General Alarm

See the following explanation to get more details about the different settings for General Alarm.



Press the button Event in the line of the settings for **Option** and select the events, that should be considered for an alarm.


To activate a relay if the defined conditions are reached press the input field in the line for the settings of **Relay**. Only relay 1 can be assigned to general alarm. For general alarms the operation mode of the assigned relay is always inverted.

Enter the **Delay** time in seconds. A time delay requires the setpoint to be exceeded continuously for the specified length of time before activating the relay. If the condition disappears before the delay period is over, the relay will not be activated.

7.7 ISM / Sensor Alarm

PATH:  \ CONFIG \ ISM / Sensor Alarm

See the following explanation to get more details about the different settings for ISM / Sensor Alarm.

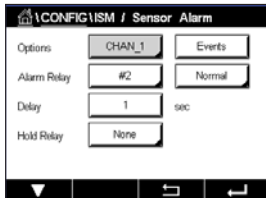

Depending on the assigned sensor the **Events** that will be considered for generating an alarm can be selected. Some alarms will be considered in any case and not have to be selected or deactivated.

To select the desired relay that will be activated if an event has taken place press the input field in the line for the settings for **Relay**.

The operation mode of the relay can be defined.


Relay contacts are in normal mode until one of the selected events has taken place. Then the relay is activated and the contact states change. Select Inverted to reverse the normal operating state of the relay (i.e. normally open contacts are in a open state, and normally closed contacts are in a closed state if an event has taken place).

Enter the **Delay** time in seconds. A time delay requires the event to be occurred continuously for the specified length of time before activating the relay. If the condition disappears before the delay period is over, the relay will not be activated.

7.8 Clean

PATH:  \ CONFIG \ Clean

See the following explanation to get more details about the different settings for Clean

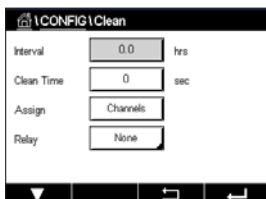

Enter the cleaning **Interval** time in hours. The cleaning interval can be set from 0.000 to 99999 hours. Setting it to 0 turns the clean cycle off.

Enter the **Clean Time** in seconds. The clean time can be 0 to 9999 seconds and must be smaller than the cleaning interval.


Assign the channel(s) for cleaning cycles. The assigned channels will be in HOLD state during the cleaning cycle.

Choose a **Relay**. Relay contacts are in normal mode until the cleaning cycle starts, then the relay is activated and the contact states change. Select Inverted to reverse the normal operating state of the relay (i.e. normally open contacts are in a open state, and normally closed contacts are in a closed state when the cleaning cycle starts).

7.9 Display Setup

PATH:  \ CONFIG \ Display Setup

See the following explanation to get more details about the different settings for Display Setup

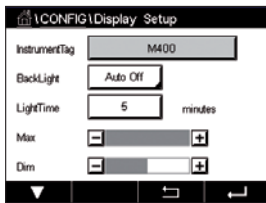

Enter the name for the M400 Transmitter (**Instrument Tag**). The instrument tag will also be displayed on the line at the top of the Start Screen and Menu Screen.

Use **BackLight** to switch off or dim the transmitter screen after a defined time period without interaction. The transmitter screen will automatically come back after pressing the display.

Enter the **Light Time** in minutes. The light time is the time period without interaction before the transmitter screen will be dimmed or switched off.

NOTE: In case of an unacknowledged warning or alarm the transmitter screen will not be dimmed or switched off even if the light time has been elapsed


The parameter **Max** allows the setting of the backlight during operation. With the parameter **Dim** the backlight of the transmitter screen during the dimmed state can be adjusted. Press the + or – buttons in the corresponding line to adjust the parameters.



7.10 Digital Inputs

PATH:  \ CONFIG \ Digital Inputs

See the following explanation to get more details about the different settings for the digital inputs

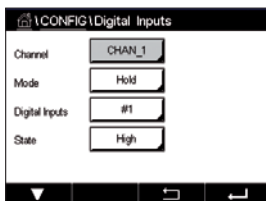

Press the input field in the line of the setting for **Mode** and select the impact of an active digital input signal. Choose 'HOLD' to lead the assigned channel in HOLD state.

Press the related button for the assignment of the **Digital Inputs** (#1 for DI1, #2 for DI2 etc.) and select the digital input signal, which has to be linked to the channel.


An additional setting can be done, if a digital input signal has been selected.

Press the input field in the line for the setting of the **State** and select if the digital input is active at high or low level of the voltage input signal.

7.11 System

PATH:  \ CONFIG \ System

See the following explanation to get more details about the different settings for the System.

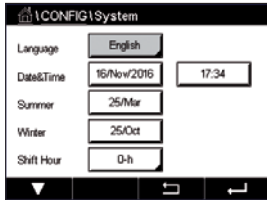

Select the desired **Language**. The following languages are available: English, French, German, Italian, Spanish, Portuguese, Russian, Chinese, Korean or Japanese.

Enter **Date&Time**.

The automatic change-over from summertime to wintertime and vice-versa frees the users from having to correct the time twice a year.

The winter to summer time-change is carried out automatically using the 12-month clock integrated in the transmitter. The date for the time-change can be set with the parameter **Summer**.

Provided it is a Sunday, the time-change would take place on the day that equates with the value, otherwise on the following Sunday. The winter/summer time-change takes place at 02:00 h.

The summer to winter time-change is carried out automatically using the 12-month clock integrated in the transmitter. The date for the time-change can be set through the parameter **Winter**.

Provided it is a Sunday, the time-change would take place on the day that equates with the value, otherwise on the following Sunday. The winter/summer time-change takes place at 03:00 h.

The number of hours, the clock will be shifted through the winter to summer and summer to winter time-change can be chosen. Press the related button for the setting of the **Shift Hour**.

7.12 PID Controller

PATH:  \ CONFIG \ PID Controller

PID control is proportional, integral and derivative control action that can provide smooth regulation of a process. Before configuring the transmitter, the following process characteristics must be identified.

Identify the **control direction** of the process

– **Conductivity:**

Dilution – direct acting where increasing measurement produces increasing control output such as controlling the feed of low conductivity diluting water to rinse tanks, cooling towers or boilers

Concentrating – reverse acting where increasing measurement produces decreasing control output, such as controlling chemical feed to attain a desired concentration

– **Dissolved Oxygen:**

Deaeration – direct acting where increasing Dissolved Oxygen concentration produces increasing control output such as controlling the feed of a reducing agent to remove oxygen from boiler feedwater

Aeration – reverse acting where increasing Dissolved Oxygen concentration produces decreasing control output, such as controlling an aerator blower speed to maintain a desired Dissolved Oxygen concentration in fermentation or wastewater treatment

– **pH/ORP:**

Acid feed only – direct acting where increasing pH produces increasing control output, also for ORP reducing reagent feed

Base feed only – reverse acting where increasing pH produces decreasing control output, also for ORP oxidizing reagent feed

Both acid and base feed – direct and reverse acting

Identify the **control output type** based on the control device to be used:

Pulse frequency – used with pulse input metering pump


Pulse length – used with solenoid valve

Analog – used with current input device such as electric drive unit, analog input metering pump or current-to-pneumatic (I/P) converter for pneumatic control valve


Default control settings provide linear control, which is appropriate for conductivity, dissolved oxygen. Therefore, when configuring PID for these parameters (or simple pH control) ignore settings of deadband and corner points in the tuning parameter section below. The non-linear control settings are used for more difficult pH/ORP control situations.

If desired, identify the non-linearity of the pH/ORP process. Improved control can be obtained if the non-linearity is accommodated with an opposing non-linearity in the controller. A titration curve (graph of pH or ORP vs. reagent volume) made on a process sample provides the best information. There is often a very high process gain or sensitivity near the setpoint and decreasing gain further away from the setpoint. To counteract this, the instrument allows for adjustable non-linear control with settings of a deadband around the setpoint, corner points further out and proportional limits at the ends of control as shown in the figure below.

Determine the appropriate settings for each of these control parameters based on the shape of the pH process titration curve.

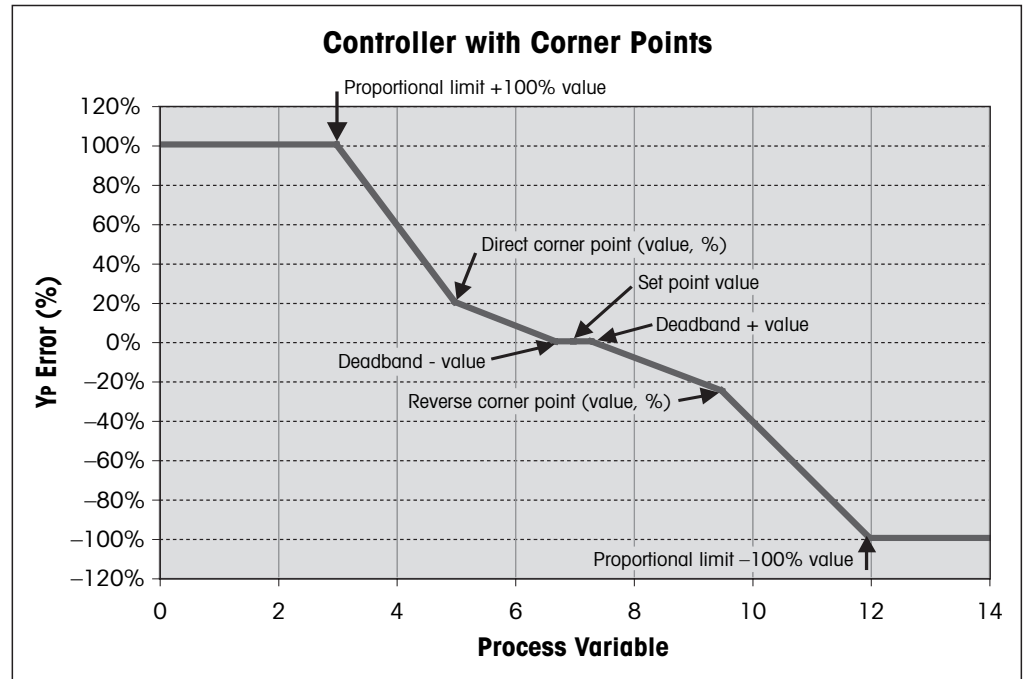

See the following explanation to get more details about the different settings for PID Controller.

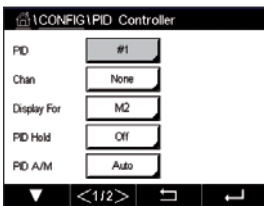

The M400 provides to one PID controller.

Press the related button for the assignment of the channel (**Chan**). Select the channel, which has to be linked to the PID Controller. To deactivate the PID controller press None.

Press the button for the assignment of the measuring parameter – based on the selected channel – that has been linked to the PID controller. Choose the measuring parameter by pressing the according field. Mx in the display indicates the measurement assigned to the PID Controller. (see chapter 7.1.1 "Channel Setup").

The M400 offers the display of control output (%PID) of the PID controller in the Start Screen and Menu Screen. Press the related button for **Display For** and select the line, the control output should be displayed by pressing the corresponding field.



NOTE: The control output of the PID controller will be displayed instead of the measurement, that has been defined to be shown in the corresponding line (see chapter "7.1.1 Channel Setup" on page 69).


Select with the parameter **PID HOLD** the state of the control output for the PID controller if the M400 Transmitter is in HOLD mode. Off means that the control output will be 0%PID if the transmitter is in HOLD mode. If Last Value has been chosen, the value for the control output signal before the transmitter went into HOLD mode will be used.

The parameter **PID A/M** allows selection of auto or manual operation for the PID controller. If auto has been chosen, the transmitter calculates the output signal based on the measured value and the settings of the parameters for the PID controller. In the case of manual operation, the transmitter shows in the Menu Screen at the line where the output signal is displayed two additional arrow buttons. Press the arrows buttons to increase or decrease the PID output signal.



NOTE: If Manual has been chosen the values for the time constants, gain, corner points, proportional limits, setpoint and deadband do not have any influence on the output signal.

Additional settings can be done by navigating to the next page of the menu.

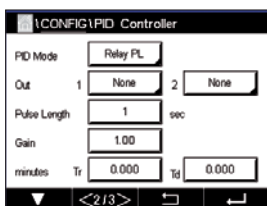

The **PID Mode** assigns a relay or analog output for PID control action. Based on the control device being used, select one of the three options Relay PL, Relay PF and Aout through pressing the corresponding field

- Relay PL: If using a solenoid valve, select Relays PL (Pulse Length).
- Relay PF: If using a pulse input metering pump, select Relays PF (Pulse Frequency)
- Aout: For using an analog control select Aout.

Link the output signal **Out1,2** of the PID controller to the desired output of the transmitter. Press the related button for Out 1 and Out 2 and select the corresponding number for the output through pressing the according field. #1 means relay 1 or Aout 1, #2 means relay 2 our Aout 2 etc.



NOTE: Take care if reed type relays are linked to the controlling function. The reed type relays could be used for pulse frequency control devices and light duty applications. The current is limited to 0.5 amps and 10 watts (see also chapter "14.2 Electrical specifications" on page 122). Do not connect to this relays higher current devices.

If the PID Mode is set to Relay PL, the Puls Length for the output signal of the transmitter can be adjusted. Press the button for **Pulse Length** and the M400 displays a keypad for modifying the value. Enter the new value in the unit seconds according to the table below and press **↵**.



NOTE: A longer pulse length will reduce wear on the solenoid valve. The % "on" time in the cycle is proportional to the control output.

| | 1 st Relay Position (Out 1) | 2 nd Relay Position (Out 2) | Pulse Length (PL) |
|------------------|--|--|--|
| Conductivity | Controlling concentrating reagent feed | Controlling dilution water | Short (PL) provides more uniform feed. Suggested start point = 30 sec |
| pH/ORP | Feeding base | Feeding acid | Reagent addition cycle: short PL provides more uniform addition of reagent. Suggested start point = 10 sec |
| Dissolved Oxygen | Reverse control action | Direct acting control action | Feed cycle time: short PL provides more uniform feed. Suggested start point = 30 sec |

If the PID Mode is set to Relay PF, the Pulse Frequency for the output signal of the transmitter can be adjusted. Press the button for **Pulse Freq** and enter the new value in the unit pulse / minute according to the table below.



NOTE: Set the pulse frequency to the maximum frequency allowed for the particular pump being used, typically 60 to 100 pulses/minute. Control action will produce this frequency at 100% output.



CAUTION: Setting the pulse frequency too high may cause the pump to overheat.

| | 1 st Relay Position = #3 | 2 nd Relay Position = #4 | Pulse Frequency (PF) |
|------------------|---|-------------------------------------|--|
| Conductivity | Controlling concentrating chemical feed | Controlling dilution water | Max allowed for the pump used (typically 60–100 pulses/minute) |
| pH/ORP | Feeding base | Feeding acid | Max allowed for the pump used (typically 60–100 pulses/minute) |
| Dissolved Oxygen | Reverse control action | Direct acting control action | Max allowed for the pump used (typically 60–100 pulses/minute) |

If the PID Mode is set to **Aout**, the type for the analog output signal of the transmitter can be selected. Press the corresponding button and choose between 4 to 20 mA and 0 to 20 mA for the output signal by pressing the according field.

For the assignment of the analog output signal consider the table below.

| | 1 st Analogout Position = Out 1 | 2 nd Analogout Position = Out 2 |
|------------------|--|--|
| Conductivity | Controlling concentrating chemical feed | Controlling dilution water |
| pH/ORP | Feeding base | Feeding acid |
| Dissolved Oxygen | Reverse control action | Direct acting control action |


Press the input field for the parameter **Gain** to enter the gain of the PID controller as a unitless value. Gain represents the maximum value of the output signal of the PID controller in per cent (value 1 corresponds to 100%).

Press the corresponding input field in the line of **min** to adjust the Parameter integral or reset time **Tr** (left button) and/or rate of derivate time **Td** (right button).



NOTE: Gain, integral and derivate time are usually adjusted later by trial and error on process response. It is recommended to start with the value $Td = 0$.

Further settings can be done by navigating to the next page of the menu.

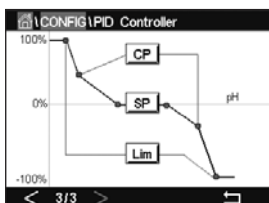

The display shows PID controller curve with input buttons for the corner points, setpoint and proportional limit for 100%.

Press the button **CP** to enter the menu for adjusting the corner points.

Page 1 shows the Corner Limit Low settings. Press the corresponding button to modify the value for the process parameter and the related output signal in %.

Browse to page 2 and the Corner Limit High settings are displayed. Press the corresponding button to modify the value for the process parameter and the related output signal in %.


Press the button **SP** to enter the menu for adjusting the setpoint and the dead band.

Press the button **Lim** to enter the menu for adjusting the proportional limit high and the proportional limit low, the range over which control action is required.

7.13 Service

PATH:  \ CONFIG \ Service

This menu is a valuable tool for troubleshooting and provides diagnostic functionality for the following items: Calibrate TouchPad, Set Analog Outputs, Read Analog Outputs, Read Analog Inputs, Set Relays, Read Relays, Read Digital Inputs, Memory and Display.

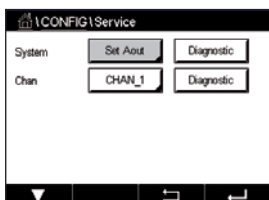

Select through the parameter **System** the desired item for diagnostic by pressing the according field.

Select through **Chan** the channel for diagnostic information of the sensor. This menu is only displayed if a sensor is connected.

The provided diagnostic functionality can now be called up through pressing the button **Diagnostic**.

7.13.1 Set Analog Outputs

The menu enables the user to set all analog outputs to any mA value within the 0–22 mA range. Use the + and – button to adjust the mA output signal. The transmitter will adjust the output signals according to the measurement and configuration of the analog output signals.

7.13.2 Read Analog Outputs

The menu shows the mA value of the analog outputs.

7.13.3 Set Relay

The menu allows the user to open or close each relay manually. If the menu is exited, the transmitter will switch the relay according to configuration.

7.13.4 Read Relay

The menu shows the state of every relay. On indicates the relay is closed, Off indicates that the relay is open.

7.13.5 Read Digital Inputs

The menu shows the state of the digital input signals.

7.13.6 Memory

If Memory is selected the transmitter will perform a memory test of all connected transmitter boards and ISM sensors.

7.13.7 Display

The transmitter shows every 5 seconds red, green, blue, grey and dark grey display and returns afterwards to the menu Service. If within the 5 seconds for every color the screen is pressed the transmitter will go to the next step.

7.13.8 Calibrate TouchPad

During the 4 calibrations steps, always press the center of the circle shown circle in the 4 corners of the display. The transmitter will show the calibration result.


7.13.9 Channel Diagnostic

If an error has occurred with the sensor, the corresponding messages are displayed.

7.14 User Management

PATH:  \ CONFIG \ User Management

This menu allows for the configuration of different user and administrator passwords, as well as setting up a list of allowed menus for the different users. The administrator has rights to access all menus. All default passwords for new transmitters are "00000000".

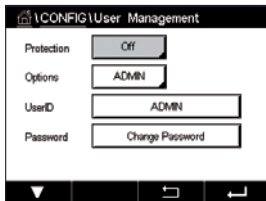

Press the input field in the line of **Protection** and select the desired kind of protection. The following options are available:

Off: No protection

Active: Activation of the Menu Screen (see chapter 3.3 "Display") has to be confirmed

Password: Activation of the Menu Screen is only possible with a password

Press the according button for **Option** to select the profile for the administrator (Admin) or one of the users.

NOTE: The administrator always has the rights to access all menus. For different users the access rights can be defined.

Press the input button for **UserID** to enter the name for the user or administrator. The name for the user or administrator will be displayed if the protection via password is selected for activation of the Menu Screen.

For changing the password of the selected user or administrator press the input field for **Password**. Enter the old password in the field Old PW, the new one in the field New PW and confirm it in the field confirm PW. The default password is "00000000" for the administrator and all users.

If the profile for a user has been selected an additional input field to define the access rights will be displayed.

To assign access rights the according button for the menu has to be pressed. In case of an assignment of the access rights, is displayed in the related button.



7.15 Reset

PATH:  \ CONFIG \ Reset

Depending on the transmitter version and configuration different options for a reset are available.

See the following explanation to get more details about the different option to reset data and / or configurations.

7.15.1 System Reset

This menu option allows the reset of the M400 Transmitter to the factory default settings (set-points off, analog outputs off, passwords, etc.). Furthermore the calibration factors for analog in- and outputs, meter etc. can be set to the last factory values.

Press the input field for **Options** and select System.


Press the input field for **Items** (Configure button) and select the different parts of the configuration that will be reset.

If an item has been selected the Action menu is displayed. Press the Reset button.

7.15.2 Reset Sensor Calibration for UniCond 2-e Sensors

For UniCond 2-e sensors, the SensorCal (sensor calibration) and ElecCal (sensor electronics calibration) can be restored to factory settings.

Press the input field for **Options** and select the channel the UniCond 2-e sensor is connected to.

Press the input field for **Item** (Configure button). Select SensorCal to Factory and/or ElecCal to Factory by checking the adjacent box. Press  enter to accept the value.


If an item has been selected the Action menu is displayed. Press the Reset button.

The M400 will bring up the confirmation dialog. Select Yes and the reset will be executed. Press No to go back to menu Reset without performing the reset.

7.16 USB

PATH:  \ CONFIG \ USB

This menu allows to output measurement values by a printer or to output measurement values for data log by USB communication.

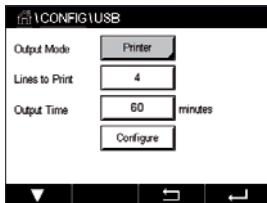


Select the Output Mode, Off or Printer or Data Log.

7.16.1 Printer Output Configuration

The Printer menu option allows configuring the M400 USB output to send data to a suitable printer. The printer output may be configured to print up to 4 configure measurements on separate lines, for each available sensor input, including pulsed input channels. At each print cycle, the output will include a header line with data and time based on the M400 internal clock, and one line for each configured measurement including channel, measurement descriptor, measurement value and unit of measure.

The output will appear as follows:


```
11/May/2012 15:36
Ch   Label Measurement
1    CHAN_1 4.01 pH
2    CHAN_1 25 centigrade
3    CHAN_1 200 DLI
```

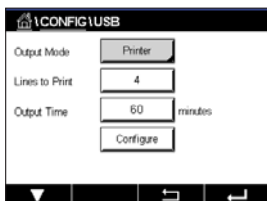

To configure the printer output, select option Printer for Output Mode. Configure the following options:

Lines to Print will configure the number of measurements that will be printed for each print cycle. Enter the total number of measurements to be configured for output. Lines to Print may be set from 1 to 4.

Output Time defines the time in minutes between each print cycle. Output time may be set from 1 to 1000 minutes.

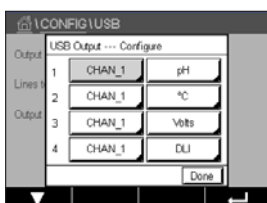


Once the output time and print lines have been established, press the Configure button to format the printer output. The number at the left of the window shows the order in which the lines will appear on the printer output. From the first dropdown, select the channel which the desired sensor is connected. This dropdown will list the labels associated with each channel as configured under Channel Setup. Using the second dropdown, select the unit associated with the measurement to be displayed.

7.16.2 USB data logging

The data log option allows configuring the M400 USB output to send data to a compatible USB memory stick. The data log may be configured to print up to 4 configure measurements on separate lines, for each available sensor input, including pulsed input channels. At each logging cycle, the output will include a header line with data and time based on the M400 internal clock, and one line for each configured measurement including channel, measurement descriptor, measurement value and unit of measure.

The output will appear as follows:

```
11/May/2012 15:36
Ch   Label Measurement
1    CHAN_1 4.01 pH
2    CHAN_1 25 °C
3    CHAN_1 200 DLI
```

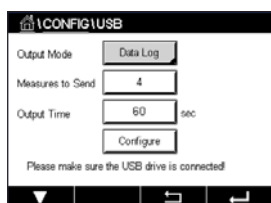

To configure the data log, select option Data log for Output Mode. Configure the following options:

Measures to Send will configure the number of measurements that will be sent for each print cycle.


Enter the total number of measurements to be configured for output. Lines to Print may be set from 1 to 4.

Output Time defines the time in minutes between each print cycle. Output time may be set from 1 to 1000 minutes.


Once the output time and print lines have been established, press the Configure button to format the data log. The number at the left of the window shows the order in which the lines will appear on the printer output. From the first dropdown, select the channel on which the desired sensor is connected. This dropdown will list the labels associated with each channel as configured under Channel Setup. Using the second dropdown, select the unit associated with the measurement to be displayed.




Start or Stop the data log with PATH: \ Config \ USB data logging after the data log is configured. Or you can setup the custom key to start or stop the data log. (See chapter "9 Custom Key" on page 109). The default setup of USB data logging is "Stop".

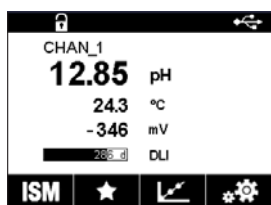


NOTE: Please make sure the USB memory stick is connected before starting the data log. A USB symbol will be displayed on the top of menu screen when a USB memory stick is connected. The compatible USB file system formats are FAT and FAT32.

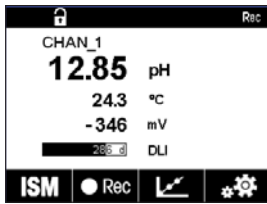

 **NOTE:** A "Rec" symbol will blink while the data is logging to a USB memory stick on the top of the menu screen.

7.17 Configuration via USB

PATH:  \ Config \ Configuration via USB

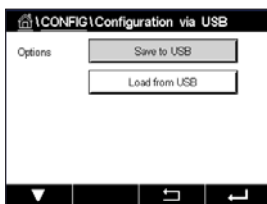

This menu allows the configuration of the current transmitter to be saved to a USB memory stick as a file or upload the configuration from a USB memory stick.



NOTE: The compatible USB file system formats are FAT and FAT32.

The naming of the configuration file must be MT_CFG_x. While x is 1 to 8. Please do not rename the saved configuration file. The configuration file saved by Transmitter Configuration Tool (TCT) can be used for uploading configuration to the transmitter.

Note: Configuration file cannot be used for M400 transmitter series across different transmitter series. e.g. M300 or M800.

8 ISM


For the menu structure refer to chapter 3.9 "Graphic Trend Measurement".

PATH:  \ ISM

8.1 iMonitor

PATH:  \ ISM \ iMonitor

The iMonitor gives an overview of the current state of the complete loop at a glance.

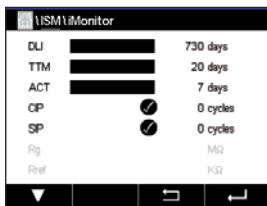

The iMonitor of the first channel is displayed on the screen. To browse through the iMonitor for the different channels press > at the bottom of the display.

The values DLI, TTM and ACT as well as TTCal in combination with UniCond 2-e sensors are shown as bar graph. If the values falls below 20% of the initial value the bar graph changes from green to yellow color. If the value falls below 10% the color changes to red.

For Cond4e sensors the days in operation of the sensor are displayed.

Furthermore SIP-, CIP-, AutoClave-, SAN-cycles as well as the values for Rg and Rref can be displayed and assigned to a colored button if the values are provided by the sensor.

The color for the related button of SIP-, CIP-, AutoClave- and SAN-cycles will change from green to yellow if less then 20% of the defined maximum quantity for the cycle remain and to red if less then 10% remain. For configuration of the maximum quantity see chapter "7.5 ISM Setup (ISM Sensors only)" on page 83.

The buttons for Rg and Rref change to yellow if the conditions for a warning messages are fulfilled and to red if the conditions for a alarm message are fulfilled. The buttons remain grey if the corresponding ISM alarm is not configured (see chapter "7.7 ISM / Sensor Alarm" on page 89).

Depending on the measured parameter (connected sensor) the following data are available in the menu iMonitor:


pH: DLI, TTM, ACT, CIP, AutoClave, SIP*, Rg**, Rref**
 Amperometric O₂: DLI, TTM, ACT, CIP, AutoClave, SIP*, Electrolyte***
 O₃: DLI, TTM, ACT, SAN
 Conductivity: Days in operation, TTCal****, CIP, SIP

- * if AutoClave has not been activated (see chapter "7.7 ISM / Sensor Alarm" on page 89)
- ** if the alarm for Rg and/or Rref has been activated (see chapter "7.7 ISM / Sensor Alarm" on page 89)
- *** if the alarm for Electrolyte Level Error has been activated (see chapter "7.7 ISM / Sensor Alarm" on page 89)
- **** if UniCond 2-e sensor is connected

8.2 Messages

PATH: \ ISM \ Messages

The messages for occurred warnings and alarms are listed in this menu. Up to 100 entries will be listed.

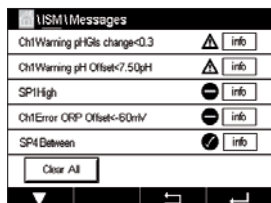

5 messages per page are listed. If more than 5 messages are available additional pages can be accessed.

Unacknowledged alarms or warning will be listed at the beginning. Then the acknowledged but still existing alarm or warning are listed. At the end of the list the already solved warning and alarms are described. Between these groups the messages are listed chronologically.

The state of the warning or alarm is indicated through the following signs:

| Symbol | Description | Meaning |
|--------|--------------------------------|--|
| | Alarm symbol is blinking | Alarm exists and has not been acknowledged |
| | Alarm symbol is not blinking | Alarm exists and has been acknowledged |
| | Warning symbol blinking | Warning exists and has not been acknowledged |
| | Warning symbol is not blinking | Warning exists and has been acknowledged |
| | OK symbol is not blinking | Warning or alarm has been solved |

An unacknowledged warning or alarm will be acknowledged by pressing the **Info** button in the corresponding line.

For every message the corresponding **Info** button can be pressed. Message information, date and time the warning or alarm has been occurred and the status of the alarm or message are displayed.

If warning or alarm has already been solved the pull up window for the message shows an additional button to clear the message i.e. to delete it from the message list.


8.3 ISM Diagnostics

PATH: \ ISM \ ISM Diagnostics

The M400 Transmitter provides for all ISM sensors a diagnostic menu. Access the menu Channel and select the channel by pressing the related input field.

Depending on the selected channel and assigned sensor different diagnostic menus are displayed. See the following explanation to get more details about the different diagnostic menus.

8.3.1 pH/ORP, Oxygen, O₃, Cond4e Sensors and TDL

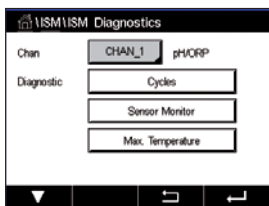

If an pH/ORP, oxygen, O₃ or Cond4e sensor is connected, the diagnostic menu cycles, sensor monitor and max. temperature are available.

Press the **Cycle** button and the information for CIP, SIP and Autoclave cycles of the connected sensor are displayed. The displayed information shows the amount of cycles the sensor has been exposed and the max. limitation for the corresponding cycle as defined in the menu ISM Setup (see chapter "7.5 ISM Setup (ISM Sensors only)" on page 83).

NOTE: For Cond4e, which are not autoclavable the menu AutoClave Cycles is not displayed.

NOTE: For O₃ sensors the SAN cycles are displayed.

NOTE: For TDL the cycles are not displayed.


Press the **Sensor Monitor** button and the information for DLI, TTM and ACT of the connected sensor are displayed. The values DLI, TTM and ACT are shown as bar graph. If the values falls below 20% of the initial value the bar graph changes from green to yellow color. If the value falls below 10% the color changes to red.

NOTE: For Cond4e sensors the operating hours are displayed.

Press the **Max. Temperature** button and the information about the maximum temperature, that the connected sensor has ever seen, together with a time stamp of this maximum is displayed. This value is stored on the sensor and cannot be changed. During autoclaving the max. temperature is not recorded.



8.3.2 UniCond 2-e and UniCond 4-e Sensors

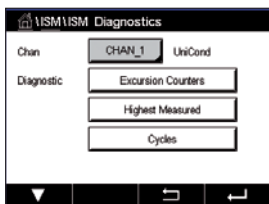

For UniCond 2-e and UniCond 4-e sensors, the following diagnostic items can be viewed: Excursion Counters including High Temp and High Conductivity, Highest Measured including Highest Temp and Highest Cond, Cycles including CIP cycles and SIP Cycles.


8.4 Calibration Data

PATH:  \ ISM \ Calibration Data

The M400 Transmitter provides a calibration history for all ISM sensors. Depending on the assigned sensor different data is available for the calibration history.

See the following explanation to get more details about the different data available for the calibration history.

8.4.1 Calibration Data for All ISM Sensors excluding UniCond 2-e and UniCond 4-e

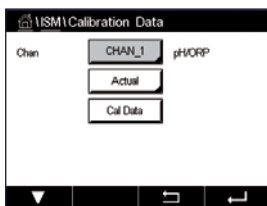

If an ISM sensor – excluding UniCond 2-e and UniCond 2-e – is connected to between the calibration data set of

Actual (Actual adjustment): This is the actual calibration dataset which is used for the measurement. This dataset moves to Cal1 position after the next adjustment.

Factory (Factory calibration): This is the original dataset, determined in the factory. This dataset remains stored in the sensor for reference and cannot be overwritten.

1.Adjust (First adjustment): This is the first adjustment after the factory calibration. This dataset remains stored in the sensor for reference and cannot be overwritten.

Cal1 (last calibration/adjustment): This is the last executed calibration/adjustment data set. This dataset moves to Cal2 and then to Cal3 when a new calibration/adjustment is performed. Afterwards, the dataset is not available anymore. Cal2 and Cal3 acting in the same way as Cal1.

Cal2 and **Cal3** can be chosen. For the selection of the calibration data set press the corresponding field.




NOTE: The amperometric oxygen sensor of THORNTON and the O₃ sensor do not provide the data set Cal1, Cal2, Cal3 and 1.Adjust.

Press the **Cal Data** button and the corresponding calibration data set is displayed. Furthermore the time stamp for the calibration and the User ID is listed.



NOTE: This function requires the correct setting of date and time during calibration and / or adjustment tasks (see chapter "7.11 System" on page 91).

8.4.2 Calibration Data for UniCond 2-e and UniCond 4-e Sensors

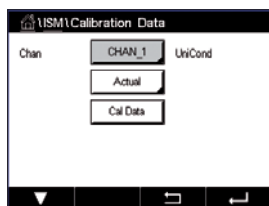

For UniCond 2-e and UniCond 4-e sensors the following three sets of calibration data may be selected:

Actual (Actual calibration): This is the actual calibration dataset which is used for the measurement.

Factory (Factory calibration): This is the original dataset, determined in the factory. This dataset remains stored in the sensor for reference and cannot be overwritten.

Cal1 (last calibration/adjustment): This is the last executed calibration/adjustment data set.

Press the Cal Data button and the corresponding calibration data set is displayed.

If the data set of the actual calibration has been chosen, on page 1, the date and time of the calibration, User ID, conductivity calibration constants, and reference conductivity values to calibrate are displayed. On page 2 the As-found conductivity values and the deviation from the reference are shown. On page 3 and 4 the same information for temperature is displayed. On page 5 the calibration cycles applied to the sensor and the next calibration date for conductivity (C) and temperature (T) are displayed.

If the dataset of the factory calibration has been chosen, on page 1, the date and time of the calibration, the conductivity calibration constants, and reference conductivity values used to calibrate are displayed. On page 2, the same values for temperature are shown.

Press ← to exit the menu Cal Data.



NOTE: This function requires the correct setting of date and time during calibration and / or adjustment tasks (see chapter "7.11 System" on page 91).


8.5 Sensor Info

PATH: \ ISM \ Sensor Info

The model, hardware and software version, last calibration date as well as the product and serial number of the ISM sensors, that are connected to the M400 Transmitter can be displayed on the screen.

Enter Sensor Info.

The data of the channel, a sensor is connected, are displayed on the screen.

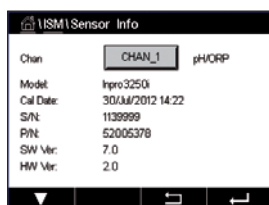

The data Model, Cal Date (date of last adjustment), S/N (serial number), P/N (product number), SW Ver (software version) and HW Ver (hardware version) of the select sensor are displayed.




NOTE: If a UniCond 2-e sensor is connected the following data is also displayed, Temp Sens. (temperature sensor) Electrode (electrode material), Body/Ins Mat: (body and/or insulator material), Inner: (inner electrode material), Outer (outer electrode material) Fitting: (fitting material), Class VI (FDA Class VI material).

To exit the menu Sensor Info press \leftarrow . To return to the Menu Screen press Home .

8.6 HW / SW Version

PATH: Home \ ISM \ HW/SW Version

The hardware and software version as well as the product number and serial number of the M400 Transmitter itself or the different boards, that are plugged in can be displayed on the screen.

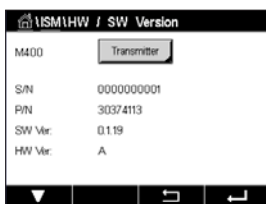

The data of the transmitter is displayed on the screen. Press the input field in the line of **M400**. To select the data of the desired board or the transmitter itself press the corresponding field.

The data S/N (serial number), P/N (product number), SW Ver (software version) and HW Ver (hardware version) of the select board or transmitter are displayed.

9 Custom Key


PATH:  \ Config \ Custom Key Setup



This menu allows the setting of a customize menu to the second left button on the menu screen as a shortcut. The custom key is a convenient option for soft key operation especially when touch-screen is not used.

Options: The "FAV" favorite is the default option. See chapter „Set Favorite“ for favorite setup.

- "Lock screen" can be selected for locking the screen
- "Trend" can be selected for graphic trend display
- "Messages" can be selected for the shortcut to access messages menu.
- "PID" can be selected for manual PID adjustment
- "Data log" can be selected for starting or stopping USB data logging.

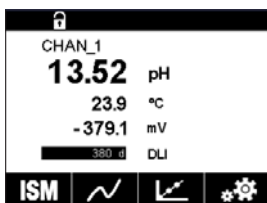

After the custom key setup, the selected custom key will be displayed at the second left button on the menu screen.




NOTE: The option "Data log" will be displayed only if "USB data logging" is selected. The option "PID" will be displayed only if manual PID controller is set.

9.1 Set Favorite

PATH:  \ FAVORITE \ Set Favorite

The M400 Transmitter allows set up of up to 4 favorites to ensure a quick access for frequently used functions.

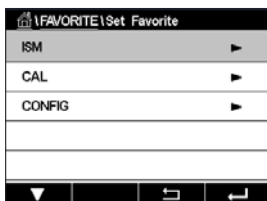

The main menus are displayed. Choose the menu, that contains the function, which should be defined as a favorite, e.g. ISM through pressing the corresponding arrow ► in the same line.

Choose the function, that should be set as a favorite by activating the option. A function, which is set as a favorite shows ★ icon.



NOTE: Deactivate the option by pressing on the icon again. The favorite ★ icon is not shown any more.

Access the menu Set Favorites. The favorites defined are listed on this page. Press the corresponding arrow ► for the function in the same line.

10 Maintenance

10.1 Front panel cleaning

Clean the surfaces with a soft damp cloth and dry the surfaces with a cloth carefully.

11 Software History

11.1 M400 Type 1

| Software version | Release date | Software changes | Documentation / Issue |
|------------------|--------------|------------------|---|
| V1.0.0 | March 2016 | – | 30 413 330 E M400 Transmitter 05/2018 |

11.2 M400 Type 2

| Software version | Release date | Software changes | Documentation / Issue |
|------------------|--------------|------------------|---|
| V1.0.0 | March 2016 | – | 30 413 330 E M400 Transmitter 05/2018 |

11.3 M400 Type 3

| Software version | Release date | Software changes | Documentation / Issue |
|------------------|--------------|------------------|---|
| V1.0.0 | March 2016 | – | 30 413 330 E M400 Transmitter 05/2018 |

11.4 M400 4-wire FF

| Software version | Release date | Software changes | Documentation / Issue |
|------------------|--------------|------------------|---|
| V1.0.0 | May 2018 | – | 30 413 330 E M400 Transmitter 05/2018 |

12 Troubleshooting

If the equipment is used in a manner not specified by Mettler-Toledo, the protection provided by the equipment may be void.

Review the table below for possible causes of common problems:

| Problem | Possible Cause |
|----------------------------------|--|
| Display is blank. | <ul style="list-style-type: none"> – No power to M400. – Hardware failure. |
| Incorrect measurement readings. | <ul style="list-style-type: none"> – Sensor improperly installed. – Incorrect units multiplier entered. – Temperature compensation incorrectly set or disabled. – Sensor or transmitter needs calibration. – Sensor or patch cord defective or exceeds recommended maximum length. – Hardware failure. |
| Measurement readings not stable. | <ul style="list-style-type: none"> – Sensors or cables installed too close to equipment that generates high level of electrical noise. – Recommended cable length exceeded. – Averaging set too low. – Sensor or patch cord defective. |
| Alarm symbol is shown. | <ul style="list-style-type: none"> – Setpoint is in alarm condition (setpoint exceeded). – Alarm has been selected (see chapter "7.7 ISM / Sensor Alarm" on page 89) and occurred. |
| Cannot change menu settings. | <ul style="list-style-type: none"> – User locked out for security reasons. |

12.1 Conductivity (resistive) Error messages/ Warning- and Alarm list for analog sensors

| Alarms | Description |
|--------------------|--|
| Watchdog time-out* | SW/System fault |
| Cond Cell open* | Cell running dry (no measurement solution) or wires are broken |
| Cond Cell shorted* | Short circuit caused by sensor or cable |

* Activate this function in the transmitter settings (see chapter 7.6 "General Alarm" PATH: Menu / General Alarm).

12.2 Conductivity (resistive) Error messages/ Warning- and Alarm list for ISM sensors

| Alarms | Description |
|--------------------|--|
| Watchdog time-out* | SW/System fault |
| Dry Cond sensor* | Cell running dry (no measurement solution) |
| Cell deviation* | Multiplier out of tolerance** (depends on sensor model). |

* Activate this function in the transmitter settings (see chapter "7.7 ISM / Sensor Alarm" on page 89

PATH: Menu/ISM/Sensor Alarm).

** For further information refer to the sensor documentation

12.3 pH Error messages/Warning- and Alarm list

12.3.1 pH, pH/pNa and Dissolved Carbon Dioxide sensors

| Warnings | Description |
|---------------------------------------|--|
| Warning pH Slope too high | Slope >102% |
| Warning pH Slope too low | Slope < 90% |
| Warning pH offset too high | pH ZeroPt > mmmptH |
| Warning pH offset too low | pH ZeroPt < nnnptH |
| Warning glass resistance low** | Glass electrode resistance changed by less than factor 0.3 |
| Warning glass resistance high** | Glass electrode resistance changed by more than factor 3 |
| Warning pH reference resistance low** | Reference electrode resistance changed by less than factor 0.3 |
| Warning reference resistance high** | Reference electrode resistance changed by more than factor 3 |

| Alarms | Description |
|--------------------------------------|--|
| Watchdog time-out* | SW/System fault |
| Error pH Slope too high | Slope >103% |
| Error pH Slope too low | Slope < 80% |
| Error pH offset too high | pH ZeroPt > xxxpH |
| Error pH offset too low | pH ZeroPt < yyypH |
| Error pH reference resistance high** | Reference electrode resistance >150 K Ω (break) |
| Error pH reference resistance low** | Reference electrode resistance < 1000 Ω (short) |
| Error pH glass resistance high ** | Glass electrode resistance > 2000 M Ω (break) |
| Error pH glass resistance low ** | Glass electrode resistance < 5 M Ω (short) |

* ISM sensors only

** Activate this function in the transmitter settings (see chapter "7.7 ISM / Sensor Alarm" on page 89 PATH: Menu/ISM/Sensor Alarm).

12.3.2 ORP messages

| Warnings* | Description |
|-----------------------------|-----------------------|
| Warning ORP ZeroPt > 30 mV | Zero offset too big |
| Warning ORP ZeroPt < -30 mV | Zero offset too small |

| Alarms* | Description |
|---------------------------|-----------------------|
| Watchdog time-out | SW/System fault |
| Error ORP ZeroPt > 60 mV | Zero offset too big |
| Error ORP ZeroPt < -60 mV | Zero offset too small |

* ISM sensors only

12.4 Amperometric O₂ Error messages/ Warning- and Alarm list

12.4.1 High level oxygen sensors

| Warnings | Description |
|---|-----------------------|
| Warning O ₂ Slope < -90 nA | Slope too big |
| Warning O ₂ Slope > -35 nA | Slope too small |
| Warning O ₂ ZeroPt > 0.3 nA | Zero offset too big |
| Warning O ₂ ZeroPt < -0.3 nA | Zero offset too small |

| Alarms | Description |
|---------------------------------------|------------------------------|
| Watchdog time-out* | SW/System fault |
| Error O ₂ Slope < -110 nA | Slope too big |
| Error O ₂ Slope > -30 nA | Slope too small |
| Error O ₂ ZeroPt > 0.6 nA | Zero offset too big |
| Error O ₂ ZeroPt < -0.6 nA | Zero offset too small |
| Electrolyte Low* | Too low level of electrolyte |

* ISM sensors only

12.4.2 Low level oxygen sensors

| Warnings | Description |
|---|-----------------------|
| Warning O ₂ Slope < -460 nA | Slope too big |
| Warning O ₂ Slope > -250 nA | Slope too small |
| Warning O ₂ ZeroPt > 0.5 nA | Zero offset too big |
| Warning O ₂ ZeroPt < -0.5 nA | Zero offset too small |

| Alarms | Description |
|---------------------------------------|---|
| Watchdog time-out* | SW/System fault |
| Error Install O ₂ Jumper | In case of using Hi Performance Oxygen a jumper has to be installed. See chapter "4.4.3 TB3 Terminal definition – Analog sensors" on page 32. |
| Error O ₂ Slope < -525 nA | Slope too big |
| Error O ₂ Slope > -220 nA | Slope too small |
| Error O ₂ ZeroPt > 1.0 nA | Zero offset too big |
| Error O ₂ ZeroPt < -1.0 nA | Zero offset too small |
| Electrolyte Low* | Too low level of electrolyte |

* ISM sensors only

12.4.3 Trace oxygen sensors

| Warnings | Description |
|---|-----------------------|
| Warning O ₂ Slope < -5000 nA | Slope too big |
| Warning O ₂ Slope > -3000 nA | Slope too small |
| Warning O ₂ ZeroPt > 0.5 nA | Zero offset too big |
| Warning O ₂ ZeroPt < -0.5 nA | Zero offset too small |

| Alarms | Description |
|---------------------------------------|------------------------------|
| Watchdog time-out | SW/System fault |
| Error O ₂ Slope < -6000 nA | Slope too big |
| Error O ₂ Slope > -2000 nA | Slope too small |
| Error O ₂ ZeroPt > 1.0 nA | Zero offset too big |
| Error O ₂ ZeroPt < -1.0 nA | Zero offset too small |
| Electrolyte Low* | Too low level of electrolyte |

* ISM sensors only

12.5 Warning- and Alarm Indication

12.5.1 Warning Indication



Warnings are indicated by a warning symbol in the head line of the display.


A warning message will be recorded and can be selected through the menu Messages (PATH: Home \ ISM \ Messages; see also chapter 8.2 "Messages").



NOTE: If the warning has not been acknowledged, the head line of the display will blink. If the warning has already been acknowledged, the head line will be displayed continuously. See also chapter 8.2 "Messages". In the case of an unacknowledged warning or alarm the transmitter screen will not be dimmed or switched off even if the light time has been elapsed (see chapter 7.9 "Display Setup").



NOTE: If at the same time a channel has both an alarm and a warning indicated, the indication of the alarm will have higher priority. The alarm will be indicated (see chapter 12.5 "Warning- and Alarm Indication") on the Menu Screen or Start Screen, while the warning will not be shown.

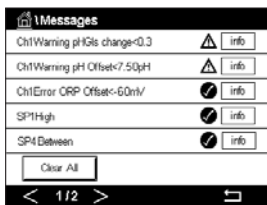

Pressing the head line on the Menu Screen will lead to the Messages. Refer to chapter 8.2 "Messages" for the description of the functionality for this menu.



NOTE: The detection of some warnings can be activated/deactivated through (de)activating the corresponding alarm. Refer to chapter 7.7 "ISM / Sensor Alarm".

12.5.2 Alarm Indication



Alarms are indicated by an alarm symbol in the head line of the display. An alarm message will be recorded and can be selected through the menu Messages (PATH: \ ISM \ Messages; see also chapter 8.2 "Messages").



NOTE: If the alarm has not been acknowledged, the head line of the display will blink. If the alarm has already been acknowledged, the head line will be displayed continuously. See also chapter 8.2 "Messages". In the case of an unacknowledged warning or alarm the transmitter screen will not be dimmed or switched off even if the light time has been elapsed (see chapter 7.9 "Display Setup").



NOTE: If at the same time a channel has born an alarm and a warning indicated, the indication of the alarm will have higher priority. The alarm will be indicated (see chapter 12.5 "Warning- and Alarm Indication") on the Menu Screen or Start Screen, while the warning will not be shown.



Pressing the head line on the Menu Screen will lead to the Messages. Refer to chapter 8.2 "Messages" for the description of the functionality for this menu.



NOTE: The detection of some alarms can be activated/deactivated. Refer therefore to chapter 7.7 "ISM / Sensor Alarm".



NOTE: Alarms which are caused by a violation of the limitation of a setpoint or the range (PATH: \ CONFIG \ Set Points; see also chapter 7.4 "Set Points") will also be indicated on the display and recorded through the menu Messages (PATH: \ ISM \ Messages; see also chapter 8.2 "Messages").

13 Ordering Information, Accessories and Spare Parts

Please contact your local Mettler-Toledo sales office or representative for details on additional accessories and spare parts.

| Transmitter | Order no. |
|--------------------|------------------|
| M400 Type 1 | 30 374 111 |
| M400 Type 2 | 30 374 112 |
| M400 Type 3 | 30 374 113 |
| M400 4-wire FF | 30 374 121 |

1) Included: 1 piece M25×1.5 cable gland, 4 pieces M20×1.5 cable glands

| Description | Order no. |
|--|------------------|
| Pipe mounting kit for ½ DIN for pipe diameter 40 to 60 mm (1.57 to 2.36") | 30 300 480 |
| Panel mounting kit for ½ DIN | 30 300 481 |
| Wall mounting kit for ½ DIN | 30 300 482 |
| Protective hood for ½ DIN | 30 073 328 |

14 Specifications

14.1 General specifications

pH/ORP (incl. pH/pNa)

| | |
|------------------------------------|--|
| Measurement parameters | pH, mV and temperature |
| pH display range | -2.00 to +16.00 pH |
| pH resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| pH accuracy ¹⁾ | Analog: ± 0.02 pH |
| mV range | -1500 to +1500 mV |
| mV resolution | Auto/0.001/0.01/0.1/1 mV (can be selected) |
| mV accuracy ¹⁾ | Analog: ± 1 mV |
| Temperature input ²⁾ | Pt 1000/Pt 100/NTC 22k |
| Temperature measuring range | -30 to +130 °C (-22 to +266 °F) |
| Temperature resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| Temperature accuracy ¹⁾ | Analog: ± 0.25 °C (± 0.45 °F) |
| Temperature compensation | Automatic/Manual |
| Max. sensor cable length | Analog: 10 to 20 m (33 to 65 ft) depending on sensor ISM: 80 m (260 ft) |
| Calibration | 1-point, 2-point or process |

1) ISM input signal causes no additional error.

2) Not required on ISM sensors

Amperometric oxygen

| | |
|------------------------------------|---|
| Measurement parameters | Dissolved oxygen (DO): Saturation or concentration and temperature Oxygen in gas: Concentration and temperature |
| Measuring current range | Analog: 0 to -7000 nA |
| Oxygen display ranges | <ul style="list-style-type: none"> • Dissolved Oxygen Saturation: 0 to 500 % air, 0 to 200 % O₂-sat Concentration: 0 ppb (µg/L) to 50.00 ppm (mg/L) • In gas Saturation: 0-100% O₂ gas Concentration: 0 to 9999 ppm O₂ gas |
| Oxygen accuracy ¹⁾ | <ul style="list-style-type: none"> • Dissolved Oxygen: Saturation ± 0.5 % of the measured value or ± 0.5 %, depending on which is larger. • Concentration at high values: ± 0.5 % of the measured value or ± 0.050 ppm/± 0.050 mg/L, depending on which is larger. • Concentration at low values: ± 0.5 % of the measured value or ± 0.001 ppm/± 0.001 mg/L, depending on which is larger • In gas: ± 0.5 % of the measured value or ± 5 ppb, depending on which is larger for ppm O₂ gas. • ± 0.5 % of the measured value or ± 0.01 %, depending on which is larger for vol-% O₂. |
| DO resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| Polarization voltage | <ul style="list-style-type: none"> • O₂ High: Cal/Meas: -675 mV (Configurable) • O₂ Low: Cal: -675 mV, Meas: -500 mV (Configurable) |
| Temperature input | Pt 1000/Pt 100/NTC 22k |
| Temperature compensation | Automatic |
| Temperature measuring range | -10 to +80 °C (+14 to +176 °F) |
| Temperature resolution | Auto/0.001/0.01/0.1/1 °C (°F) (can be selected) |
| Temperature accuracy ¹⁾ | ± 0.25 °C (± 0.45 °F) |
| Max. sensor cable length | <ul style="list-style-type: none"> • Analog: 20 m (65 ft) • ISM: 80 m (260 ft) |
| Calibration | 1-point (slope and offset) or process (slope and offset) |

1) ISM input signal causes no additional error.

Optical oxygen

| | |
|-----------------------------|---|
| Measurement parameters | Dissolved Oxygen (DO): Saturation or concentration and temperature Oxygen in gas: Concentration and temperature |
| Oxygen display ranges | <ul style="list-style-type: none"> • Dissolved Oxygen Saturation: 0 to 500 % air, 0 to 200 % O₂-sat <li style="padding-left: 20px;">Concentration: 0 ppb (ug/L) to 50.00 ppm (mg/L) • In gas Saturation: 0 to 100 vol-% O₂ <li style="padding-left: 20px;">Concentration: 0 to 9999 ppb O₂ gas |
| Oxygen accuracy | ± 1 digit |
| Oxygen resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| Temperature compensation | Automatic |
| Temperature measuring range | -30 to +150 °C (-22 to +302 °F) |
| Temperature resolution | Auto/0.001/0.01/0.1/1 °C (°F) (can be selected) |
| Temperature accuracy | ± 1 digit |
| Max. sensor cable length | 80 m (260 ft) |
| Calibration | 1-point (depending on sensor model) 2-point or process, process scaling |

Dissolved carbon dioxide

| | |
|-------------------------------|--|
| Measurement parameters | Dissolved carbon dioxide and temperature |
| CO ₂ display range | 0 to 5000 mg/L 0 to 200 % sat 0 to 1500 mm Hg 0 to 2000 mbar 0 to 2000 hPa |
| CO ₂ accuracy | ± 1 digit |
| CO ₂ resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| mV range | -1500 to +1500 mV |
| mV resolution | Auto/0.01/0.1/1 mV (can be selected) |
| mV accuracy | ± 1 digit |
| Total pressure range | 0 to 4000 mbar |
| Temperature measuring range | -30 to +150 °C (-22 to +302 °F) |
| Temperature resolution | Auto/0.001/0.01/0.1/1 °C (°F) (can be selected) |
| Temperature accuracy | ± 1 digit |
| Max. sensor cable length | 80 m (260 ft) |
| Calibration | 1-point (offset), 2-point (slope and offset) or process (offset) |

CO₂ hi (thermal conductivity)

| | |
|----------------------------------|---|
| Measurement parameters | Dissolved carbon dioxide and temperature |
| CO ₂ display ranges | 0 to 10 bar p (CO ₂)/0 to 145 psi p (CO ₂) 0 to 15 g/L 0 to 7 V/V CO ₂ |
| Accuracy in fluids ¹⁾ | ± 1 % of reading (within ± 5 % of calibration temperature) ± 2 % of reading over temperature range 0 to 50 °C (32 to 122 °F) |
| Calibration | 1-point or process |

1) Complete loop of sensor and transmitter

GPro 500 TDL

| | |
|---|--|
| Measurement parameters | O ₂ , O ₂ and temperature, CO (ppm), CO (%) , H ₂ O, CO ₂ (%) |
| Gas display ranges | 0 to 100 % |
| Gas accuracy, resolution, repeatability and low detection limit | Depending on sensor model |
| Linearity | Better than 1 % |
| Drift | Negligible (< 2 % of measurement range between maintenance intervals) |
| Sampling rate | 1 second |
| Response time (t ₉₀) | Depending on sensor model |
| Process pressure ranges | Depending on sensor model |
| Process temperature ranges | 0 to 250 °C (32 to 482 °F) optional (for probe installation) 0 to 600 °C (32 to 1112 °F) with additional thermal barrier 0 to 150 °C (32 to 302 °F) (white cell) |
| Max. sensor cable length | 40 m (130 ft) (FM version) |
| Calibration | 1-point (offset) or process (slope or offset) |

Dissolved ozone

| | |
|------------------------------------|--|
| Measurement parameters | Concentration and temperature |
| Display range for current | Analog: 0 to -7000 nA |
| Ozone measuring range | 0 to 5000 ppb (µg/L) O ₃ |
| Ozone accuracy | ± 1 % (or 0.4 ppb) up to 2000 ppb ± 2.5 % (or 50 – 125 ppb) from 2000 to 5000 ppb |
| Resolution | ± 1 digit |
| Temperature compensation | Automatic |
| Temperature measuring range | 5 to +50 °C (+41 to +122 °F) |
| Temperature resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| Temperature accuracy ¹⁾ | Analog: ± 0.25 °C (± 0.45 °F) |
| Max. sensor cable length | 80 m |
| Calibration | 1-point (offset) or process (slope and offset) |

Conductivity 2-e/4-e

| | |
|---|--|
| Measurement parameters | Conductivity/resistivity and temperature |
| Conductivity ranges | See sensor specification |
| Chemical concentration curves (used with 4-e sensors) | NaCl: 0–26 % @ 0 °C to 0–28 % @ +100 °C NaOH: 0–12 % @ 0 °C to 0–16 % @ +40 °C to 0–6 % @ +100 °C HCl: 0–18 % @ -20 °C to 0–18 % @ 0 °C to 0–5 % @ +50 °C HNO ₃ : 0–30 % @ -20 °C to 0–30 % @ 0 °C to 0–8 % @ +50 °C H ₂ SO ₄ : 0–26 % @ -12 °C to 0–26 % @ +5 °C to 0–9 % @ +100 °C H ₃ PO ₄ : 0–35 % @ +5 °C to +80 °C |
| TDS ranges | NaCl, CaCO ₃ |
| Cond/Res accuracy ¹⁾ | Analog: ± 0.5 % of reading or 0.25 Ω, whichever is greater |
| Cond/Res repeatability ¹⁾ | Analog: ± 0.25 % of reading or 0.25 Ω, whichever is greater |
| Cond/Res resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| Temperature input | Pt 1000 |
| Temperature measuring range | -40 to +200 °C (-40 to +392 °F) |
| Temperature resolution | Auto/0.001/0.01/0.1/1 (can be selected) |
| Temperature accuracy | Analog: ± 0.25 °C (± 0.45 °F) within -30 to +150 °C (-22 to +302 °F); ± 0.50 °C (± 0.90 °F) outside |
| Max. sensor cable length | • Analog: 2-e sensors: 61 m (200 ft); 4-e sensors: 15 m (50 ft) • ISM: 2-e sensors: 90 m (300 ft); 4-e sensors: 80 m (260 ft) |
| Calibration | 1-point, 2-point or process |

1) ISM input signal causes no additional error.

14.2 Electrical specifications

| | |
|--|---|
| Supply voltage | <ul style="list-style-type: none"> • 80 to 255 VAC, 50 to 60 Hz, 10 VA • 20 to 30 VDC, 10 VA |
| Connection terminal | Detachable screw terminals, appropriate for wire cross section 0.2 to 1.5 mm ² (AWG 16–24) |
| Mains fuse | 2.0 A slow blow, type FC |
| Analog output ¹⁾ | 4 × 0/4 to 20 mA, 22 mA alarm, galvanically isolated from input and from earth/ground |
| Measurement error through analog outputs | < ±0.05 mA over 1 to 22 mA range |
| Analog output configuration | Linear, bi-linear, logarithmic, auto range |
| Load | Max. 500 Ω |
| PID process controller | 1 × PID with pulse length, pulse frequency or analog control output signal |
| Cycle time analog output | Ca. 1 s |
| Hold input/Alarm contact | Yes/Yes |
| Alarm output delay | 0 to 999 s, selectable |
| Relays | <ul style="list-style-type: none"> • 2 SPDT, mechanical, 250 VAC or 30 VDC, 3 A • 2 SPST, Reed, 250 VAC or 250 VDC, 0.5 A, 10 W |
| Digital input | 2 With switching limits 0.00 VDC to 1.00 VDC inactive, 2.30 VDC to 30.00 VDC active; galvanically isolated up to 60 V from output, analog input and ground/earth |
| Analog input ²⁾ | 1 × 0/4 to 20 mA |
| User interface | <ul style="list-style-type: none"> • TFT touch-screen 4" • Black and white • Resolution: ¼ VGA (320 pixel × 240 pixel) |
| Keypad | • 4 tactile feedback keys |
| Languages | 10 (English, German, French, Italian, Spanish, Portuguese, Russian, Japanese, Korean and Chinese) |
| Interfaces | <ul style="list-style-type: none"> • 1 USB Host: Printer connection, data logging, loading configuration from USB stick and saving configuration to USB stick • 1 USB Device: Software update interface |

1) For M400 Type 1, Type 2, Type 3 only.

2) For M400 Type 2, Type 3 and M400 4-wire FF only.

14.3 FOUNDATION fieldbus specifications

| | |
|-------------------------------------|--|
| Supply voltage for FF block | 9 to 32 V DC |
| Current | 22 mA |
| Max. current in case of fault (FDE) | < 28 mA |
| Physical interface | According to IEC 61158-2 |
| Transfer rate | 31.25 kbit / s |
| Profile | FF_H1 (Foundation fieldbus) |
| Communication protocol | FF-816 |
| ITK version | 6.1.0 |
| Manufacturer ID | (DEV_TYPE) 0x465255 |
| FF Type | (DEV_REV) 1 |
| FF communication model | <ul style="list-style-type: none"> • 1 Resource Block • 2 Transducer Blocks • 4 Analog Input Blocks • 1 Analog Output Block • 2 Discrete Input Blocks • 2 Discrete Output Blocks |

14.4 Environmental specifications

| | |
|---|---|
| Storage temperature | -40 to +70 °C (-40 to +158 °F) |
| Ambient temperature \neq operating range | -20 to +50 °C (-4 to +122 °F) |
| Relative humidity | 0 to 95 % non-condensing |
| Altitude | Max. 2000 m |
| EMC | Compliant with EN 61326-1:2013 (Industrial environment) Emission: Class A, Immunity: Class A |
| UL | Installation (overvoltage) Category II |
| CE mark | The measuring system is in conformity with the statutory requirements of the EC Directives. METTLER TOLEDO confirms successful testing of the device by affixing to it the CE mark. |
| Ex Approvals ¹⁾ | <ul style="list-style-type: none"> • cCSAus Class I, Division 2, Groups A, B, C, D T4 Class I, Zone 2, AEx nA nC IIC T4 Gc • ATEX II 3G Ex nA nC IIC T4 Gc • IECEx Ex ec ic nC IIC T4 Gc |

1) For M400 Type 1, Type 2, Type 3 only.

14.5 Mechanical specifications

| | | |
|------------------|----------------------------------|---|
| Dimensions | Housing – Height × Width × Depth | 136 × 136 × 116 mm (5.35 × 5.35 × 4.57") |
| | Front bezel – Height x Width | 150 × 150 mm (5.91 × 5.91") |
| | Max. depth – panel mounted | 116 mm (4.57") (excludes plug-in connectors) |
| Weight | 1.50 kg (3.3 lb) | |
| Material | Aluminum (ADC12) die cast | |
| Enclosure rating | IP 66/NEMA 4X | |

15 Warranty

METTLER TOLEDO warrants this product to be free from significant deviations in material and workmanship for a period of one year from the date of purchase. If repair is necessary and not the result of abuse or misuse within the warranty period, please return by pre-paid freight and an amendment will be made without any charge. METTLER TOLEDO's Customer Service Dept. will determine if the product problem is due to deviations or customer abuse. Out-of-warranty products will be repaired on an exchange basis at cost.

The above warranty is the only warranty made by METTLER TOLEDO and is in lieu of all other warranties, expressed or implied, including, without limitation, implied warranties of merchantability and fitness for a particular purpose. METTLER TOLEDO shall not be liable for any loss, claim, expense or damage caused by, contributed to or arising out of the acts or omissions of the

Buyer or Third Parties, whether negligent or otherwise. In no event shall METTLER TOLEDO's liability for any cause of action whatsoever exceed the cost of the item giving rise to the claim, whether based in contract, warranty, indemnity, or tort (including negligence).

16 Buffer tables

M400 Transmitters have the ability to do automatic pH buffer recognition. The following tables show different buffers that are automatically recognized.

16.1 Standard pH buffers

16.1.1 Mettler-9

| Temp (°C) | pH of buffer solutions | | | |
|-----------|------------------------|------|------|------|
| 0 | 2.03 | 4.01 | 7.12 | 9.52 |
| 5 | 2.02 | 4.01 | 7.09 | 9.45 |
| 10 | 2.01 | 4.00 | 7.06 | 9.38 |
| 15 | 2.00 | 4.00 | 7.04 | 9.32 |
| 20 | 2.00 | 4.00 | 7.02 | 9.26 |
| 25 | 2.00 | 4.01 | 7.00 | 9.21 |
| 30 | 1.99 | 4.01 | 6.99 | 9.16 |
| 35 | 1.99 | 4.02 | 6.98 | 9.11 |
| 40 | 1.98 | 4.03 | 6.97 | 9.06 |
| 45 | 1.98 | 4.04 | 6.97 | 9.03 |
| 50 | 1.98 | 4.06 | 6.97 | 8.99 |
| 55 | 1.98 | 4.08 | 6.98 | 8.96 |
| 60 | 1.98 | 4.10 | 6.98 | 8.93 |
| 65 | 1.98 | 4.13 | 6.99 | 8.90 |
| 70 | 1.99 | 4.16 | 7.00 | 8.88 |
| 75 | 1.99 | 4.19 | 7.02 | 8.85 |
| 80 | 2.00 | 4.22 | 7.04 | 8.83 |
| 85 | 2.00 | 4.26 | 7.06 | 8.81 |
| 90 | 2.00 | 4.30 | 7.09 | 8.79 |
| 95 | 2.00 | 4.35 | 7.12 | 8.77 |

16.1.2 Mettler-10

| Temp (°C) | pH of buffer solutions | | | | |
|-----------|------------------------|------|------|-------|--|
| 0 | 2.03 | 4.01 | 7.12 | 10.65 | |
| 5 | 2.02 | 4.01 | 7.09 | 10.52 | |
| 10 | 2.01 | 4.00 | 7.06 | 10.39 | |
| 15 | 2.00 | 4.00 | 7.04 | 10.26 | |
| 20 | 2.00 | 4.00 | 7.02 | 10.13 | |
| 25 | 2.00 | 4.01 | 7.00 | 10.00 | |
| 30 | 1.99 | 4.01 | 6.99 | 9.87 | |
| 35 | 1.99 | 4.02 | 6.98 | 9.74 | |
| 40 | 1.98 | 4.03 | 6.97 | 9.61 | |
| 45 | 1.98 | 4.04 | 6.97 | 9.48 | |
| 50 | 1.98 | 4.06 | 6.97 | 9.35 | |
| 55 | 1.98 | 4.08 | 6.98 | | |
| 60 | 1.98 | 4.10 | 6.98 | | |
| 65 | 1.99 | 4.13 | 6.99 | | |
| 70 | 1.98 | 4.16 | 7.00 | | |
| 75 | 1.99 | 4.19 | 7.02 | | |
| 80 | 2.00 | 4.22 | 7.04 | | |
| 85 | 2.00 | 4.26 | 7.06 | | |
| 90 | 2.00 | 4.30 | 7.09 | | |
| 95 | 2.00 | 4.35 | 7.12 | | |

16.1.3 NIST Technical Buffers

| Temp (°C) | pH of buffer solutions | | | | |
|-----------|------------------------|-------|-------|-------|-------|
| 0 | 1.67 | 4.00 | 7.115 | 10.32 | 13.42 |
| 5 | 1.67 | 4.00 | 7.085 | 10.25 | 13.21 |
| 10 | 1.67 | 4.00 | 7.06 | 10.18 | 13.01 |
| 15 | 1.67 | 4.00 | 7.04 | 10.12 | 12.80 |
| 20 | 1.675 | 4.00 | 7.015 | 10.07 | 12.64 |
| 25 | 1.68 | 4.005 | 7.00 | 10.01 | 12.46 |
| 30 | 1.68 | 4.015 | 6.985 | 9.97 | 12.30 |
| 35 | 1.69 | 4.025 | 6.98 | 9.93 | 12.13 |
| 40 | 1.69 | 4.03 | 6.975 | 9.89 | 11.99 |
| 45 | 1.70 | 4.045 | 6.975 | 9.86 | 11.84 |
| 50 | 1.705 | 4.06 | 6.97 | 9.83 | 11.71 |
| 55 | 1.715 | 4.075 | 6.97 | | 11.57 |
| 60 | 1.72 | 4.085 | 6.97 | | 11.45 |
| 65 | 1.73 | 4.10 | 6.98 | | |
| 70 | 1.74 | 4.13 | 6.99 | | |
| 75 | 1.75 | 4.14 | 7.01 | | |
| 80 | 1.765 | 4.16 | 7.03 | | |
| 85 | 1.78 | 4.18 | 7.05 | | |
| 90 | 1.79 | 4.21 | 7.08 | | |
| 95 | 1.805 | 4.23 | 7.11 | | |

16.1.4 NIST standard buffers (DIN and JIS 19266: 2000–01)

| Temp (°C) | pH of buffer solutions | | | |
|-----------|------------------------|-------|-------|-------|
| 0 | | | | |
| 5 | 1.668 | 4.004 | 6.950 | 9.392 |
| 10 | 1.670 | 4.001 | 6.922 | 9.331 |
| 15 | 1.672 | 4.001 | 6.900 | 9.277 |
| 20 | 1.676 | 4.003 | 6.880 | 9.228 |
| 25 | 1.680 | 4.008 | 6.865 | 9.184 |
| 30 | 1.685 | 4.015 | 6.853 | 9.144 |
| 37 | 1.694 | 4.028 | 6.841 | 9.095 |
| 40 | 1.697 | 4.036 | 6.837 | 9.076 |
| 45 | 1.704 | 4.049 | 6.834 | 9.046 |
| 50 | 1.712 | 4.064 | 6.833 | 9.018 |
| 55 | 1.715 | 4.075 | 6.834 | 8.985 |
| 60 | 1.723 | 4.091 | 6.836 | 8.962 |
| 70 | 1.743 | 4.126 | 6.845 | 8.921 |
| 80 | 1.766 | 4.164 | 6.859 | 8.885 |
| 90 | 1.792 | 4.205 | 6.877 | 8.850 |
| 95 | 1.806 | 4.227 | 6.886 | 8.833 |



NOTE: The pH(S) values of the individual charges of the secondary reference materials are documented in a certificate of an accredited laboratory. This certificate is supplied with the respective buffer materials. Only these pH(S) values shall be used as standard values for the secondary reference buffer materials. Correspondingly, this standard does not include a table with standard pH values for practical use. The table above only provides examples of pH(PS) values for orientation.

16.1.5 Hach buffers

Buffer values up to 60 °C as specified by Bergmann & Beving Process AB.

| Temp (°C) | pH of buffer solutions | | |
|-----------|------------------------|------|-------|
| 0 | 4.00 | 7.14 | 10.30 |
| 5 | 4.00 | 60 | 10.23 |
| 10 | 4.00 | 7.04 | 10.11 |
| 15 | 4.00 | 7.04 | 10.11 |
| 20 | 4.00 | 7.02 | 10.05 |
| 25 | 4.01 | 7.00 | 10.00 |
| 30 | 4.01 | 6.99 | 9.96 |
| 35 | 4.02 | 6.98 | 9.92 |
| 40 | 4.03 | 6.98 | 9.88 |
| 45 | 4.05 | 6.98 | 9.85 |
| 50 | 4.06 | 6.98 | 9.82 |
| 55 | 4.07 | 6.98 | 9.79 |
| 60 | 4.09 | 6.99 | 9.76 |

16.1.6 Ciba (94) buffers

| Temp (°C) | pH of buffer solutions | | | | |
|-----------|------------------------|-------|-------|-------|--|
| 0 | 2.04 | 4.00 | 7.10 | 10.30 | |
| 5 | 2.09 | 4.02 | 7.08 | 10.21 | |
| 10 | 2.07 | 4.00 | 7.05 | 10.14 | |
| 15 | 2.08 | 4.00 | 7.02 | 10.06 | |
| 20 | 2.09 | 4.01 | 6.98 | 9.99 | |
| 25 | 2.08 | 4.02 | 6.98 | 9.95 | |
| 30 | 2.06 | 4.00 | 6.96 | 9.89 | |
| 35 | 2.06 | 4.01 | 6.95 | 9.85 | |
| 40 | 2.07 | 4.02 | 6.94 | 9.81 | |
| 45 | 2.06 | 4.03 | 6.93 | 9.77 | |
| 50 | 2.06 | 4.04 | 6.93 | 9.73 | |
| 55 | 2.05 | 4.05 | 6.91 | 9.68 | |
| 60 | 2.08 | 4.10 | 6.93 | 9.66 | |
| 65 | 2.07* | 4.10* | 6.92* | 9.61* | |
| 70 | 2.07 | 4.11 | 6.92 | 9.57 | |
| 75 | 2.04* | 4.13* | 6.92* | 9.54* | |
| 80 | 2.02 | 4.15 | 6.93 | 9.52 | |
| 85 | 2.03* | 4.17* | 6.95* | 9.47* | |
| 90 | 2.04 | 4.20 | 6.97 | 9.43 | |
| 95 | 2.05* | 4.22* | 6.99* | 9.38* | |

* Extrapolated

16.1.7 Merck Titrisole, Riedel-de-Haën Fixanale

| Temp (°C) | pH of buffer solutions | | | | |
|-----------|------------------------|------|------|------|-------|
| 0 | 2.01 | 4.05 | 7.13 | 9.24 | 12.58 |
| 5 | 2.01 | 4.05 | 7.07 | 9.16 | 12.41 |
| 10 | 2.01 | 4.02 | 7.05 | 9.11 | 12.26 |
| 15 | 2.00 | 4.01 | 7.02 | 9.05 | 12.10 |
| 20 | 2.00 | 4.00 | 7.00 | 9.00 | 12.00 |
| 25 | 2.00 | 4.01 | 6.98 | 8.95 | 11.88 |
| 30 | 2.00 | 4.01 | 6.98 | 8.91 | 11.72 |
| 35 | 2.00 | 4.01 | 6.96 | 8.88 | 11.67 |
| 40 | 2.00 | 4.01 | 6.95 | 8.85 | 11.54 |
| 45 | 2.00 | 4.01 | 6.95 | 8.82 | 11.44 |
| 50 | 2.00 | 4.00 | 6.95 | 8.79 | 11.33 |
| 55 | 2.00 | 4.00 | 6.95 | 8.76 | 11.19 |
| 60 | 2.00 | 4.00 | 6.96 | 8.73 | 11.04 |
| 65 | 2.00 | 4.00 | 6.96 | 8.72 | 10.97 |
| 70 | 2.01 | 4.00 | 6.96 | 8.70 | 10.90 |
| 75 | 2.01 | 4.00 | 6.96 | 8.68 | 10.80 |
| 80 | 2.01 | 4.00 | 6.97 | 8.66 | 10.70 |
| 85 | 2.01 | 4.00 | 6.98 | 8.65 | 10.59 |
| 90 | 2.01 | 4.00 | 7.00 | 8.64 | 10.48 |
| 95 | 2.01 | 4.00 | 7.02 | 8.64 | 10.37 |

16.1.8 WTW buffers

| Temp (°C) | pH of buffer solutions | | | |
|-----------|------------------------|------|------|-------|
| 0 | 2.03 | 4.01 | 7.12 | 10.65 |
| 5 | 2.02 | 4.01 | 7.09 | 10.52 |
| 10 | 2.01 | 4.00 | 7.06 | 10.39 |
| 15 | 2.00 | 4.00 | 7.04 | 10.26 |
| 20 | 2.00 | 4.00 | 7.02 | 10.13 |
| 25 | 2.00 | 4.01 | 7.00 | 10.00 |
| 30 | 1.99 | 4.01 | 6.99 | 9.87 |
| 35 | 1.99 | 4.02 | 6.98 | 9.74 |
| 40 | 1.98 | 4.03 | 6.97 | 9.61 |
| 45 | 1.98 | 4.04 | 6.97 | 9.48 |
| 50 | 1.98 | 4.06 | 6.97 | 9.35 |
| 55 | 1.98 | 4.08 | 6.98 | |
| 60 | 1.98 | 4.10 | 6.98 | |
| 65 | 1.99 | 4.13 | 6.99 | |
| 70 | | 4.16 | 7.00 | |
| 75 | | 4.19 | 7.02 | |
| 80 | | 4.22 | 7.04 | |
| 85 | | 4.26 | 7.06 | |
| 90 | | 4.30 | 7.09 | |
| 95 | | 4.35 | 7.12 | |

16.1.9 JIS Z 8802 buffers

| Temp (°C) | pH of buffer solutions | | | |
|-----------|------------------------|-------|-------|-------|
| 0 | 1.666 | 4.003 | 6.984 | 9.464 |
| 5 | 1.668 | 3.999 | 6.951 | 9.395 |
| 10 | 1.670 | 3.998 | 6.923 | 9.332 |
| 15 | 1.672 | 3.999 | 6.900 | 9.276 |
| 20 | 1.675 | 4.002 | 6.881 | 9.225 |
| 25 | 1.679 | 4.008 | 6.865 | 9.180 |
| 30 | 1.683 | 4.015 | 6.853 | 9.139 |
| 35 | 1.688 | 4.024 | 6.844 | 9.102 |
| 38 | 1.691 | 4.030 | 6.840 | 9.081 |
| 40 | 1.694 | 4.035 | 6.838 | 9.068 |
| 45 | 1.700 | 4.047 | 6.834 | 9.038 |
| 50 | 1.707 | 4.060 | 6.833 | 9.011 |
| 55 | 1.715 | 4.075 | 6.834 | 8.985 |
| 60 | 1.723 | 4.091 | 6.836 | 8.962 |
| 70 | 1.743 | 4.126 | 6.845 | 8.921 |
| 80 | 1.766 | 4.164 | 6.859 | 8.885 |
| 90 | 1.792 | 4.205 | 6.877 | 8.850 |
| 95 | 1.806 | 4.227 | 6.886 | 8.833 |

16.2 Dual membrane pH electrode buffers

16.2.1 Mettler-pH/pNa buffers (Na+ 3.9M)

| Temp (°C) | pH of buffer solutions | | | |
|-----------|------------------------|------|------|------|
| 0 | 1.98 | 3.99 | 7.01 | 9.51 |
| 5 | 1.98 | 3.99 | 7.00 | 9.43 |
| 10 | 1.99 | 3.99 | 7.00 | 9.36 |
| 15 | 1.99 | 3.99 | 6.99 | 9.30 |
| 20 | 1.99 | 4.00 | 7.00 | 9.25 |
| 25 | 2.00 | 4.01 | 7.00 | 9.21 |
| 30 | 2.00 | 4.02 | 7.01 | 9.18 |
| 35 | 2.01 | 4.04 | 7.01 | 9.15 |
| 40 | 2.01 | 4.05 | 7.02 | 9.12 |
| 45 | 2.02 | 4.07 | 7.03 | 9.11 |
| 50 | 2.02 | 4.09 | 7.04 | 9.10 |

For addresses of METTLER TOLEDO
Market Organizations please go to:
www.mt.com/pro-MOs



Management System
certified according to
ISO 9001 / ISO 14001



Subject to technical changes.
© 07/2019 Mettler-Toledo GmbH, Process Analytics
Printed in Switzerland. 30 413 330 E

Mettler-Toledo GmbH, Process Analytics
Im Hackacker 15, CH-8902 Urdorf, Switzerland
Tel. +41 44 729 62 11, Fax +41 44 729 66 36

www.mt.com/pro